


Intrada Organic Acid

150 x 2 mm

Application

LC-MS: Organic acids in Energy drink (Brand-M)

エナジードリンク(Brand-M)中の有機酸のLC-MS分析

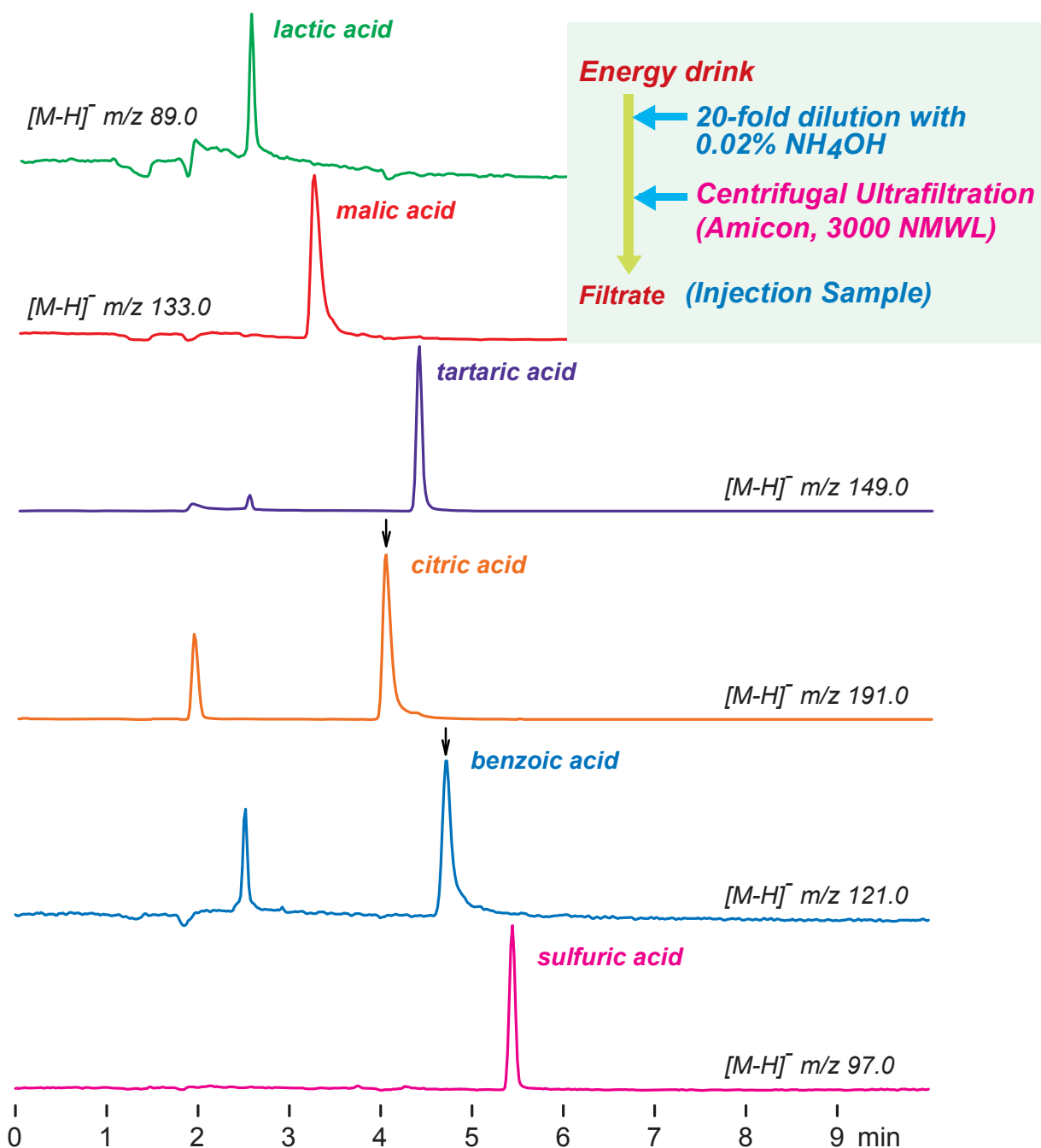

Intrada Organic Acid, 150 x 2 mm

A: acetonitrile / water / HCOOH = 10 / 90 / 0.1

B: acetonitrile / 100mM HCOONH4 = 10 / 90

0%B (0-1min), 0-100%B (1-7min), 100%B (7-10min)

0.2mL/min (12MPa), 37°C, 5µL

MS (ESI, SIM, Negative)