



Scherzo SM-C18
Unison UK-C18

150 x 3 mm

Application

Mononucleotides

モノヌクレオチド

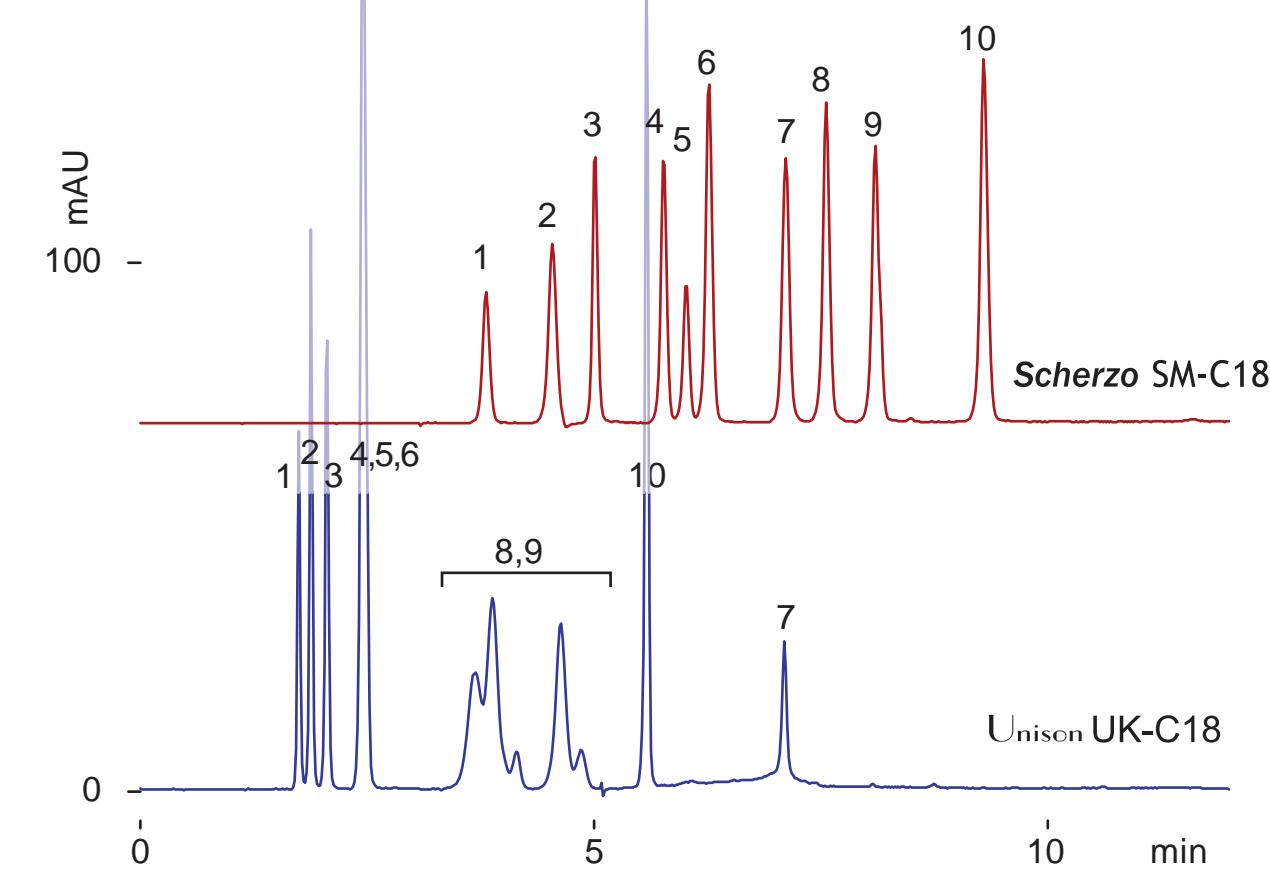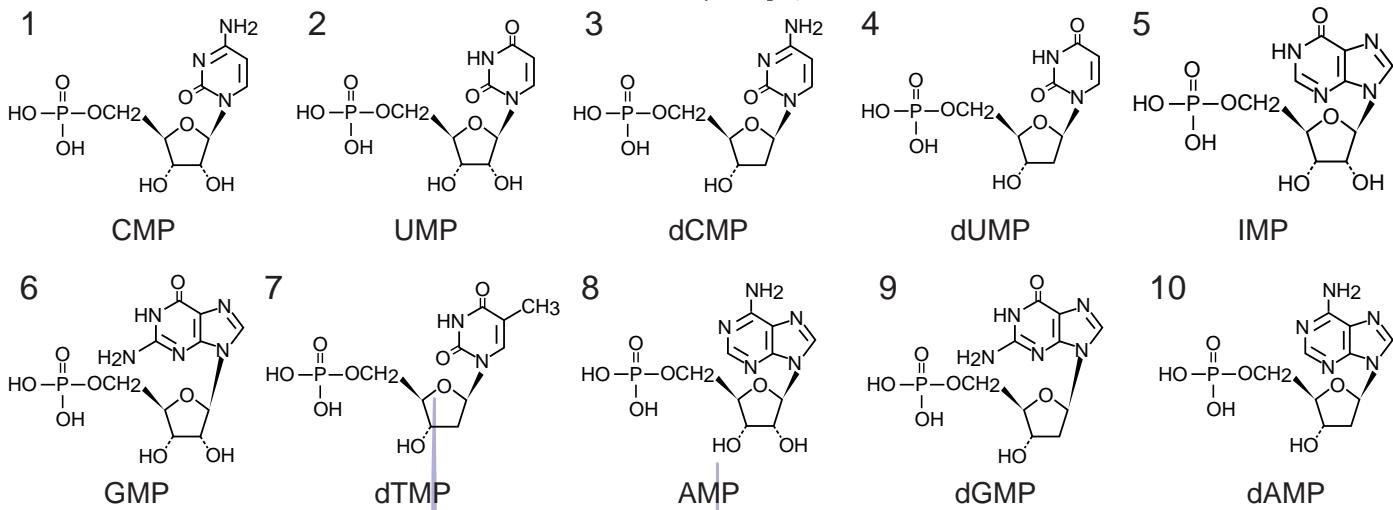

150 x 3 mm

A: 5mM NH₄H₂PO₄ : 5mM (NH₄)₂HPO₄ = 1 : 1, pH 6.8

B: (10mM NH₄H₂PO₄ : 10mM (NH₄)₂HPO₄ = 1 : 1, pH 6.8) /acetonitrile = 50 /50

0-30 %B (0-10min), 0.4 mL/min (10MPa), 37 deg.C, 260 nm, 1 uL (0.2ug)