


Presto FF-C18

150 x 4.6 mm


Application

Polyvinylpyrrolidone (PVP)

ポリビニルピロリドン

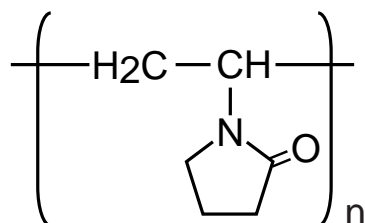

polyvinylpyrrolidone
(ca.50 kDa)

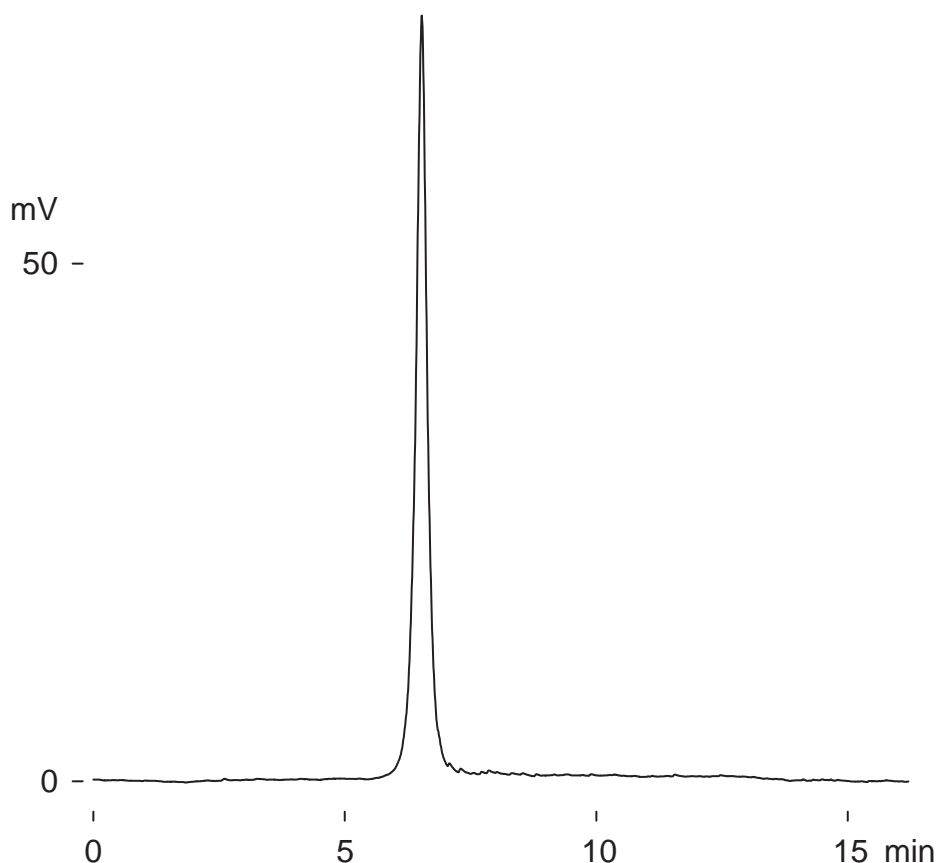

Presto FF-C18, 150 x 4.6mm

A: water / trifluoroacetic acid = 100 / 0.1

B: acetonitrile / trifluoroacetic acid = 100 / 0.1

20-60 %B (0-10min), 0.4 mL/min (23MPa), 37 deg.C

ELSD (spray chamber 50deg.C, drift tube 100deg.C), 1 uL (5ug)