








Unison UK-Amino

250 x 3 mm

Application

Unsaturated chondroitin disaccharides

不飽和コンドロイチン二糖類

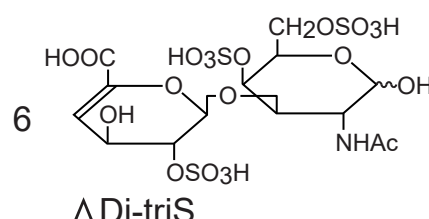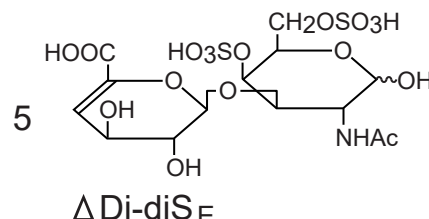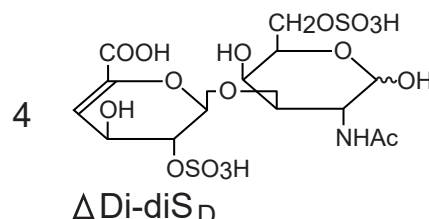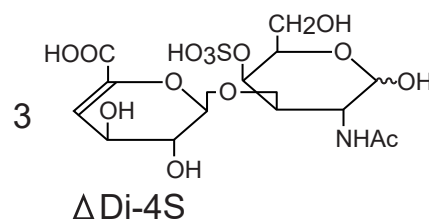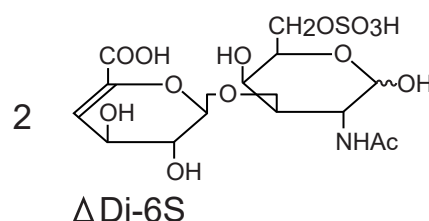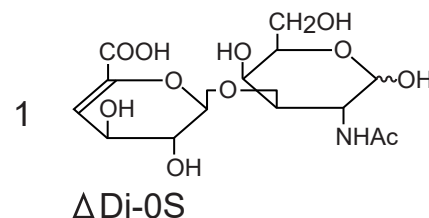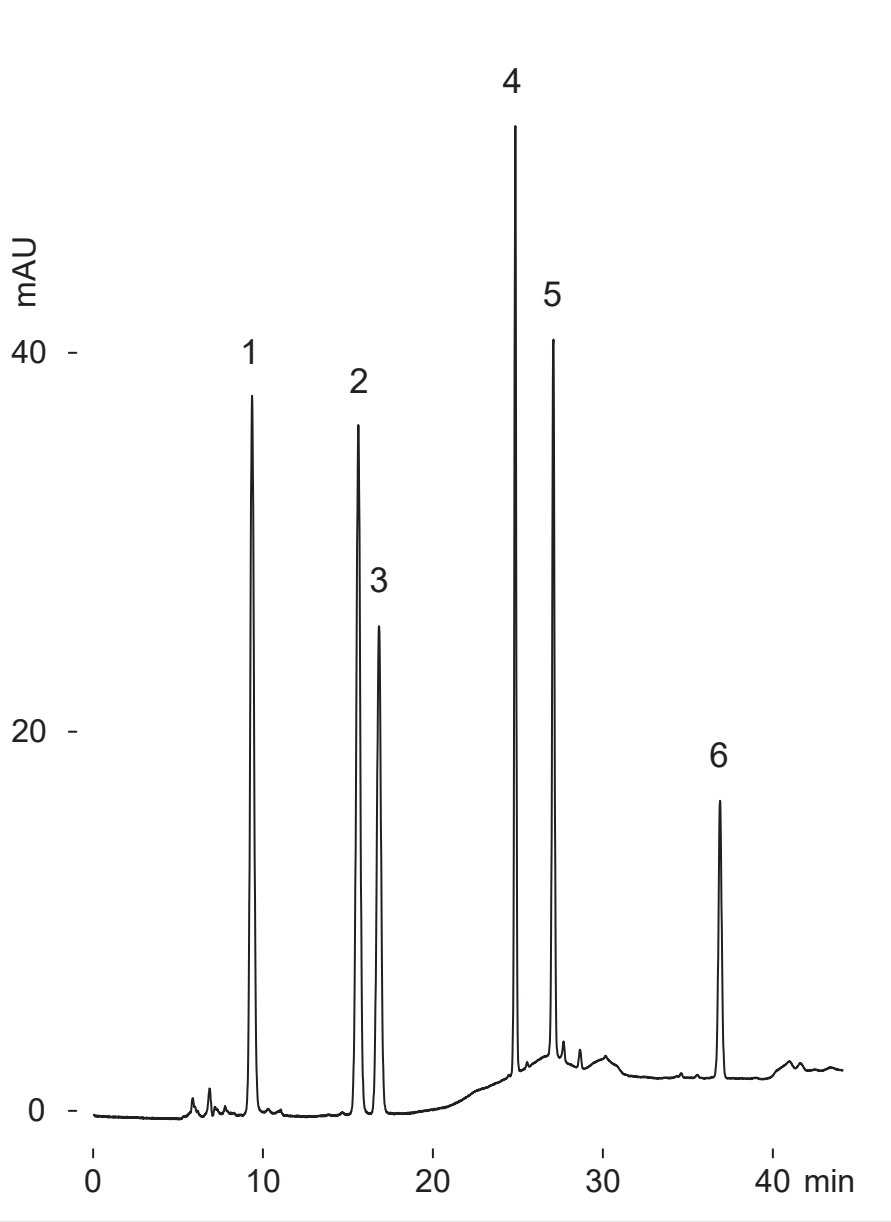

Unison UK-Amino, 250 x 3 mm  
 A: 5 mM (NH<sub>4</sub>)<sub>2</sub>HPO<sub>4</sub> / acetonitrile = 40 / 60  
 B: 100 mM (NH<sub>4</sub>)<sub>2</sub>HPO<sub>4</sub>  
 2-3 %B (0-10 min), 3-40 %B (10-40 min)  
 0.4 mL/min (5 MPa), 50 deg.C, 232 nm  
 4 uL (0.26-0.75 ug, 100 mM (NH<sub>4</sub>)<sub>2</sub>HPO<sub>4</sub>)