



Scherzo SS-C18

150 x 3 mm


Application

Histamine related compounds


ヒスタミン関連化合物

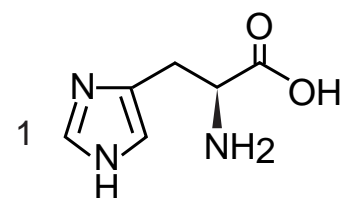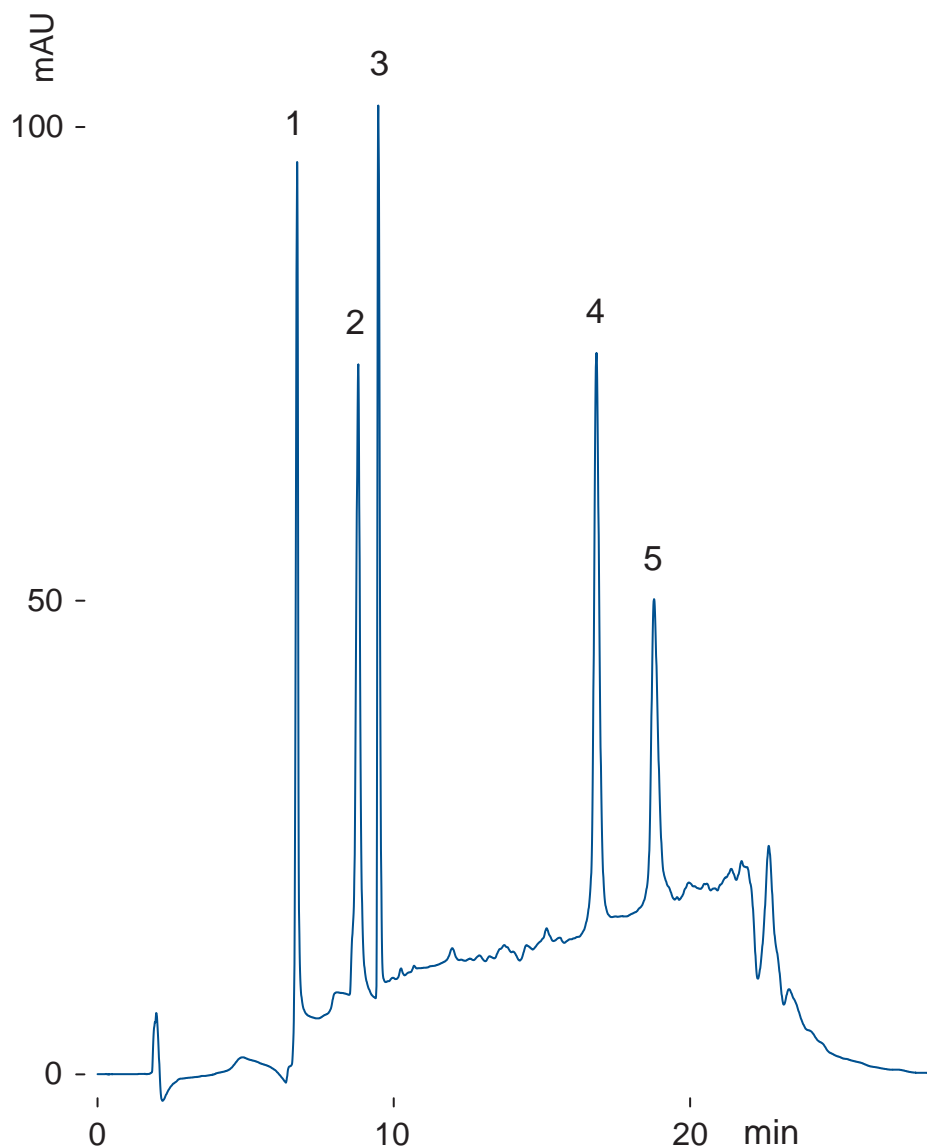


L-histidine

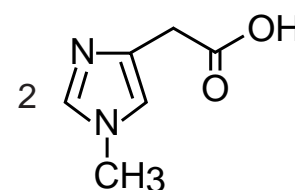


N-methyl-4-imidazoleacetic acid

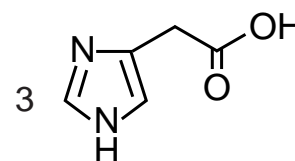

4-imidazoleacetic acid

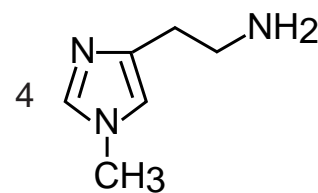

N-methylhistamine

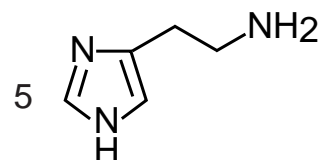

histamine

Scherzo SS-C18, 150 x 3 mm

A: methanol

B: 100 mM ammonium formate : 100 mM ammonium hydroxide = 95 : 5 (pH8.1)

0-20 %B (0-20 min)

0.4 mL/min (6 MPa), 60 deg.C, 230 nm, 1.6 uL (0.6-2.8 ug, MeOH)