


Unison UK-Phenyl

150 x 3 mm

Application

## LC-MS: Isomers separation for fluoroamphetamines and fluoromethamphetamine

LC-MS: フルオロアンフェタミンとフルオロメタンフェタミンの異性体分離

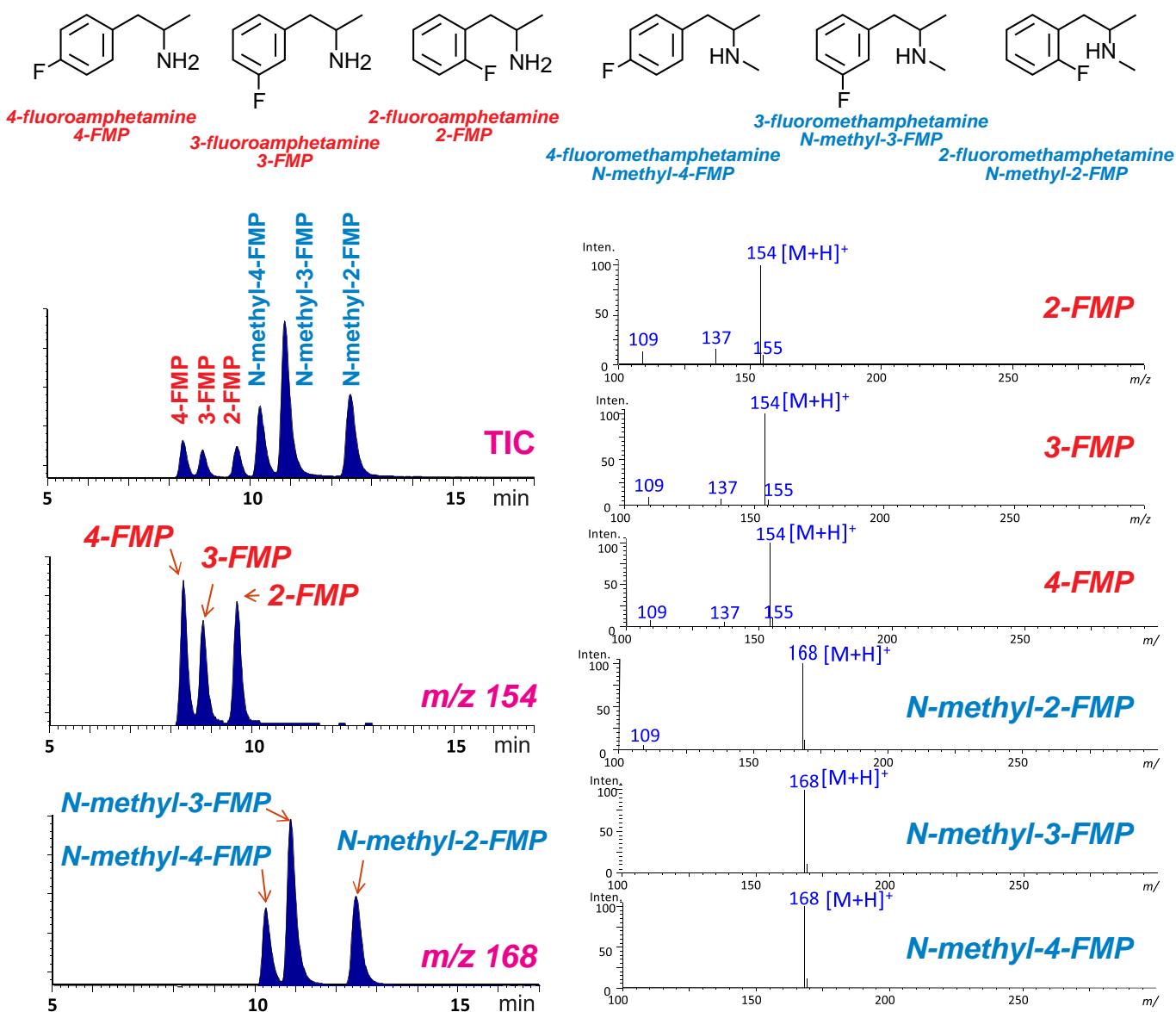

Unison UK-Phenyl, 150 x 3 mm  
 methanol /acetonitrile /10mM ammonium formate (pH 3.5) = 9 /1 /90  
 0.4 mL/min, 40 deg.C, 1-20uL injection  
 MS: ESI, positive

Ref) The 133th Annual Meeting of the Pharmaceutical Society of Japan, Yokohama, 2013