



Scherzo SM-C18

100 x 3 mm


Application

LC-MS: Photo Acid Generators



光酸発生剤のLC-MS分析

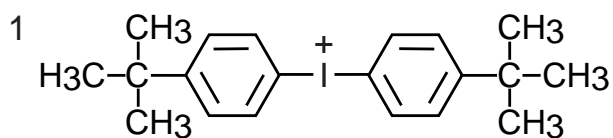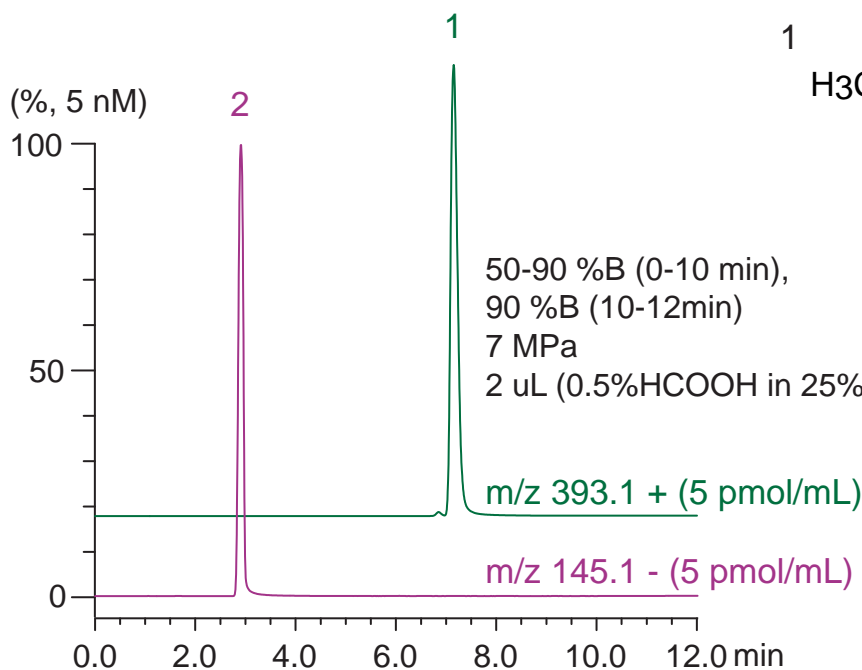


bis(4-tert-butylphenyl)iodonium

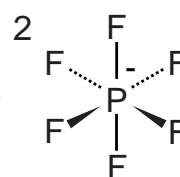

hexafluorophosphate

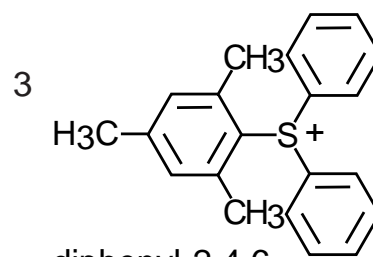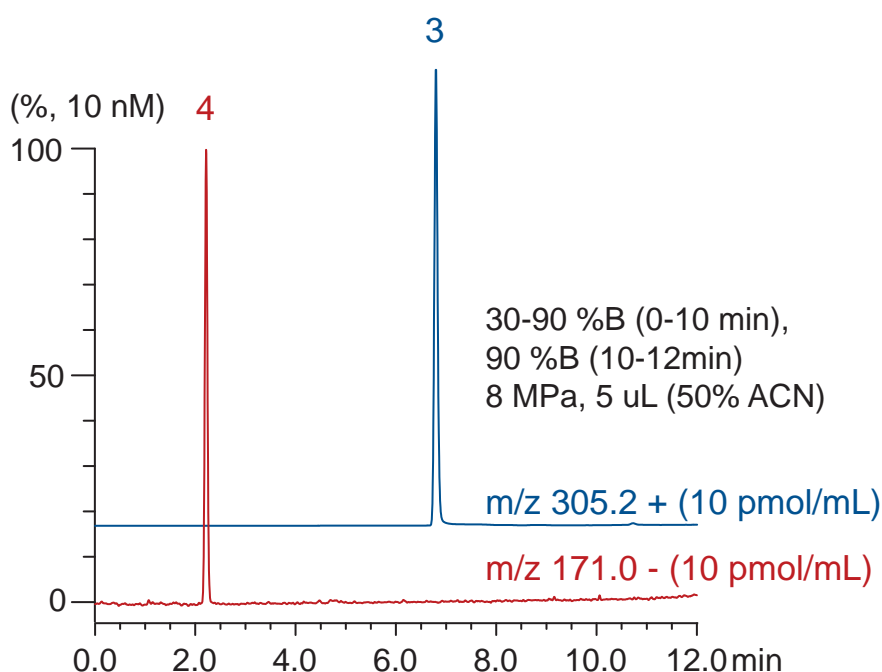

diphenyl-2,4,6-trimethylphenylsulfonium

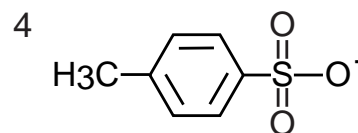

p-toluenesulfonate

Scherzo SM-C18, 100 x 3 mm  
 A: 50mM ammonium formate  
 B: acetonitrile  
 0.4 mL/min, 37°C  
 Single Quad. MS: ESI (SIM, positive / negative)