

Intrada Amino Acid


150 x 3 mm

Application

LC-MS: 20 standard amino acids

20種アミノ酸(標準試料)のLC-MS分析

150 x 3 mm

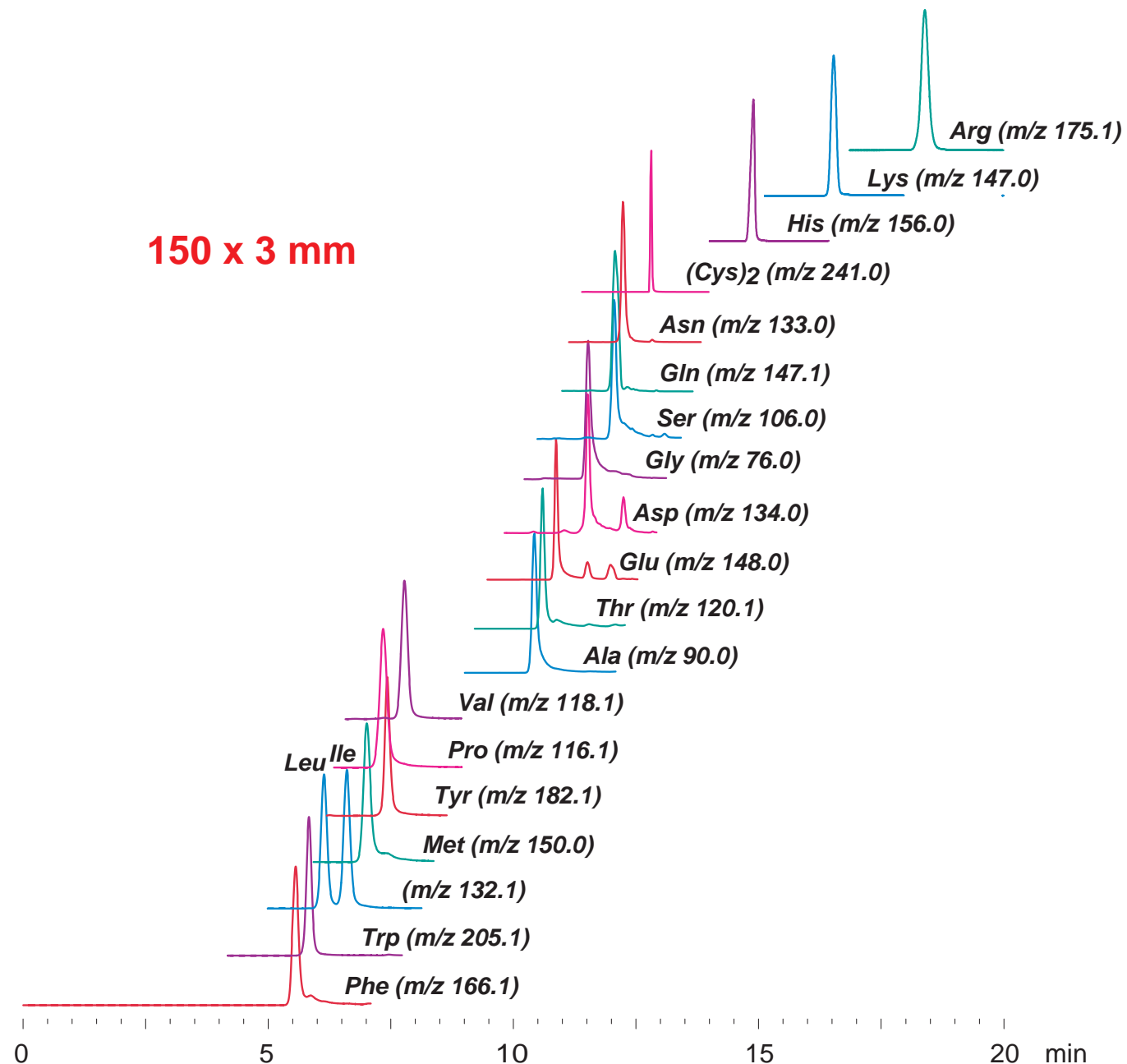

Intrada Amino Acid, 150 x 3 mm

A: acetonitrile / formic acid = 100 / 0.3

B: acetonitrile / 100mM ammonium formate = 20 / 80

20 %B (0-5 min), 20-35 %B (5-11 min), 100 %B (11-20 min)

0.6 mL/min (8 MPa), 37°C, 5 uL (1000nmol/mL, 0.1N HCl)

Single Quad. MS: ESI (SIM, positive)