



Intrada Amino Acid

100 x 3 mm

Application

LC-MS: Nicotine related compounds (3)

ニコチン関連化合物 (3)

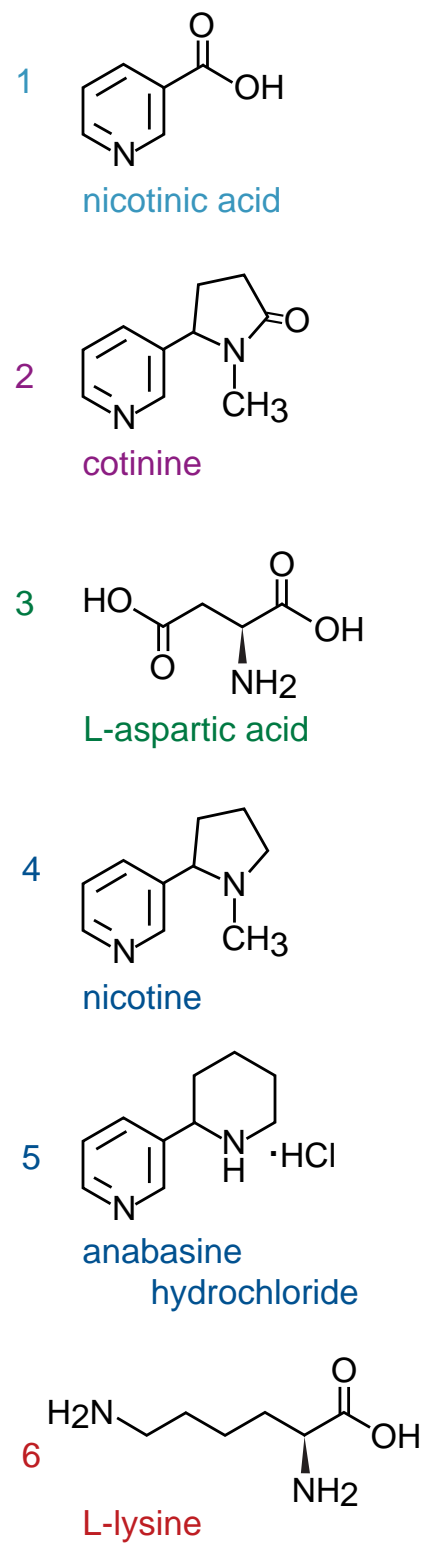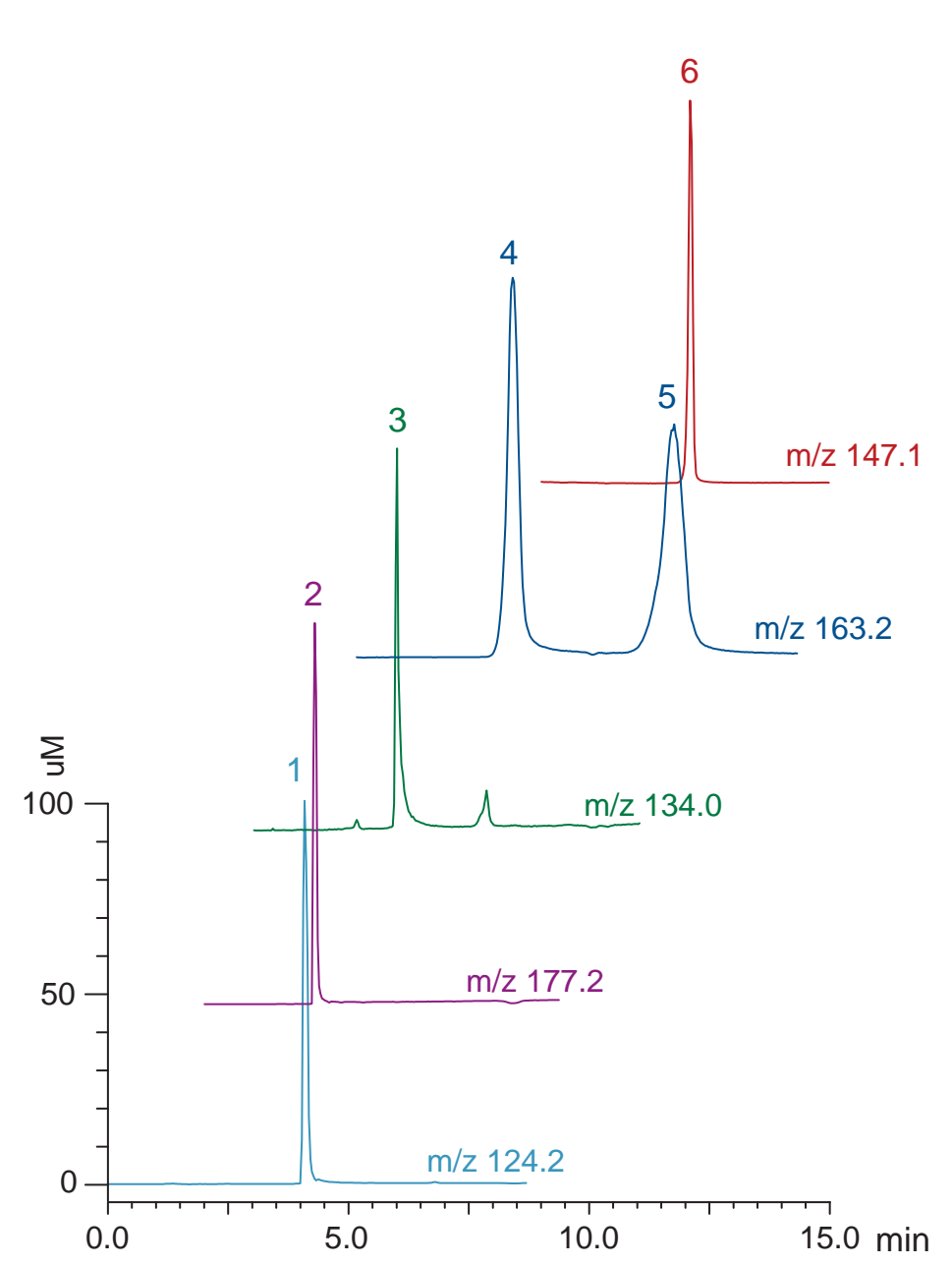

Intrada Amino Acid, 100 x 3 mm

A: acetonitrile / formic acid = 100 / 0.5

B: acetonitrile / 100mM ammonium formate = 20 / 80

0-100 %B (0-10 min), 100 %B (10-15 min)

0.5 mL/min (4 MPa), 37°C, 20 uL (100 nmol /mL, 1% HCOOH)

Single Quad. MS (ESI, positive)