


Scherzo SS-C18

100 x 3 mm

Application

LC-MS: Polyoxin D

ポリオキシン DのLC-MS分析

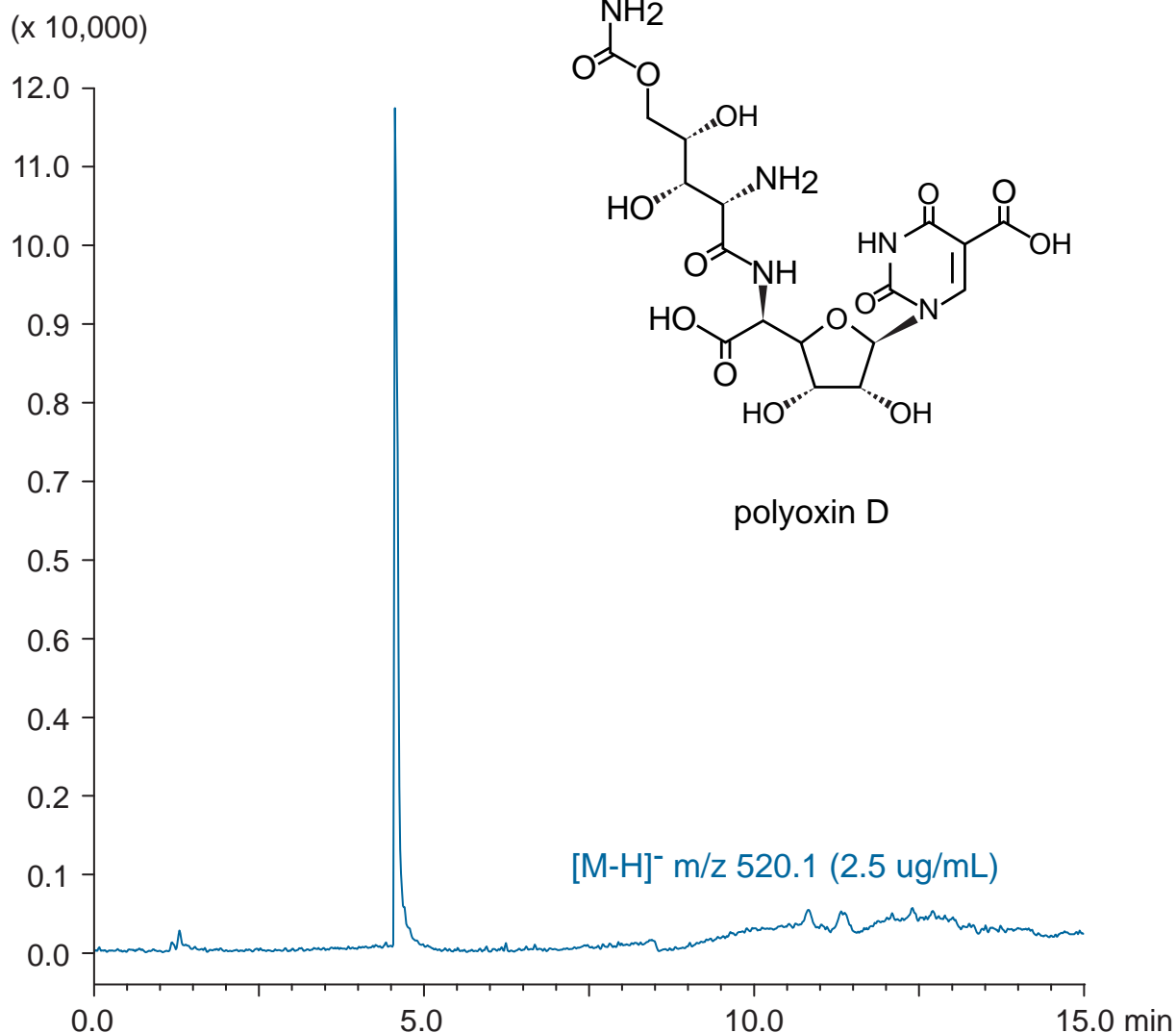

Scherzo SS-C18, 100 x 3 mm

A: water / formic acid = 100 / 0.3

B: acetonitrile / 100mM ammonium formate = 20 / 80

0-100 %B (0-10 min), 100 %B (10-15 min)

0.4 mL/min (6 MPa), 37°C, 10 uL (2.5 ug/mL, 1% HCOOH)

Single Quad. MS (ESI, negative)