


Scherzo SS-C18

100 x 3 mm

Application

LC-MS: 20 standard amino acids

20種アミノ酸(標準試料)のLC-MS分析

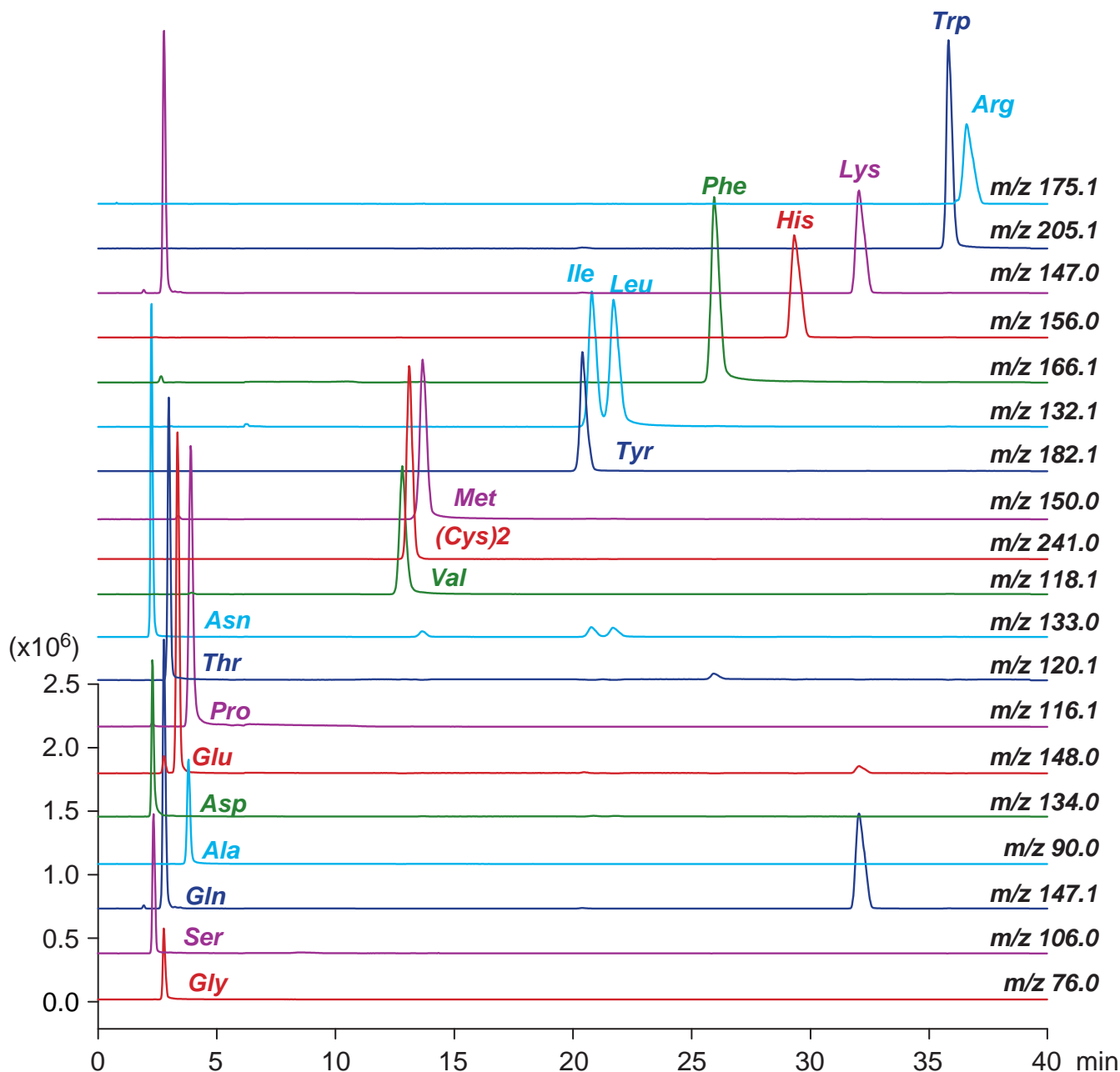

Scherzo SS-C18, 100 x 3 mm
 A: water / formic acid = 100 / 0.1
 B: acetonitrile / 100mM ammonium formate = 40 / 60
 0 %B (0-3 min), 0-10 %B (3-25 min), 10-25 %B (25-40 min)
 0.6 mL/min (9MPa), 37°C, 5 µL (1000nmol/mL, 0.1N HCl)
 Single Quad. MS: ESI (SIM, positive)