

**Intrada SEC**

250 x 4.6 mm

Application


**SEC: Polyvinyl Acetate (PVAc)**

ポリ酢酸ビニル(PVAc)のSEC分離




**Polyvinyl Acetate (PVAc)**

**Intrada SEC**

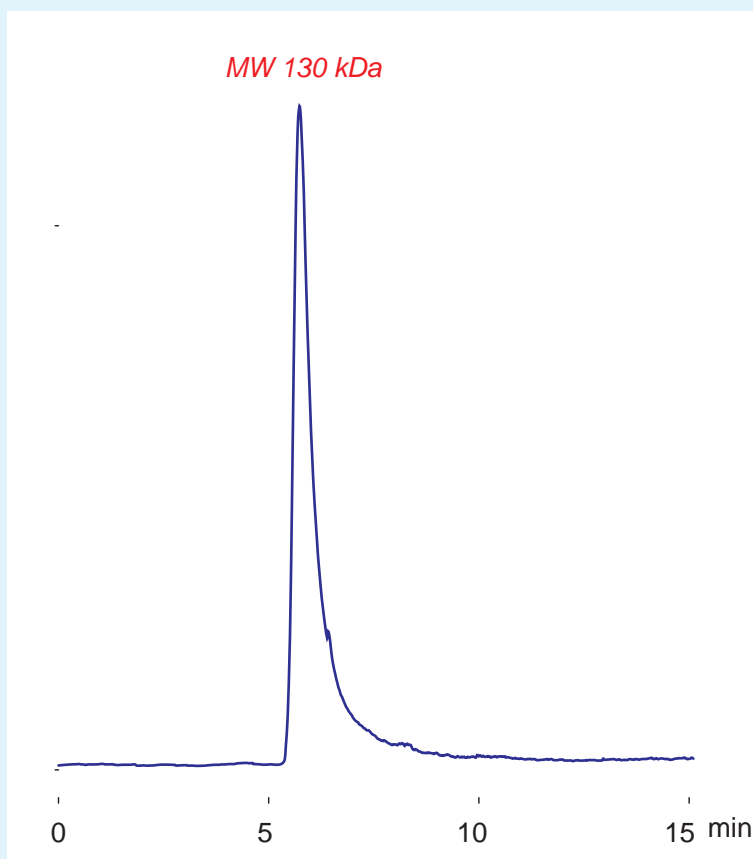

**Vinyl Acetate, Polymer**

n = 1400 - 1600

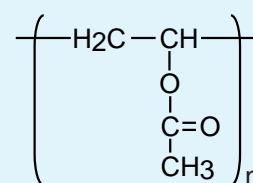

874

*Intrada SEC, 250 x 4.6 mm*

*THF*

*0.3mL/min (3MPa)*

*37 °C*

*ELSD (spray chamber 50°C, drift tube 100°C)*

*1 µL (2mg/mL, THF)*