

Intrada SEC

250 x 4.6 mm

Application

SEC-MS: Elastase, from Porcine pancreas

エラスターゼ(ブタ膵臓由来)のSEC-MS分析




SEC-MS: Elastase


Intrada SEC

Elastase, from Porcine pancreas

20170807_PLUS_IntradaSEC_100mMAF_2017... 8/7/2017 9:25:00 AM

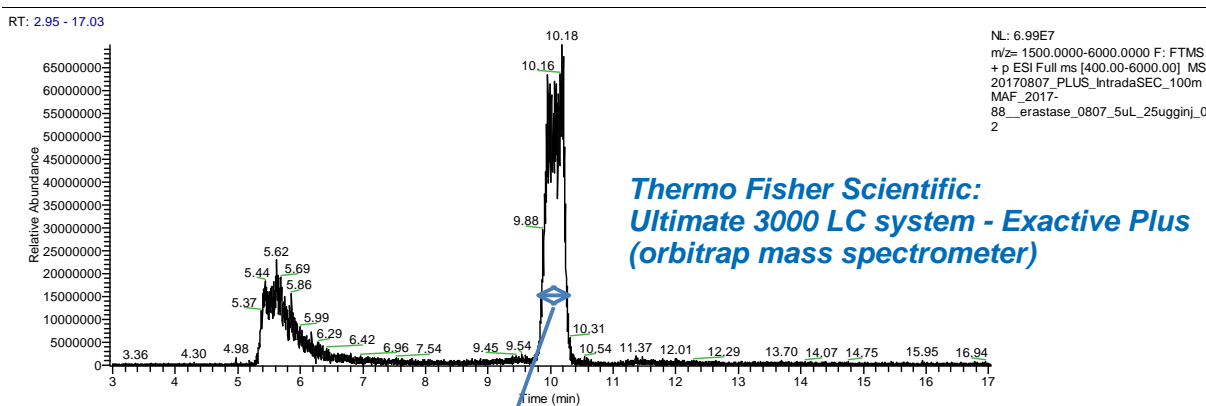

20170807_PLUS_IntradaSEC_100mMAF_2017-88_erastase_0807_5uL_25ugginj_0 #4699-4972 RT: 9.84-10.25 AV: 274 SB: 120 7.24-7.32, 7.98-8.07 NL: 2.59E6 T: FTMS + p ESI Full ms [400.00-6000.00]

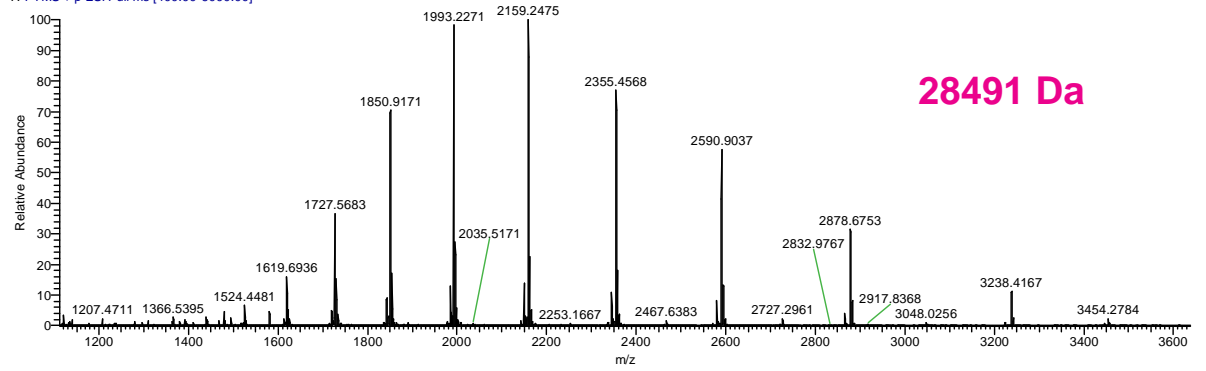

*Intrada SEC, 250x4.6mm
100mM ammonium formate
0.3 mL/min
37°C, 5µL (25µg)*

Courtesy of Dr. Keiko KUWATA, Nagoya University, JAPAN