


Scherzo SM-C18

50 x 3 mm

Application

LC-MS: Ionic Liquid (1-butyl-1-methylpyrrolidinium bis(trifluoromethylsulfonyl)imide)

イオン液体 (1-ブチル-1-メチルピロリジニウム
ビス(トリフルオロメチルスルホニル)イミド

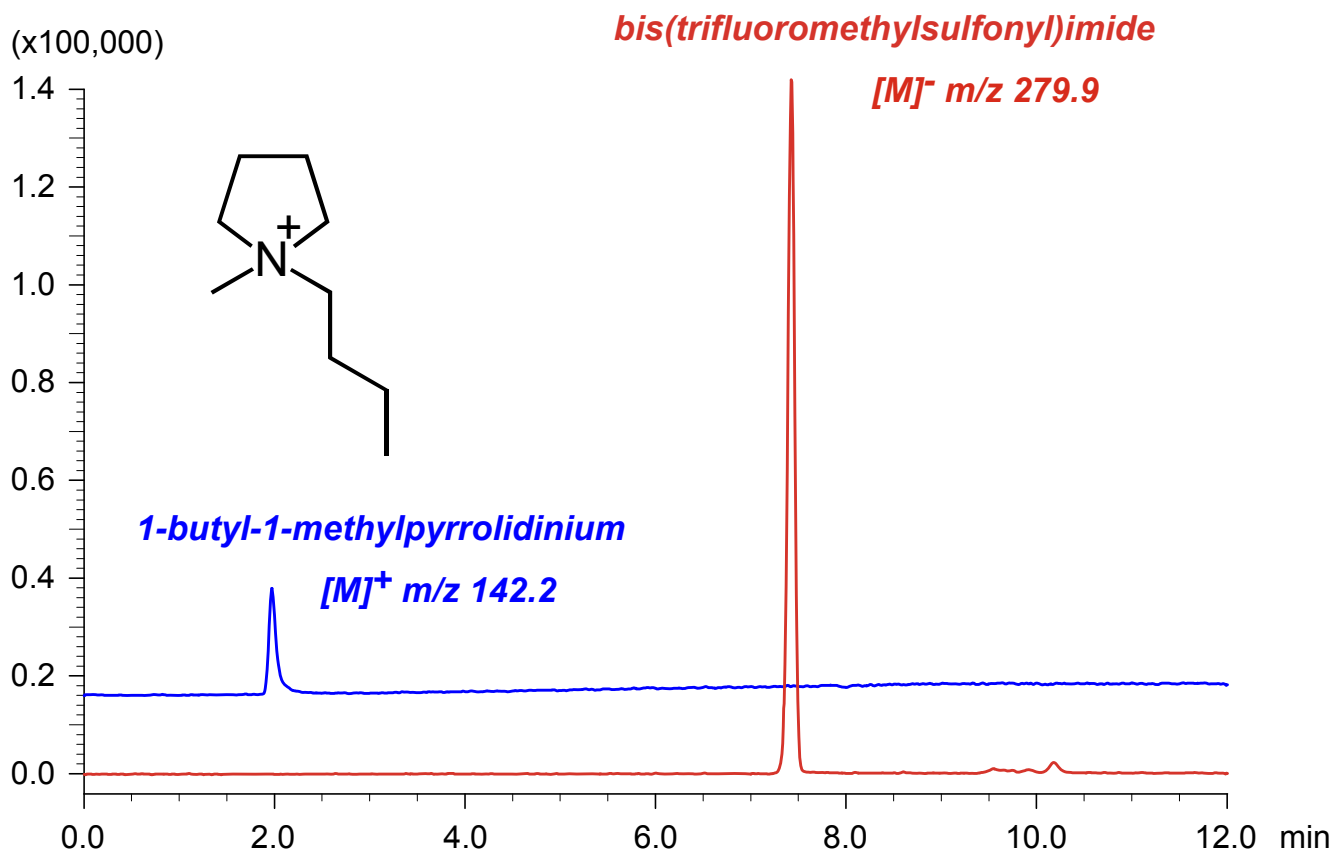

Scherzo SM-C18, 50 x 3 mm
 A: 20mM ammonium formate
 B: 100mM ammonium formate / acetonitrile = 30 / 70
 15-100%B (0-8 min), 100%B (8-12 min)
 0.4 mL/min (4MPa), 37°C
 2 μL (100ng/mL in 50mM HCOONH₄ / ACN = 50 / 50)
 Single Quad. MS (ESI)