


Unison UK-Phenyl

100 x 3 mm


Application

LC-MS: Carbonic anhydrase inhibitor (antigraaucoma drug)


炭酸脱水酵素阻害剤 (緑内障治療薬)

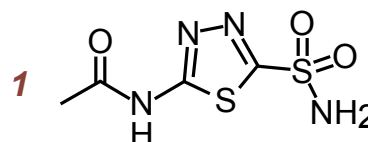

acetazolamide

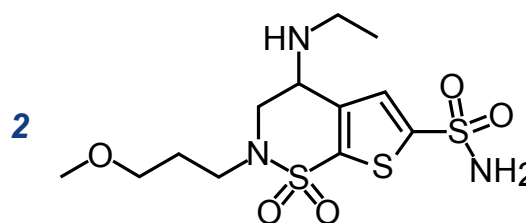

brinzolamide

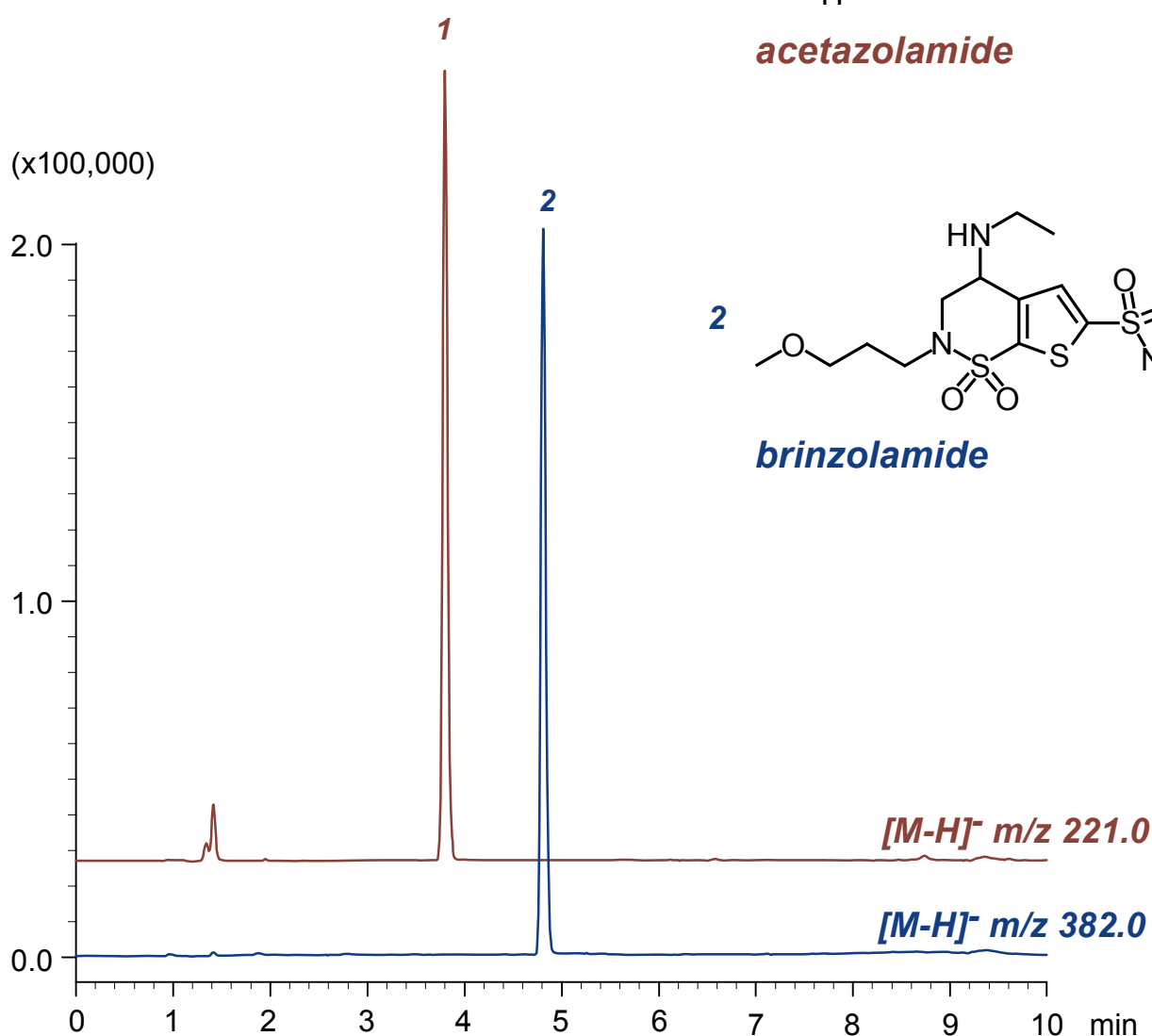

Unison UK-Phenyl , 100 x 3 mm
 A: water/ formic acid = 100 / 0.1
 B: acetonitrile / formic acid = 100 / 0.1
 0-50%B (0-8 min), 90%B (8-10 min)
 0.4mL/min (6 MPa), 37°C
 10 µL (500ng/mL, 0.1% formic acid)
 Single Quad. MS (ESI)