


**Intrada Organic Acid**

150 x 2 mm

Application

LC-MS: Polyvalent / polycarboxylic acids (Normal-Phase Mode)

多価酸, ポリカルボン酸のLC-MS分析 (順相モード)

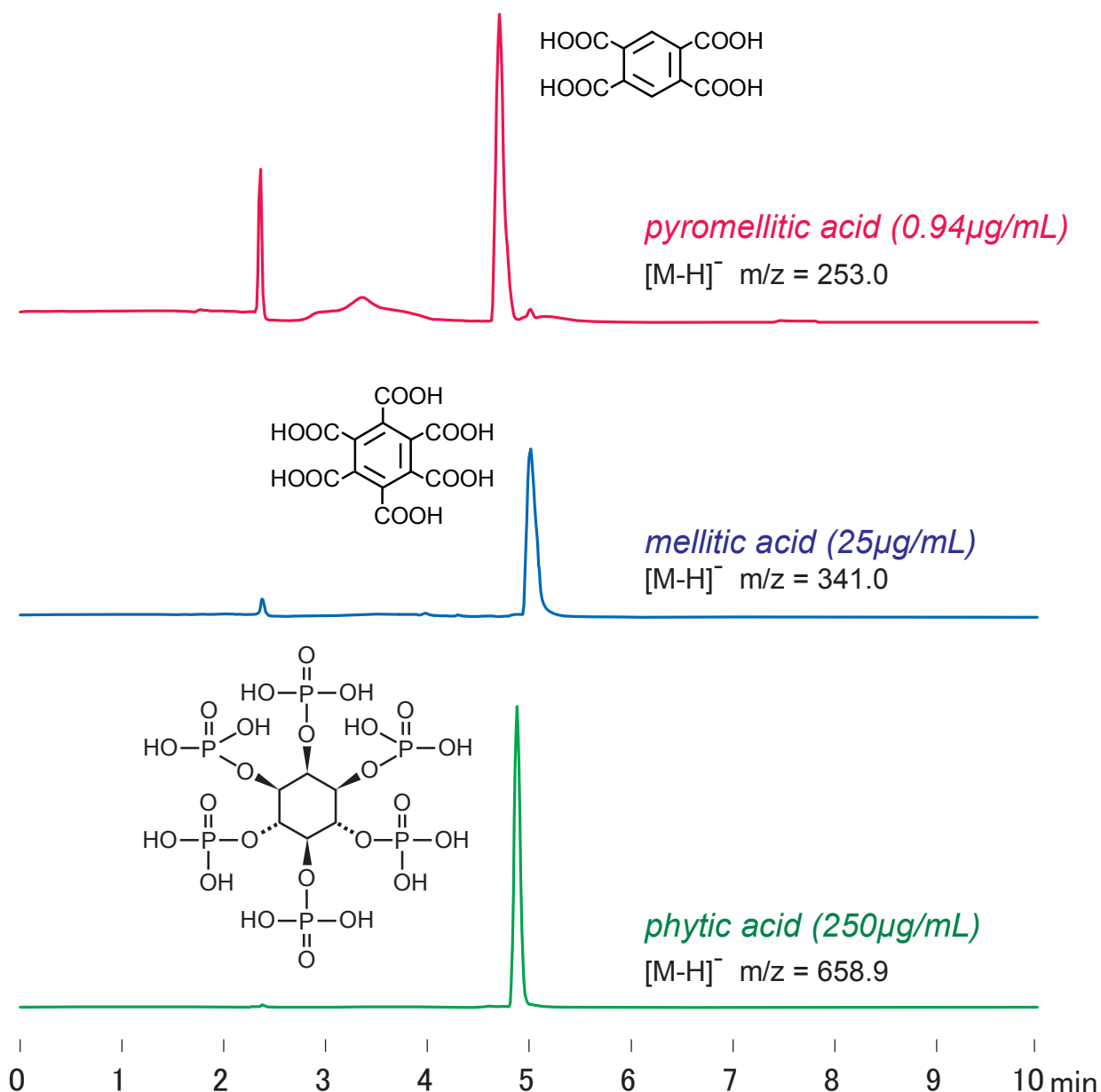

**Intrada Organic Acid, 150 x 2 mm**

**A: acetonitrile / 100mM ammonium acetate = 95 / 5**

**B: acetonitrile / 100mM ammonium acetate = 10 / 90**

**0%B (0-1 min), 0-100%B (1-4 min), 100%B (4-10 min)**

**0.2 mL/min (5-6MPa), 37 °C, MS(ESI, SIM)**

**5 μL (50mM ammonium acetate solution)**