


Intrada Organic Acid

150 x 2 mm

Application

LC-MS: Sulfonic acids and sulfuric acids

スルホン酸, 硫酸類のLC-MS分析

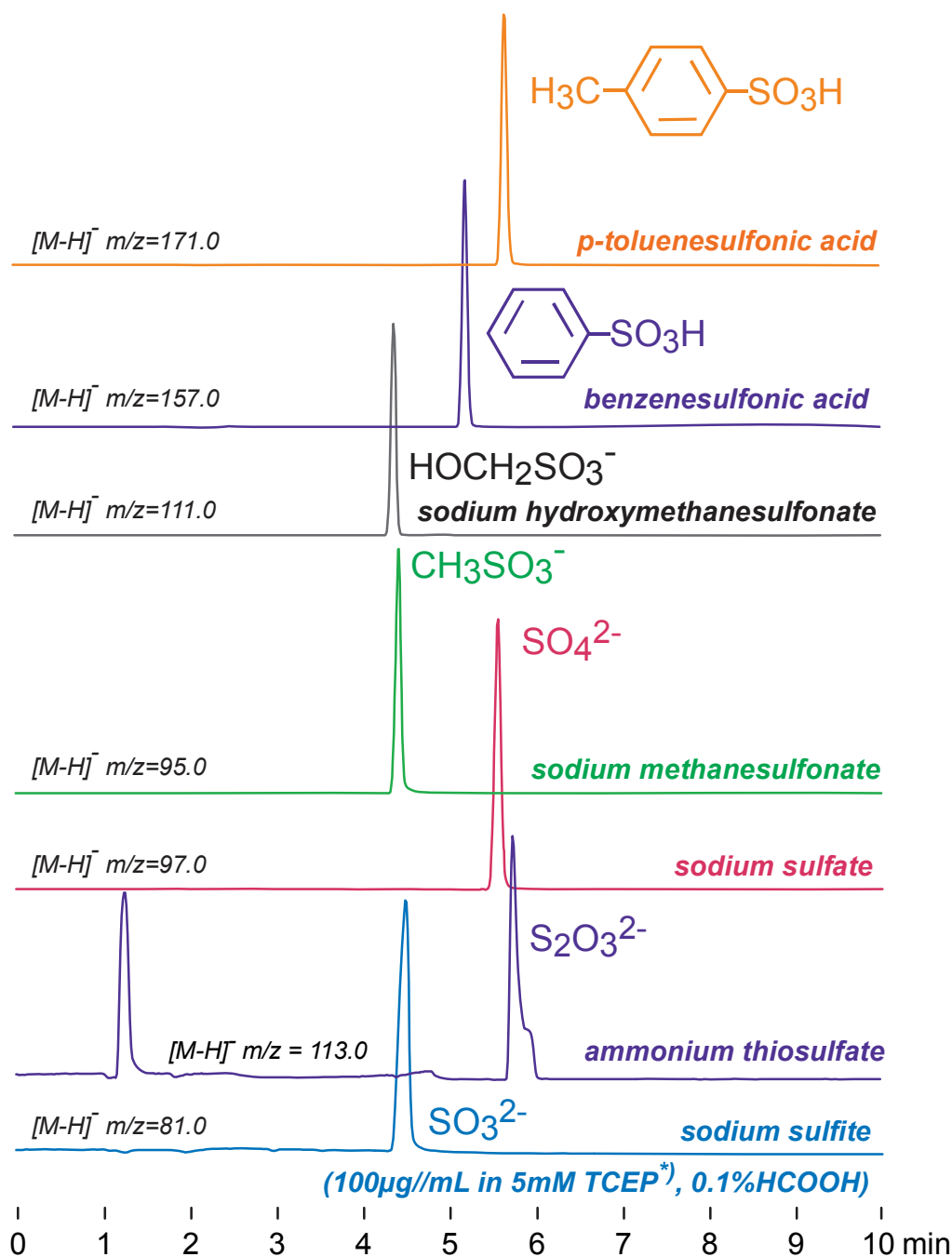

^{*}) TCEP: tris(2-carboxyethyl)phosphine)

Intrada Organic Acid, 150 x 2mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / 100mM ammonium formate = 10 / 90

0%B (0-1min), 0-100%B (1-7min), 100%B (7-10min)

0.2mL/min, 37°C, MS (ESI, SIM, Negative)

5 μ L (50-100 μ g/mL in 0.1% HCOOH)