

Cadenza CD-C18

Cadenza CX-C18


Unison UK-C18

50 x 3 mm

Application


LC-MS: Omarigliptin (Oral antidiabetic drug)

オマリグリプチン(経口糖尿病用剤)のLC-MS分析

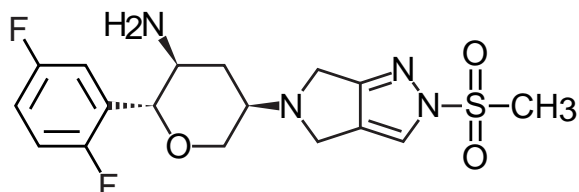

omarigliptin

$[M+H]^+$ m/z 399.1

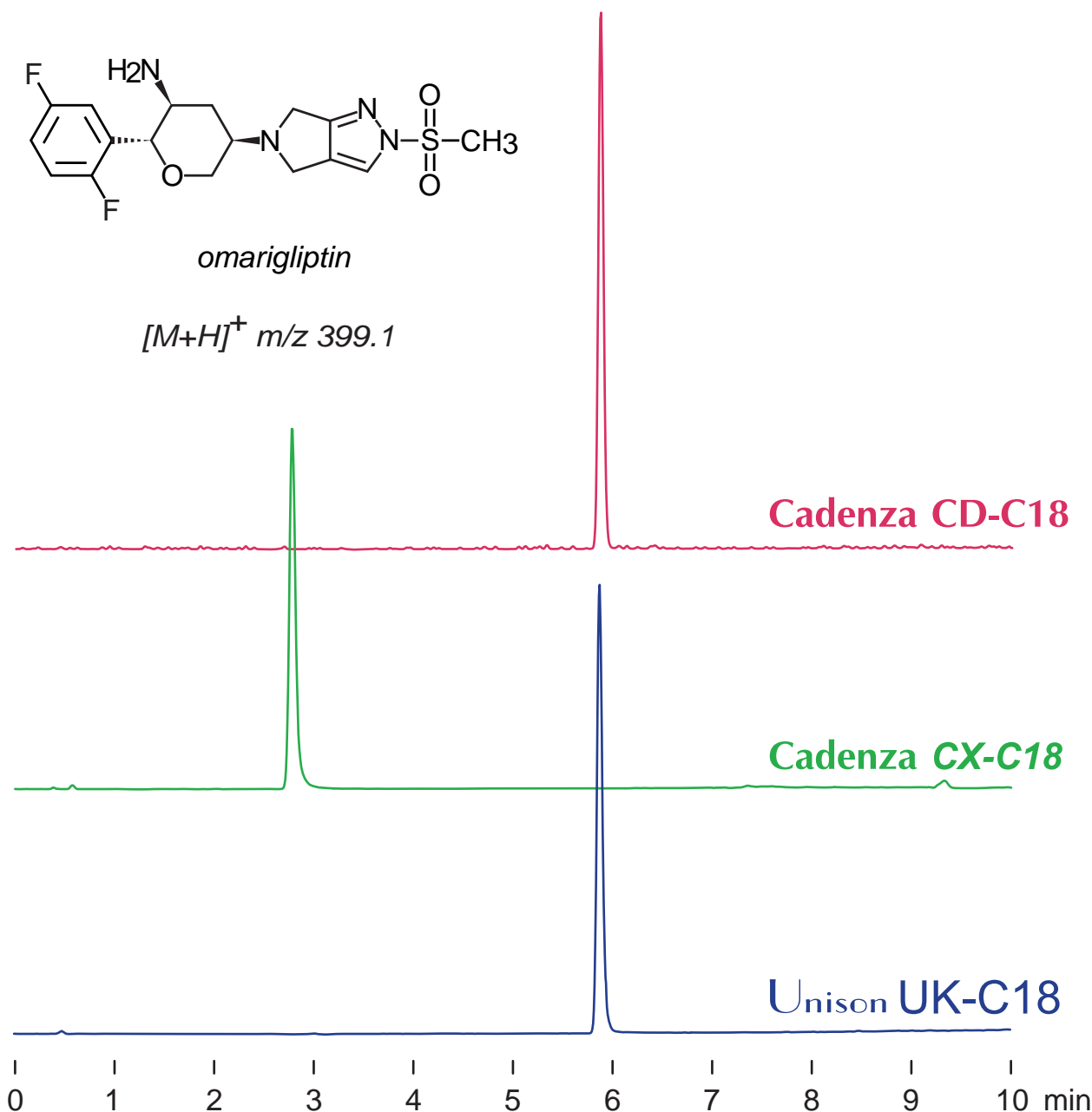

50 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-30%B (0-10min)

0.4mL/min (4MPa), 37°C, MS (ESI, SIM, positive)

5µL (1µg/mL, 0.1% HCOOH)