




Unison UK-Phenyl
Scherzo SS-C18

100 x 3 mm

Application

LC-MS: N-nitrosoamines

N-ニトロソアミンのLC-MS分析

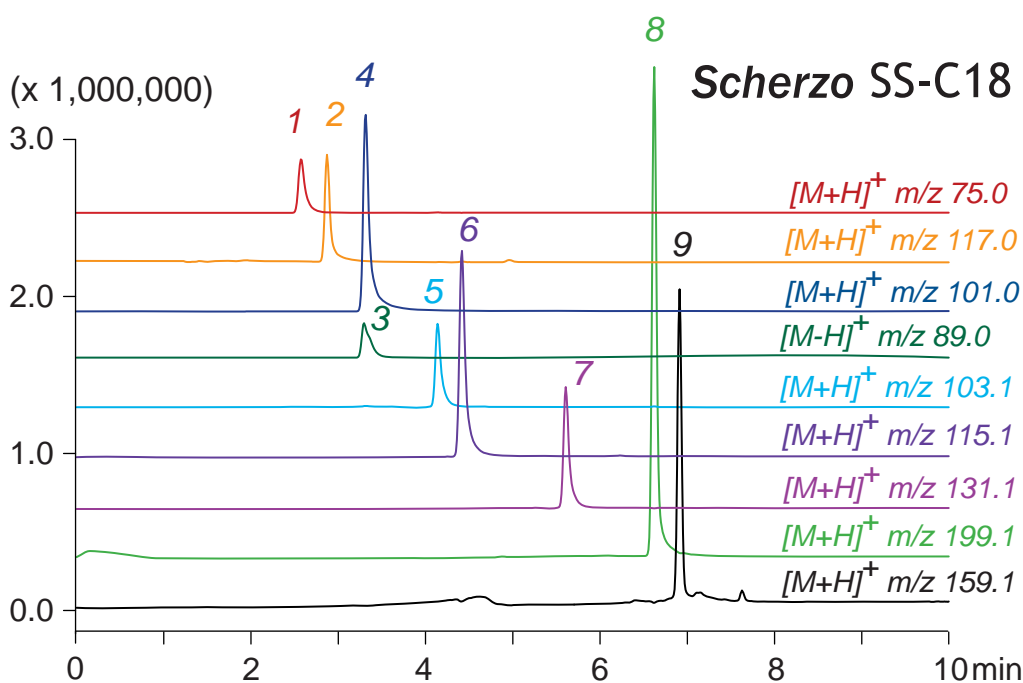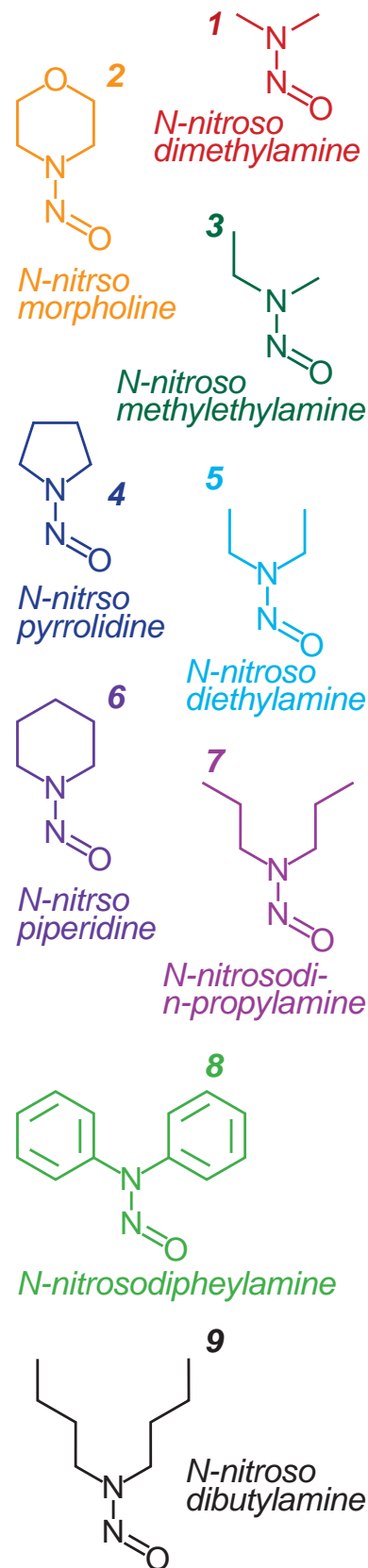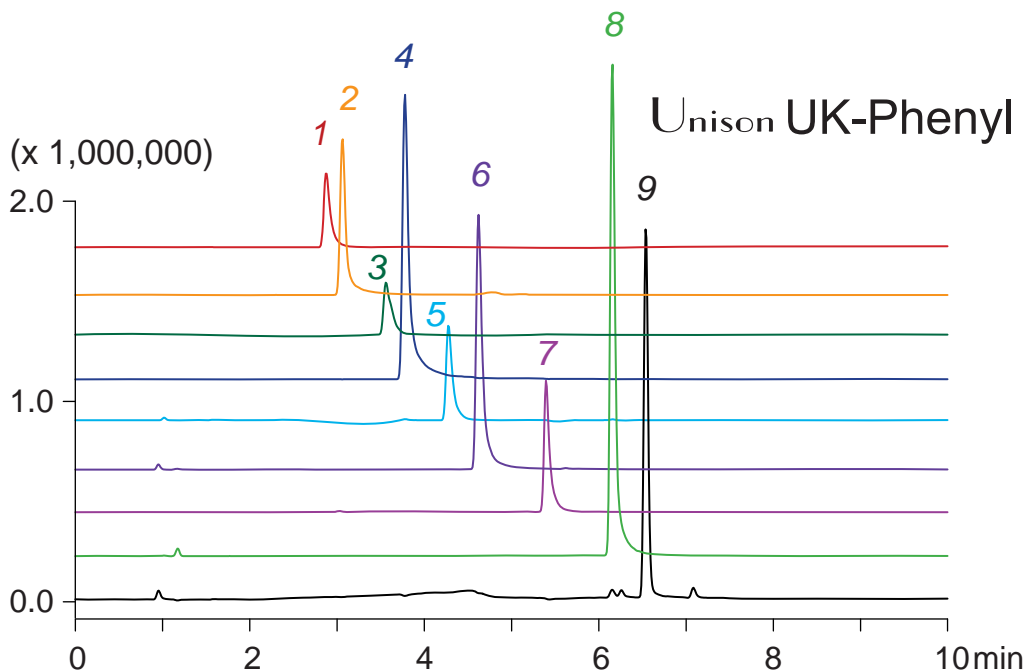

100 x 3 mm

A: water / formic acid = 100 / 0.1

B: methanol / formic acid = 100 / 0.1

0-90%B (0-6min), 90%B (6-10min)

0.4mL/min (6-7MPa), 37°C, MS (ESI, SIM, positive)

5µL (10µg/mL, MeOH/water/HCOOH = 50/50/0.1)