


Intrada Organic Acid

150 x 2 mm

Application

LC-MS: Inorganic anions

無機陰イオンのLC-MS分析

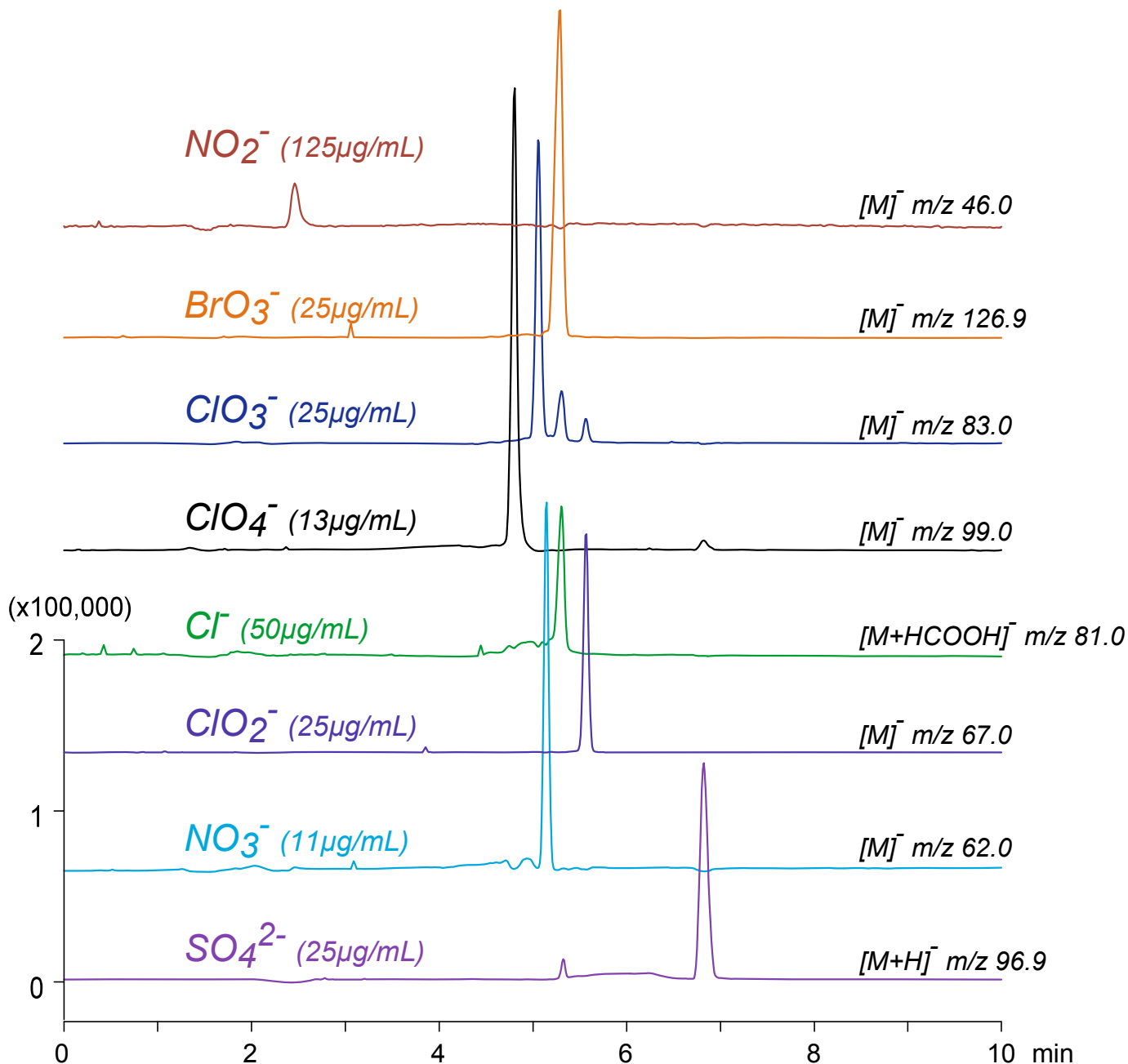

Intrada Organic Acid, 150 x 2 mm

A: acetonitrile / water / formic acid = 70 / 30 / 0.3

B: acetonitrile / 100mM ammonium formate = 10 / 90

0%B (0-1.5min), 0-100%B (1.5-8min), 100%B (8-10min)

0.2mL/min (8MPa), 37°C, MS (ESI, SIM)

5μL (water)