


Intrada Organic Acid

100 x 2 mm

Application

LC-MS: Deoxyoligonucleotides

デオキシオリゴヌクレオチドのLC-MS分析

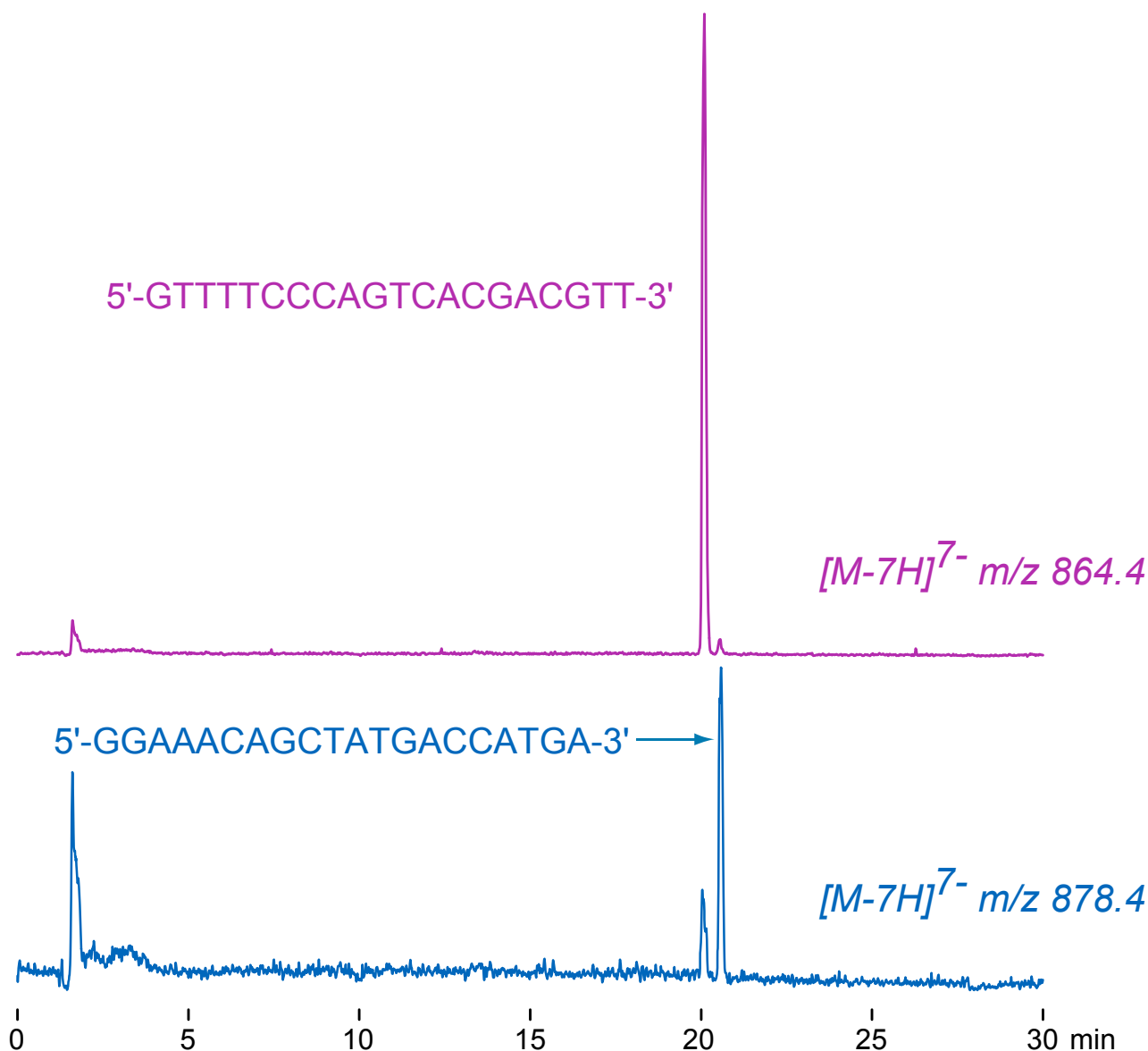

Intrada Organic Acid, 100 x 2 mm

A: acetonitrile / 100mM ammonium formate = 95 / 5

B: acetonitrile / 100mM ammonium formate = 10 / 90

10%B (0-1min), 10-70%B (1-30min)

0.2mL/min (3MPa), 60°C, MS (ESI, SIM)

100µL (0.8nmol/mL, 0.05% HCOOH)