

Scherzo SW-C18
Scherzo SS-C18
Scherzo SM-C18

100 x 3 mm


Application

LC-MS: Strong organic base catalysts


有機強塩基触媒のLC-MS分析




1,4-diazabicyclo[2.2.2]octane
diazabicyclooctane (DBO)

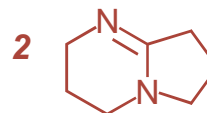

1,5-diazabicyclo[4.3.0]non-5-ene
diazabicyclononene (DBN)

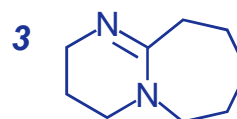

1,8-diazabicyclo[5.4.0]undec-7-ene
diazabicycloundecene (DBU)

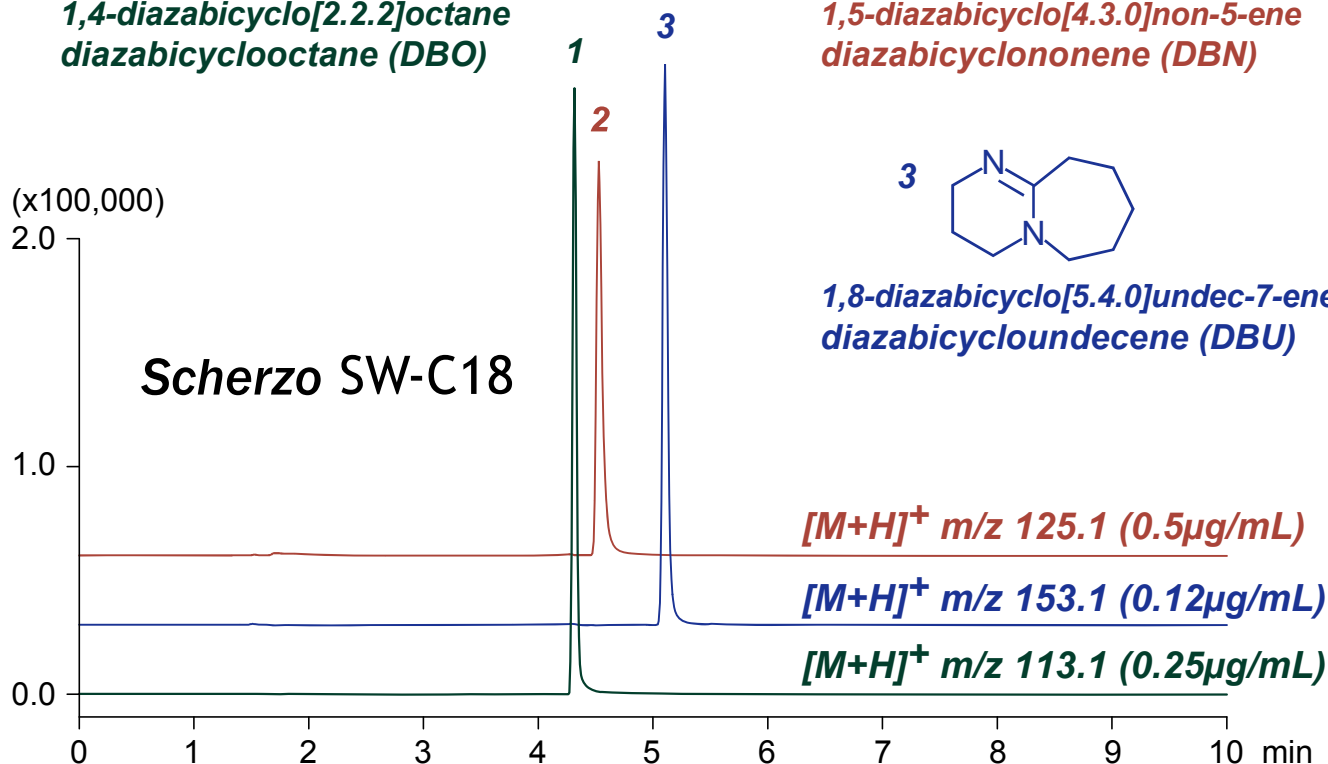

100 x 3 mm
A: water / formic acid = 100 / 0.1
B: acetonitrile / 100mM ammonium formate = 60 / 40
0%B (0-2min), 0-100%B (2-6min), 100%B (6-10min)
0.4mL/min (6MPa), 37°C, MS (ESI, SIM)
5 μ L (0.12-0.50 μ g/mL, 0.5% HCOOH)

