

Unison UK-C1

150 x 4.6 mm


Application

Lipo-Oxytocin-1 (LOT-1)

リポオキシトシン-1

Lipo-Oxytocin-1 (Code: 3422-v / Peptide Institute, Japan)

Palmitoyl-Cys-Tyr(Palmitoyl)-Ile-Gln-Asn-Cys-Pro-Leu-Gly-NH₂

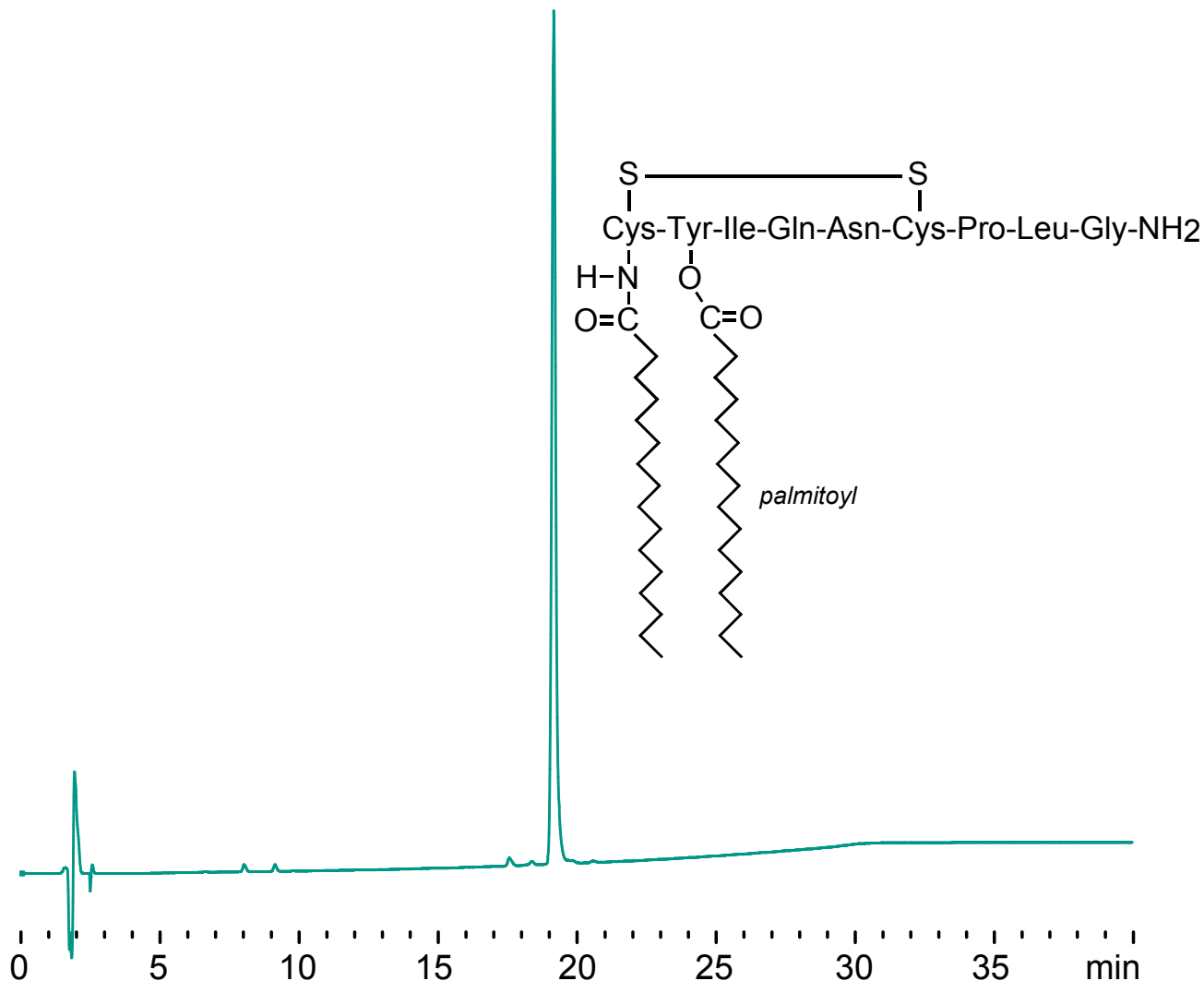

Unison UK-C1, 150 x 4.6 mm
 A: 0.05% TFA in water
 B: 0.04% TFA in acetonitrile
 50-80%B (0-25min), 80%B (25-40min)
 1mL/min (11MPa), 25°C, 220nm
 0.6µL (0.5mg/50µL in trifluoroethanol)

Courtesy of PEPTIDE INSTITUTE, INC. (JAPAN)