

Cadenza **CX-C18**



Unison **UK-C18**

150 x 3 mm

HPLC COLUMN APPLICATION

Remdesivir (RNA polymerase inhibitor)

レムデシビル (抗ウイルス薬)

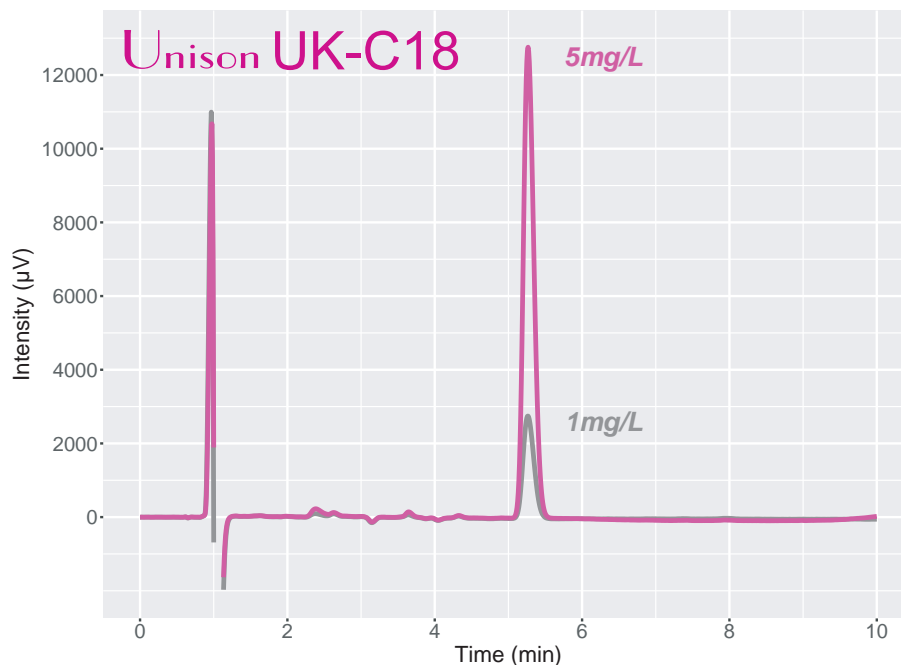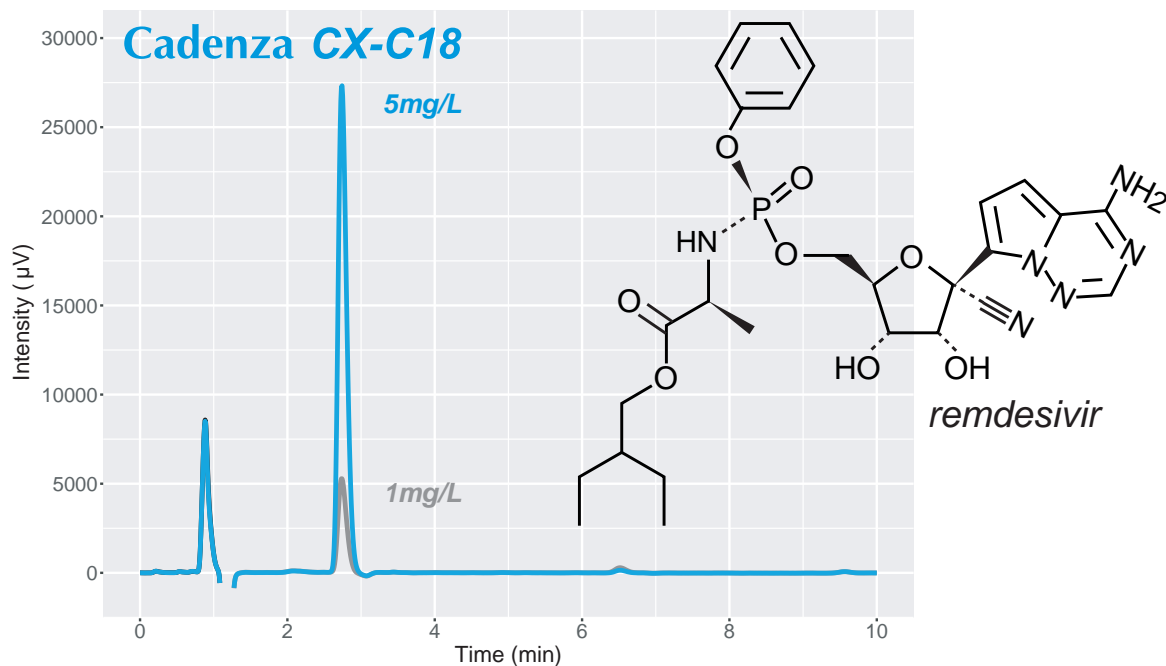

150 x 3 mm

water / acetonitrile / formic acid = 60 / 40 / 0.1 (v/v/v)

0.8mL/min, 40°C, 244nm

20µL (acetonitrile)

Courtesy of Prof. Dr.TSUJI Yasuhiro and OGAMI Chika, NIHON UNIVERSITY, JAPAN