

Unison UK-Phenyl

Unison UK-C8

Cadenza CX-C18

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm

Application

*Pranlukast (Antiasthmatic, Leukotriene receptor antagonist)*

برانلکاست (喘息治療薬, ロイコトリエン受容体拮抗薬)

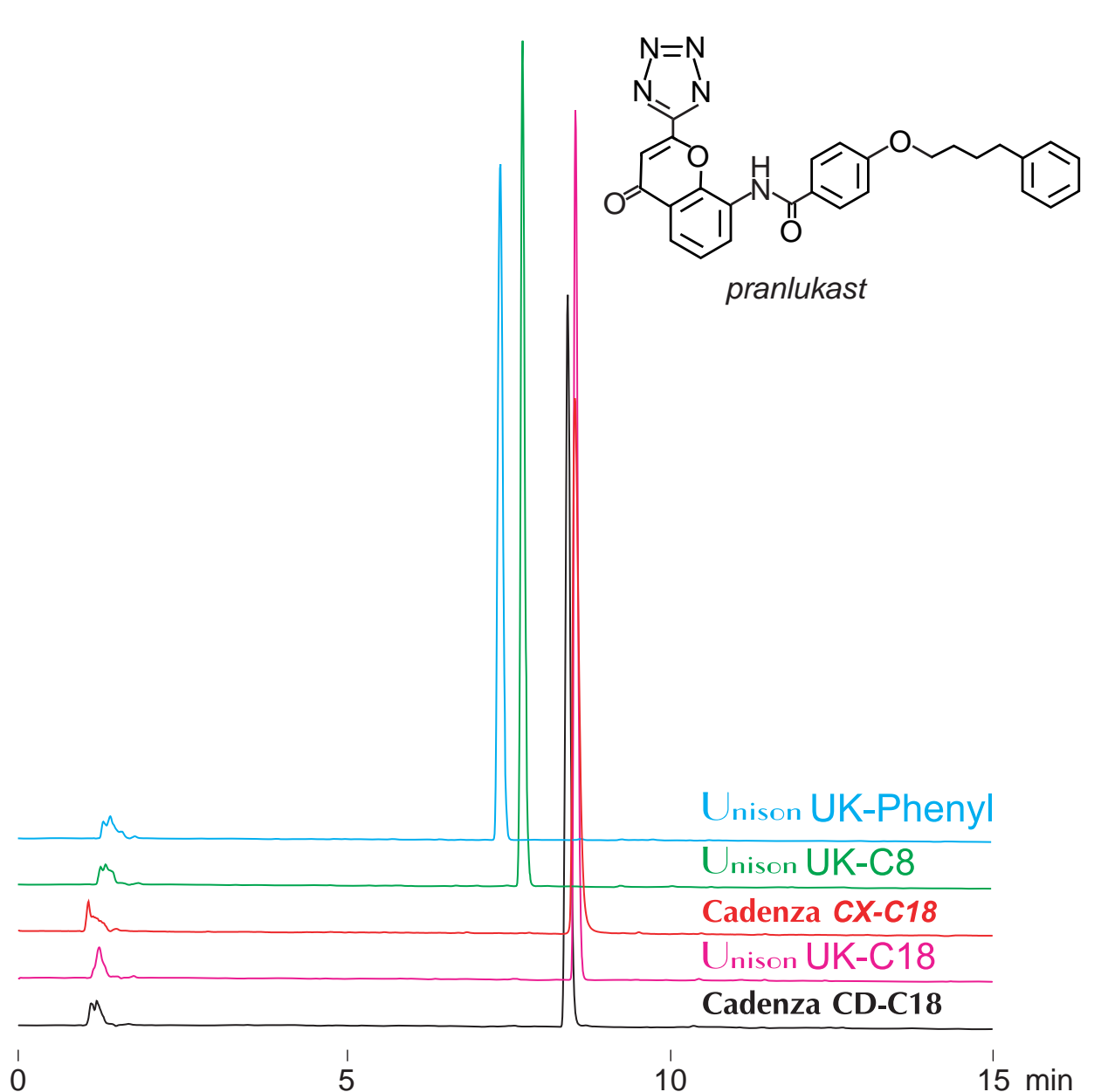

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate= 90 / 10

25-100%B (0-15min)

1.0 mL/min (6MPa), 37°C, 260nm

10 µL (1mg/mL, dimethyl sulfoxide)