

Unison UK-Phenyl

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

Application

*Diazepam (Antianxiety, benzodiazepine derivatives)*

ジアゼパム (ベンゾジアゼピン系抗不安薬)

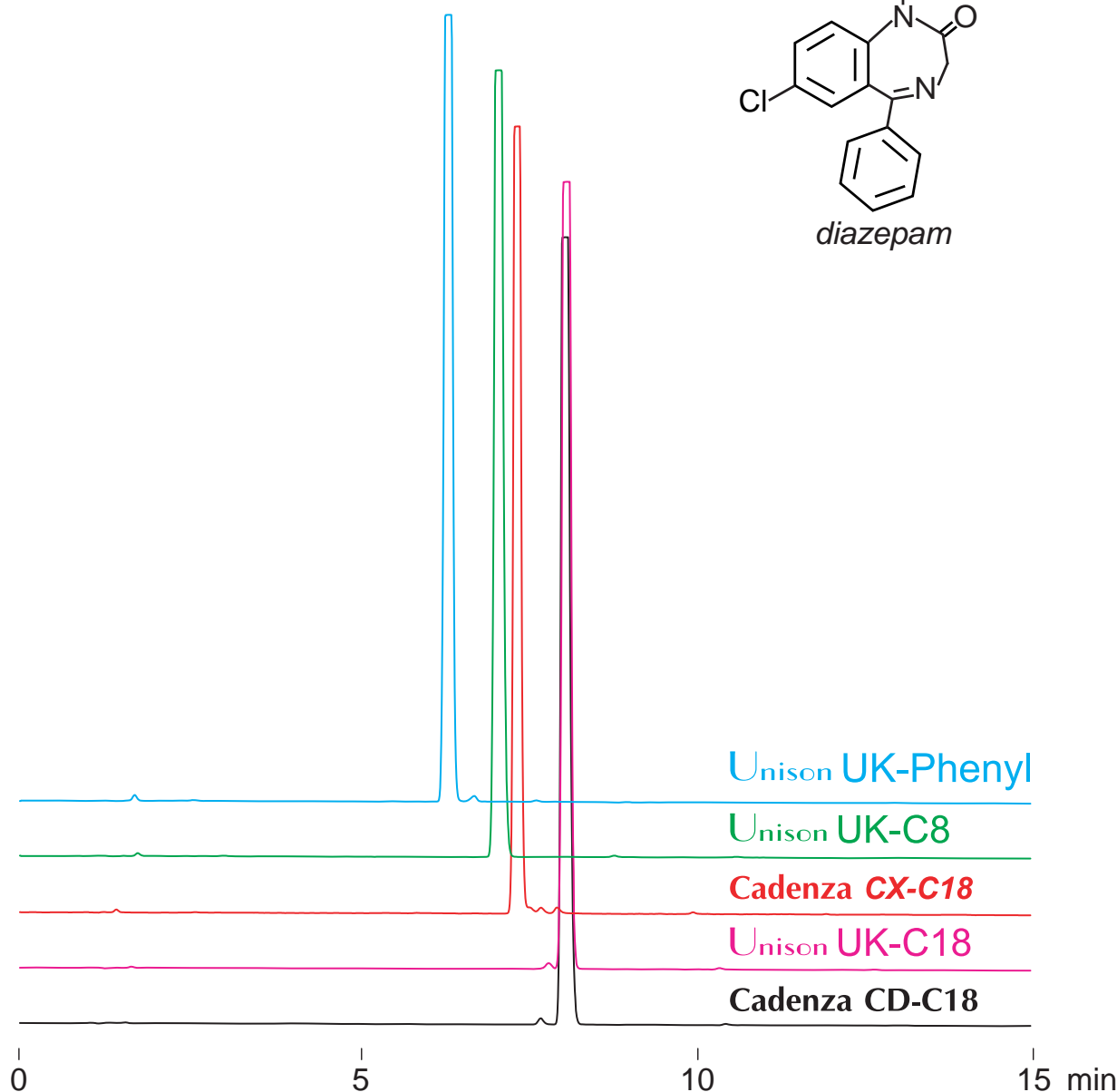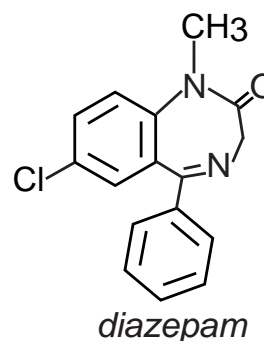

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

30-100%B (0-15min)

1.0mL/min (6MPa), 37°C, 260nm

10µL (0.5mg/mL, acetonitrile / formic acid = 100 / 1)