


Unison UK-Phenyl

150 x 3 mm


HPLC COLUMN APPLICATION

LC-MS: Bilirubin (Heme metabolite)

ビリルビン (ヘム代謝物) の LC-MS 分析

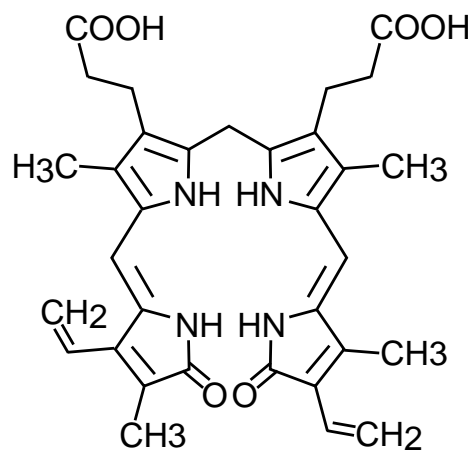

Bilirubin

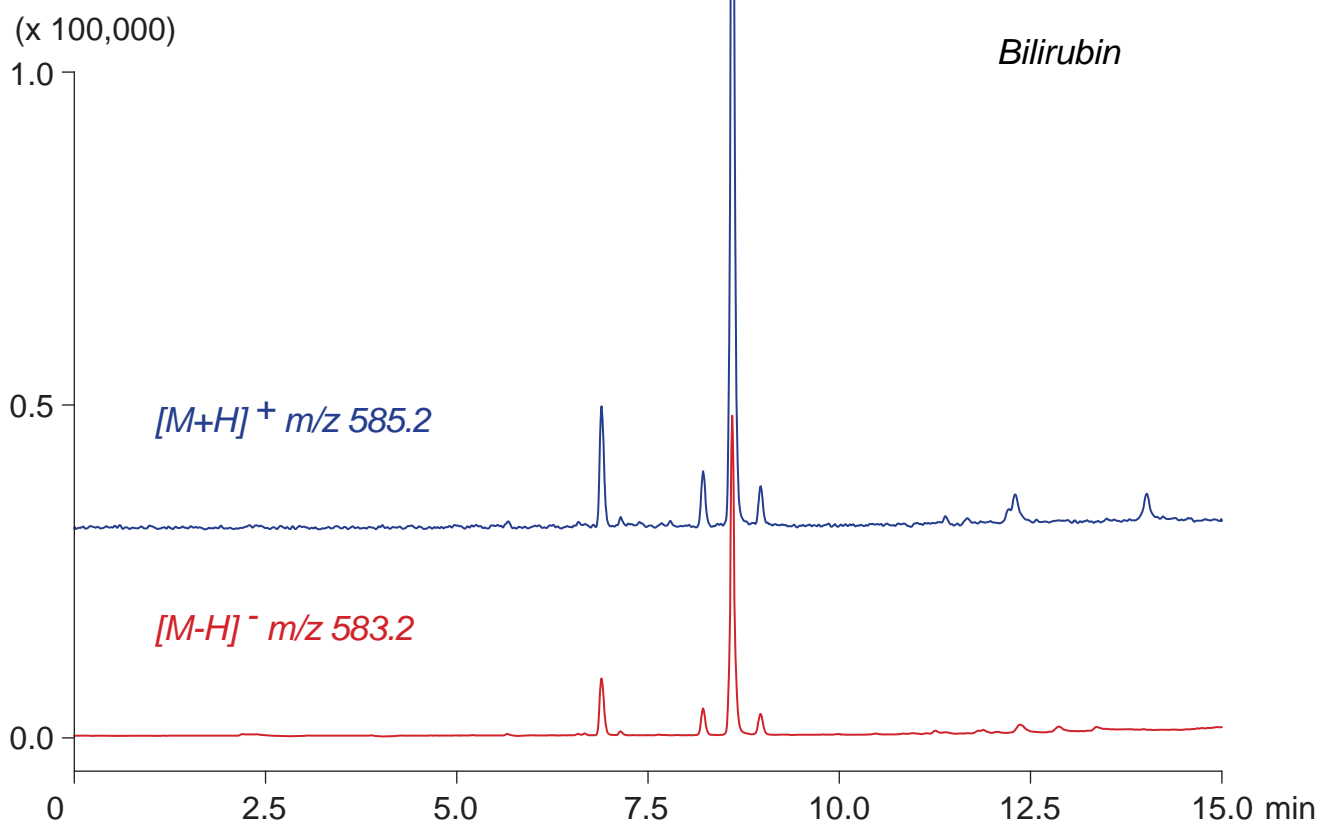

Unison UK-Phenyl, 150 x 3 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

15-100%B (0-10min)

0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)

5µL (5µg/mL, dimethylsulfoxide)