








Nardis ND-NX **MF**

75 x 2 mm

HPLC COLUMN APPLICATION

K-value (freshness index for fish) related compounds

K 値 (魚鮮度指標) 関連化合物

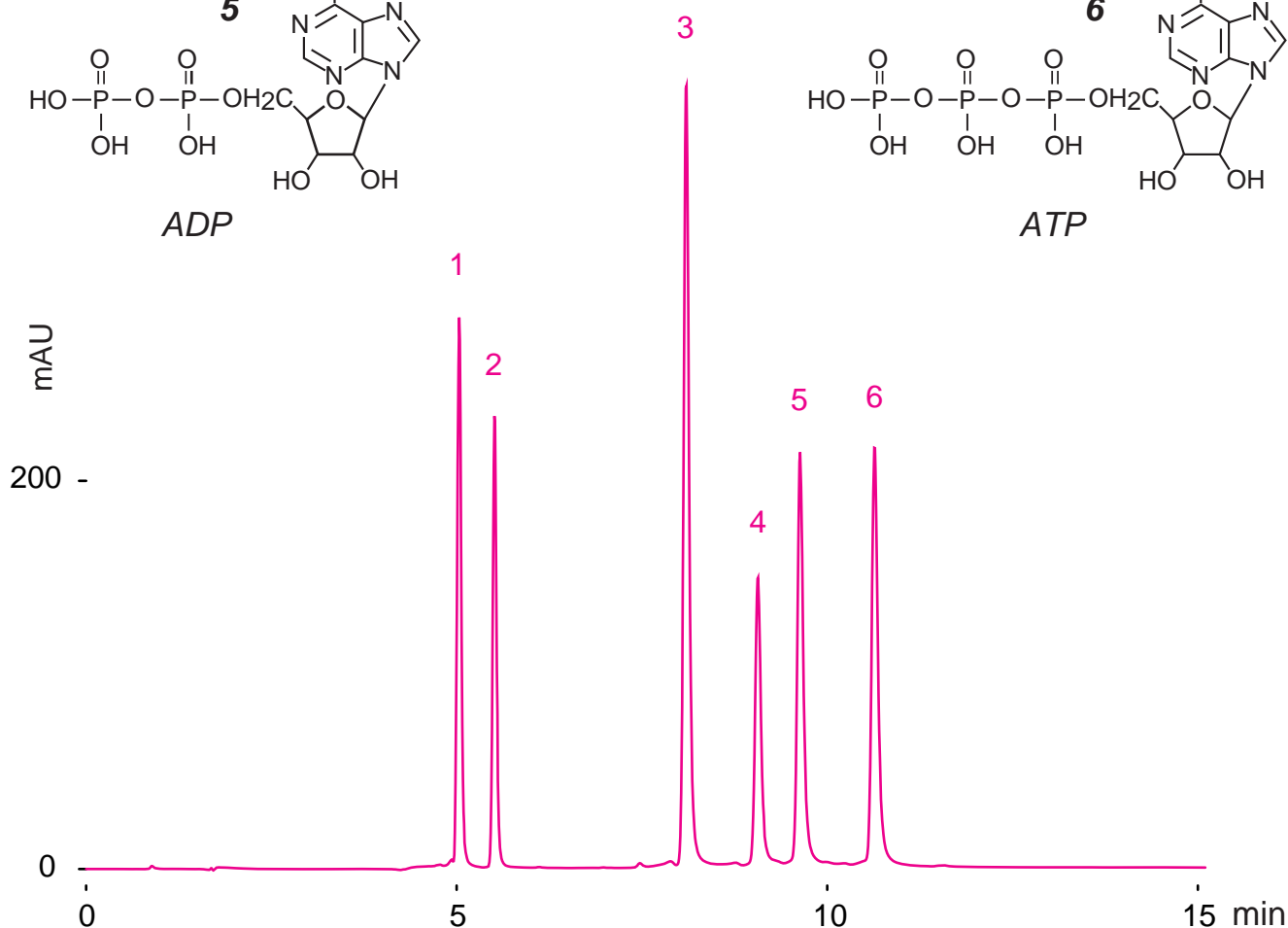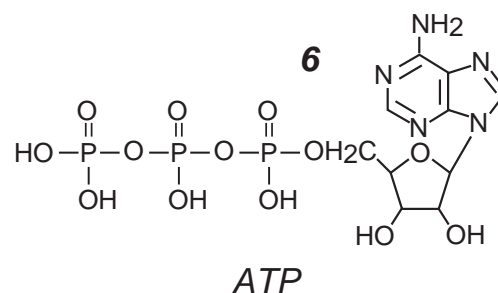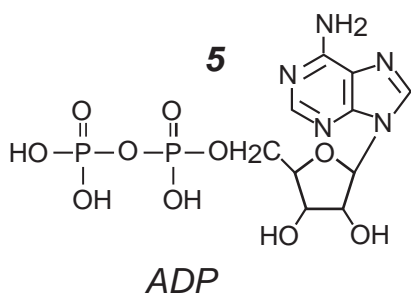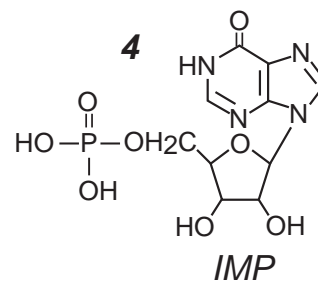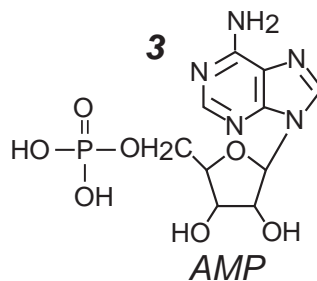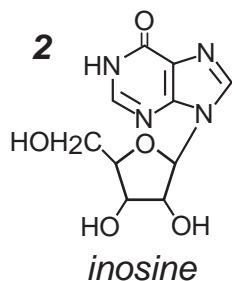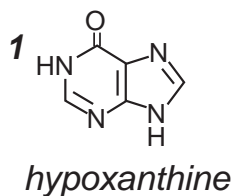

Nardis ND-NX MF, 75 x 2 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium acetate = 50 / 50

5-50%B (0-2.5min), 50-100%B (2.5-15min)

0.2mL/min (2MPa), 37°C, 260nm

3µL (0.06-0.2mg/mL, 100mM ammonium acetate)