
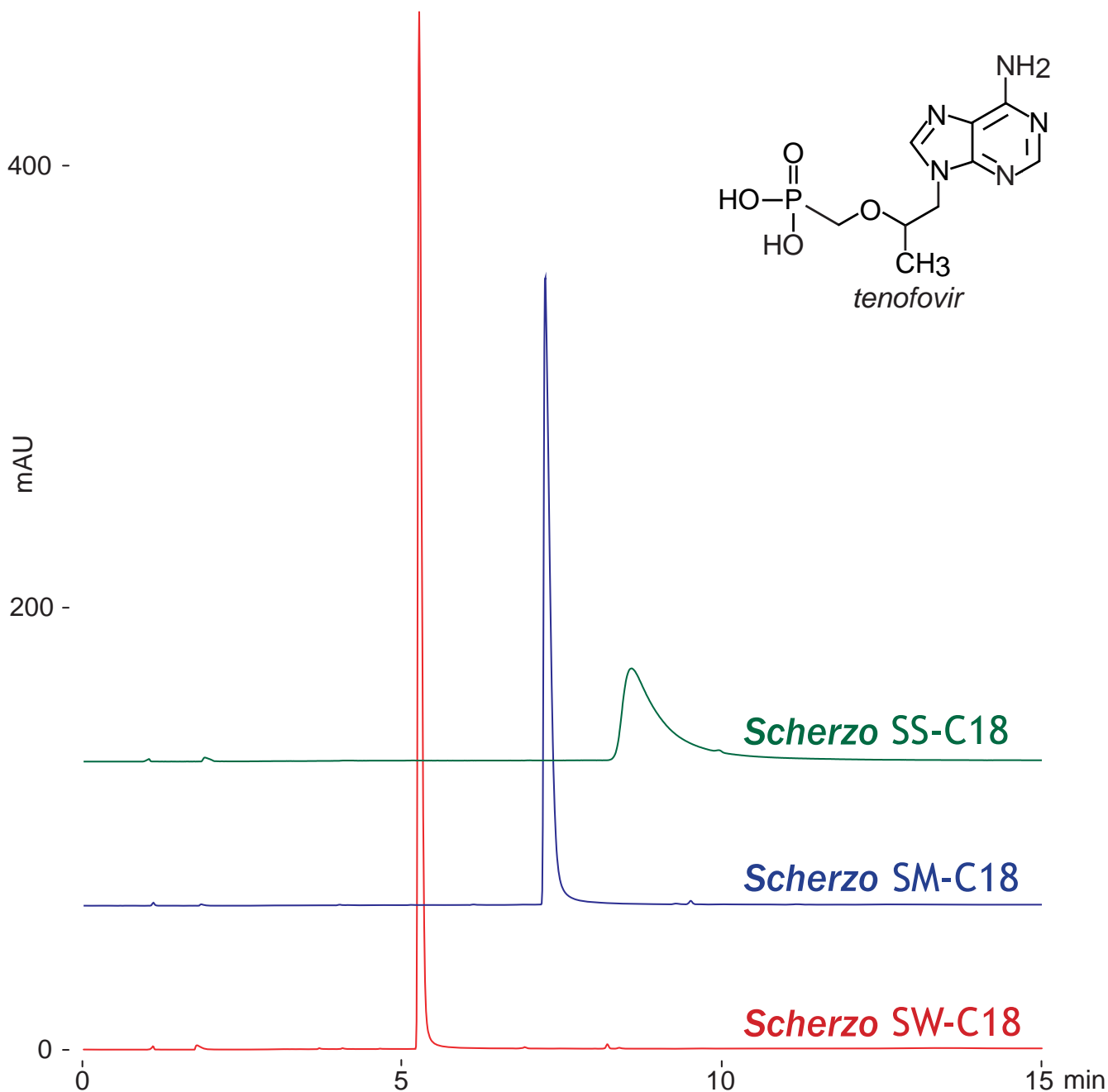


Scherzo SW-C18
 Scherzo SM-C18
 Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Tenofovir (Antiviral)
 テノホビル (抗ウイルス薬)

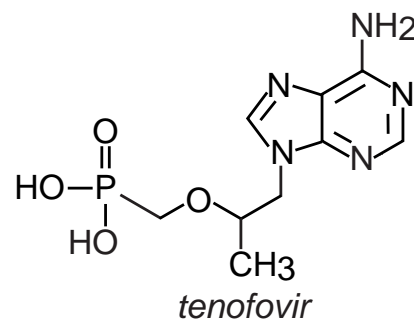

100 x 4.6 mm
 A: 10mM ammonium formate
 B: 100mM ammonium formate / acetonitrile = 50 / 50
 0-50%B (0-15min)
 1.0mL/min (8MPa), 37°C, 260nm
 10µL (0.15mg/mL, 100mM ammonium formate)