


Nardis ND-NX
Unison UK-Amino

150 x 4.6 mm

HPLC COLUMN APPLICATION

Tenofovir (Antiviral)
テノホビル (抗ウイルス薬)

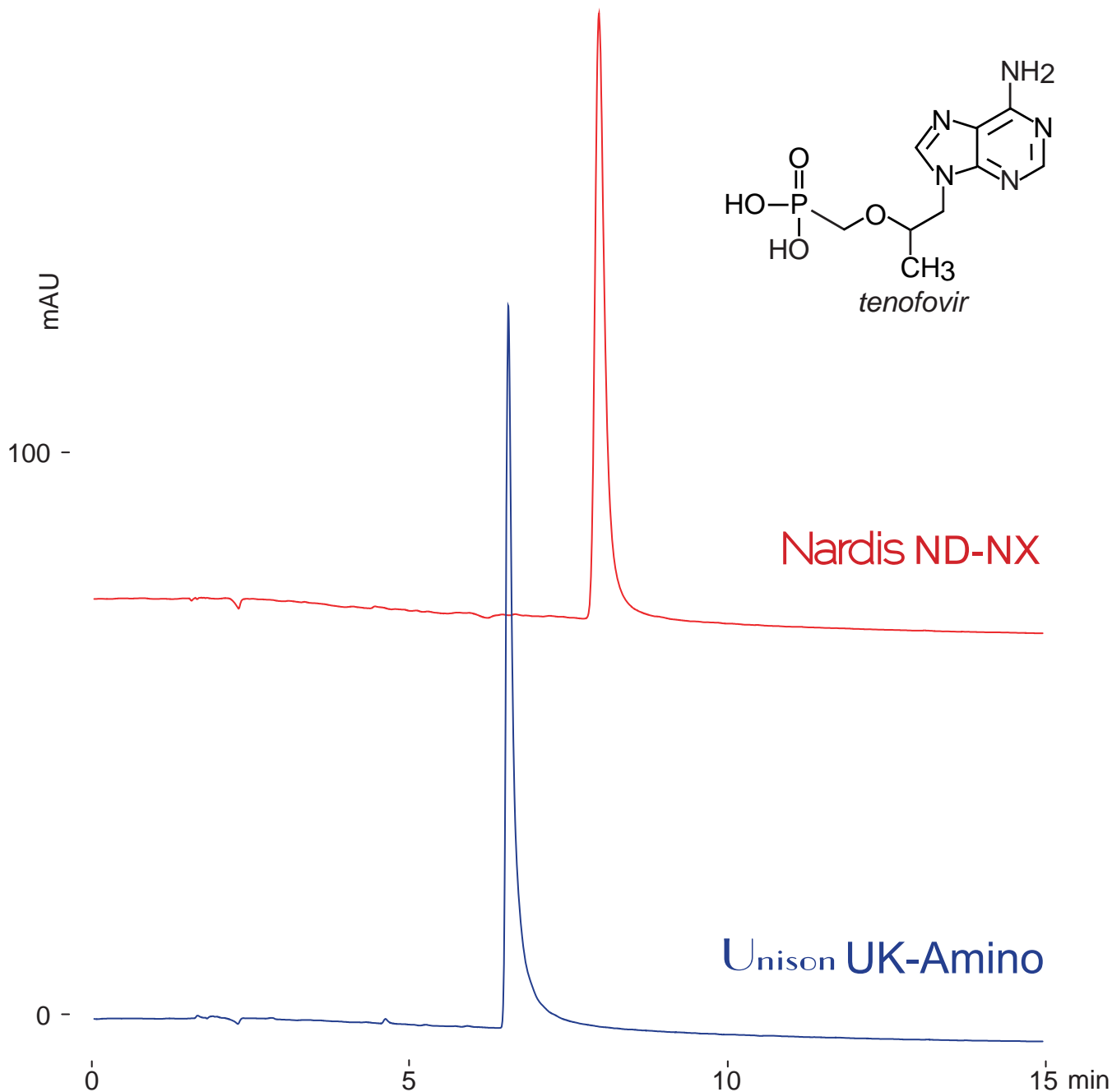

150 x 4.6 mm
 A: acetonitrile / formic acid = 100 / 0.5
 B: acetonitrile / 100mM ammonium formate = 10 / 90
 30-100%B (0-15min)
 1mL/min (6MPa), 37°C, 260nm
 10µL (0.15mg/mL, 100mM ammonium formate)