


Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

Citicoline (psychostimulant)

シチコリン (神経刺激薬) の分析

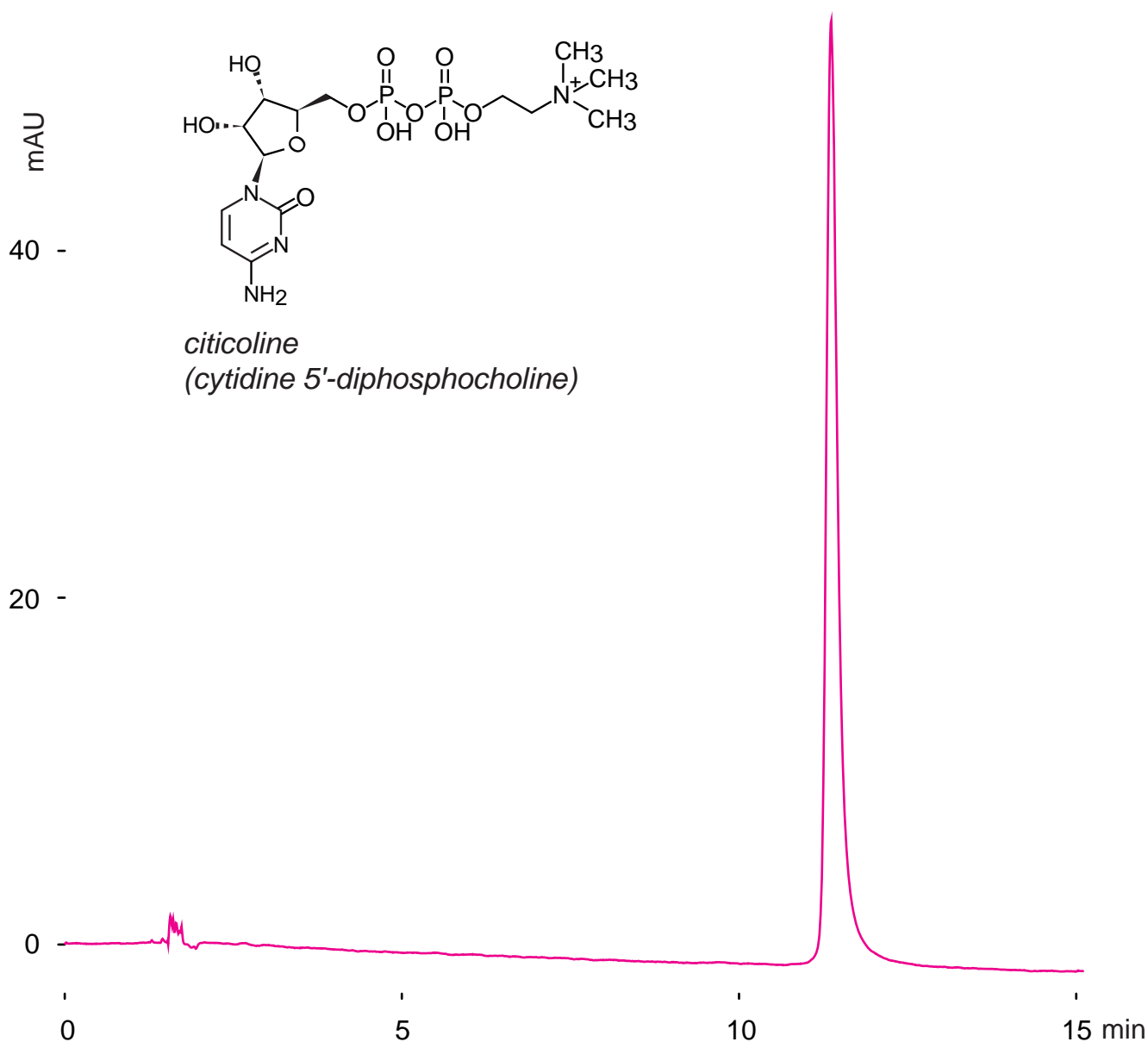

Nardis ND-NX, 150 x 4.6 mm
 A: acetonitrile / formic acid = 100 / 0.5
 B: acetonitrile / water / formic acid = 50 / 50 / 0.5
 50-100%B (0-15min)
 1.0mL/min (6MPa), 37°C, 260nm
 10µL (3mg/mL, water/HCOOH=100/1)