

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Bergapten (Furanocoumarin)

ベルガプテン (フロクマリン類)

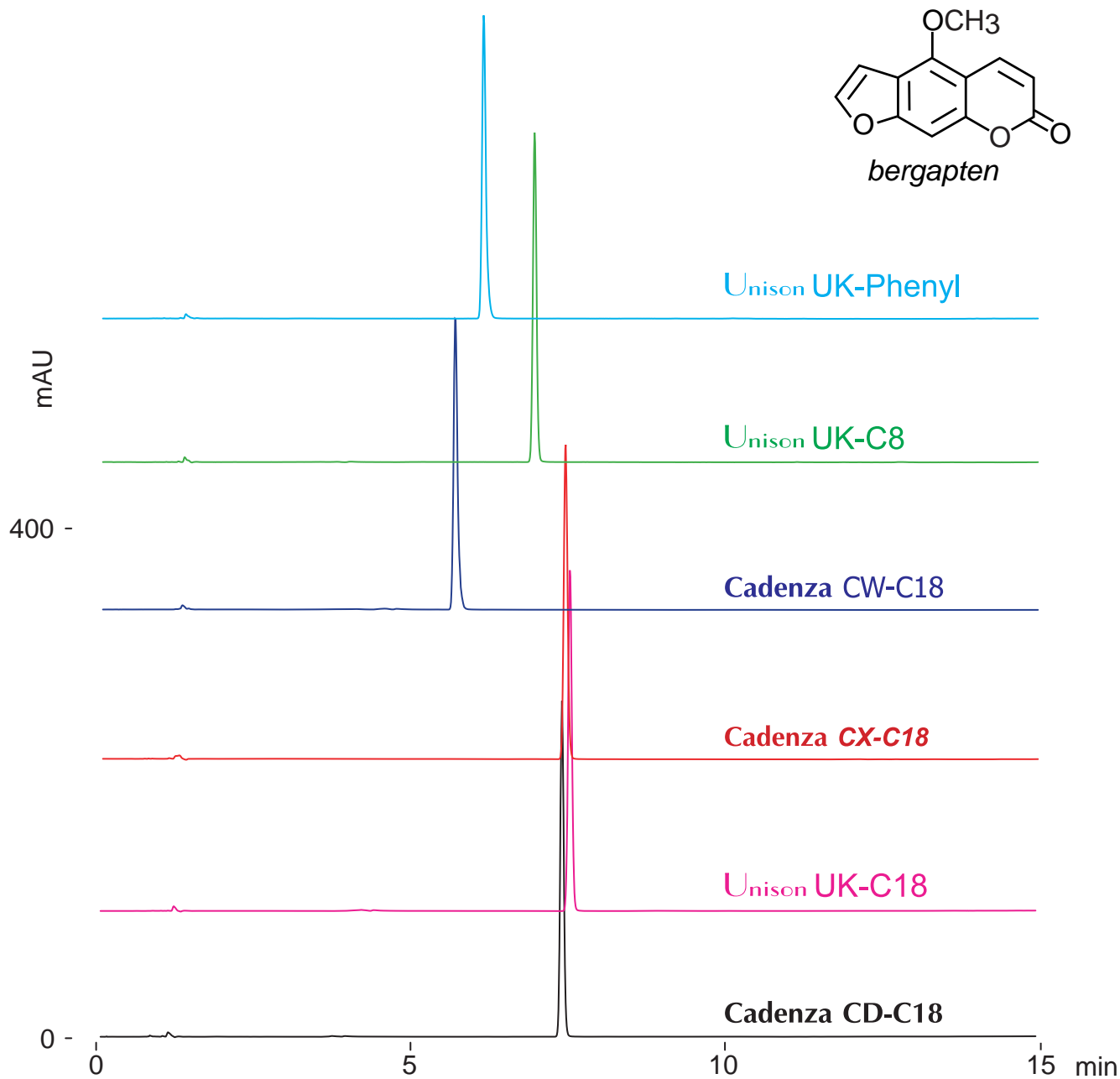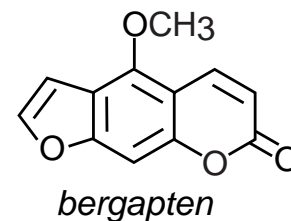

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

20-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.04mg/mL, water / acetonitrile / formic acid = 33 / 67 / 1)