


Unison UK-Phenyl Cadenza CL-C18  
 Unison UK-C8 Cadenza CW-C18  
 Cadenza CX-C18 Unison UK-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

### Helicide (Phenolic glycosides)

4-ホルミルフェニル -D-アロピラノシド (フェノール性配糖体)

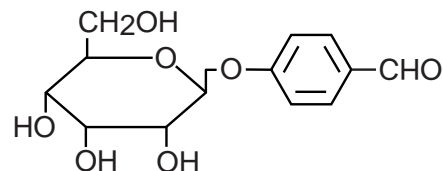

4-formylphenyl beta-D-allopyranoside  
(helicide)

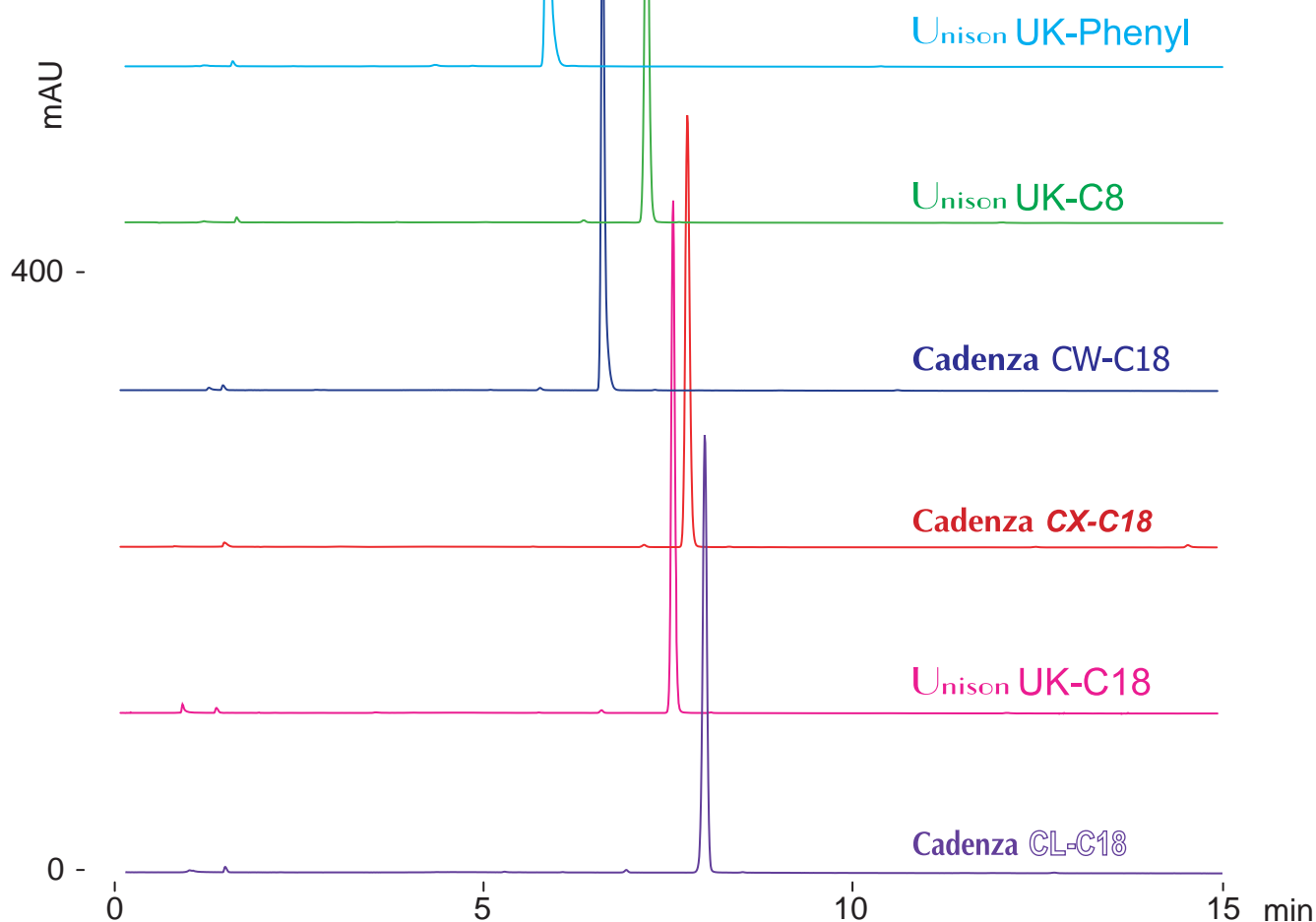

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 50 / 50 / 0.1

0%B (0-1min), 0-40%B (1-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.1mg/mL, water / formic acid = 100 / 1)