

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Camptothecin (Alkaloid)

カンプトテシン (アルカロイド)

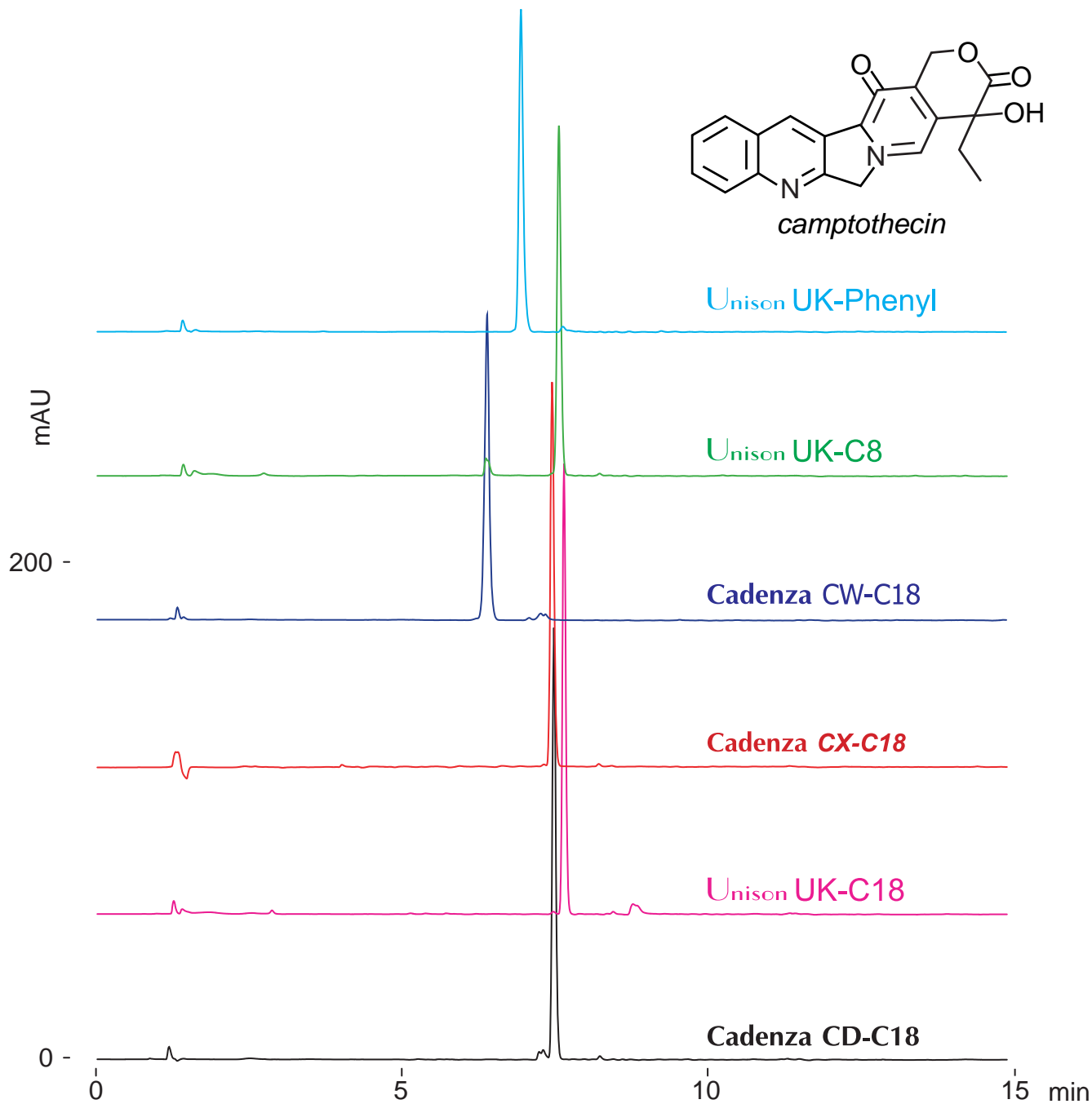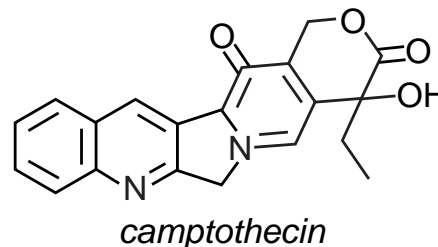

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-85%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.033mg/mL, water / acetonitrile / formic acid = 20 / 80 / 1)