


Unison UK-Phenyl Unison UK-C1

Unison UK-C18 Unison UK-C8 100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Dioctyl sulfosuccinate (Double-tail anionic surfactant)

ジオクチルスルホコハク酸 (二鎖型アニオン界面活性剤) の LC-MS 分析

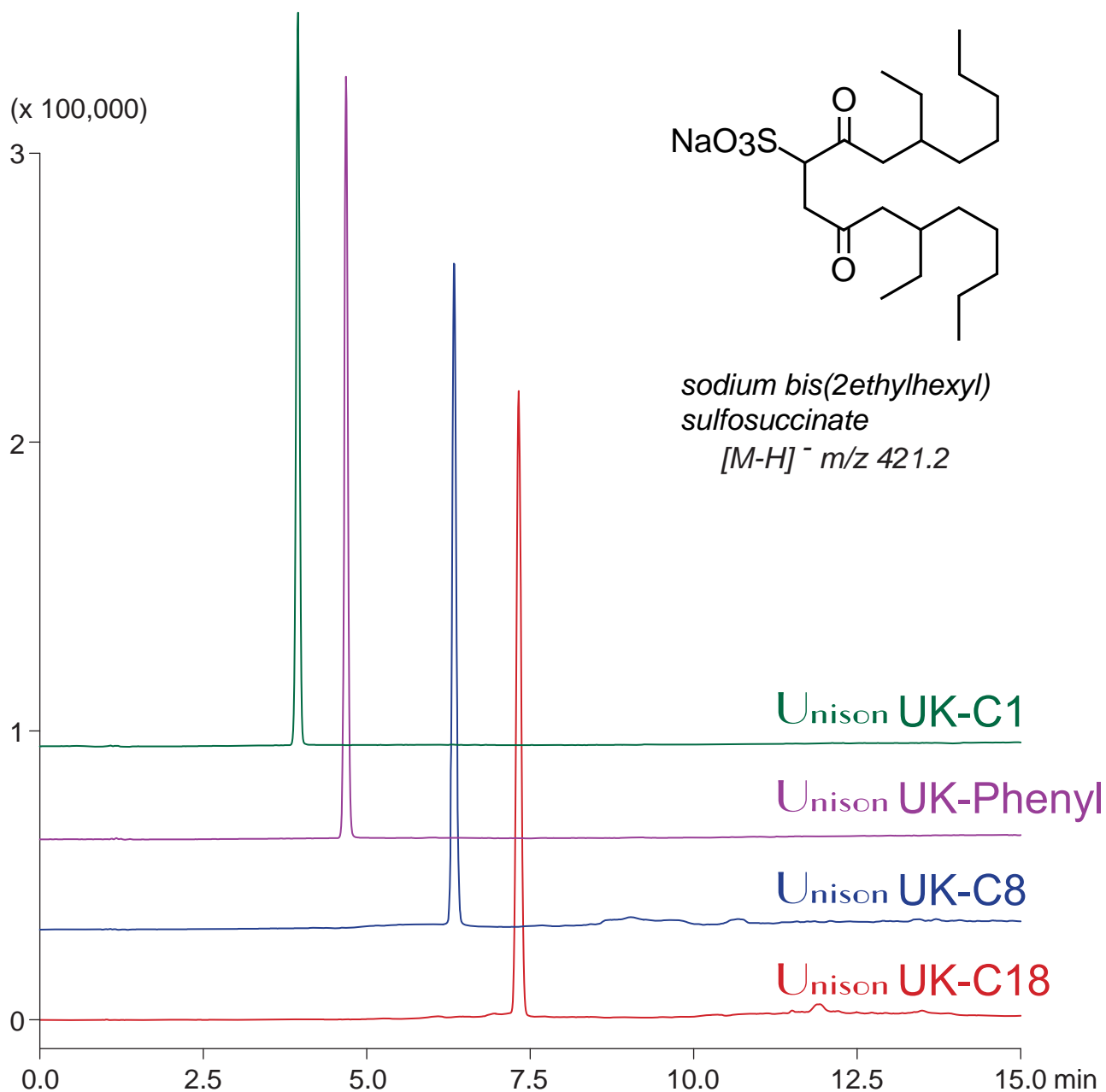

100 x 3 mm
 A: acetonitrile / 100mM ammonium formate = 40/ 60
 B: acetonitrile / 100mM ammonium formate = 90/ 10
 0-100%B (0-10min), 100%B (10-15min)
 0.4mL/min (6-7MPa), 37°C, Single Quad. MS (ESI)
 5uL (50ng/mL, ACN / water / HCOOH = 50 / 50 / 1)