


Cadenza CX-C18

150 x 3 mm


HPLC COLUMN APPLICATION

Mirtazapine (Antidepressant)

ミルタザピン (抗うつ薬)

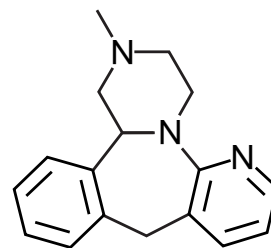

mirtazapine

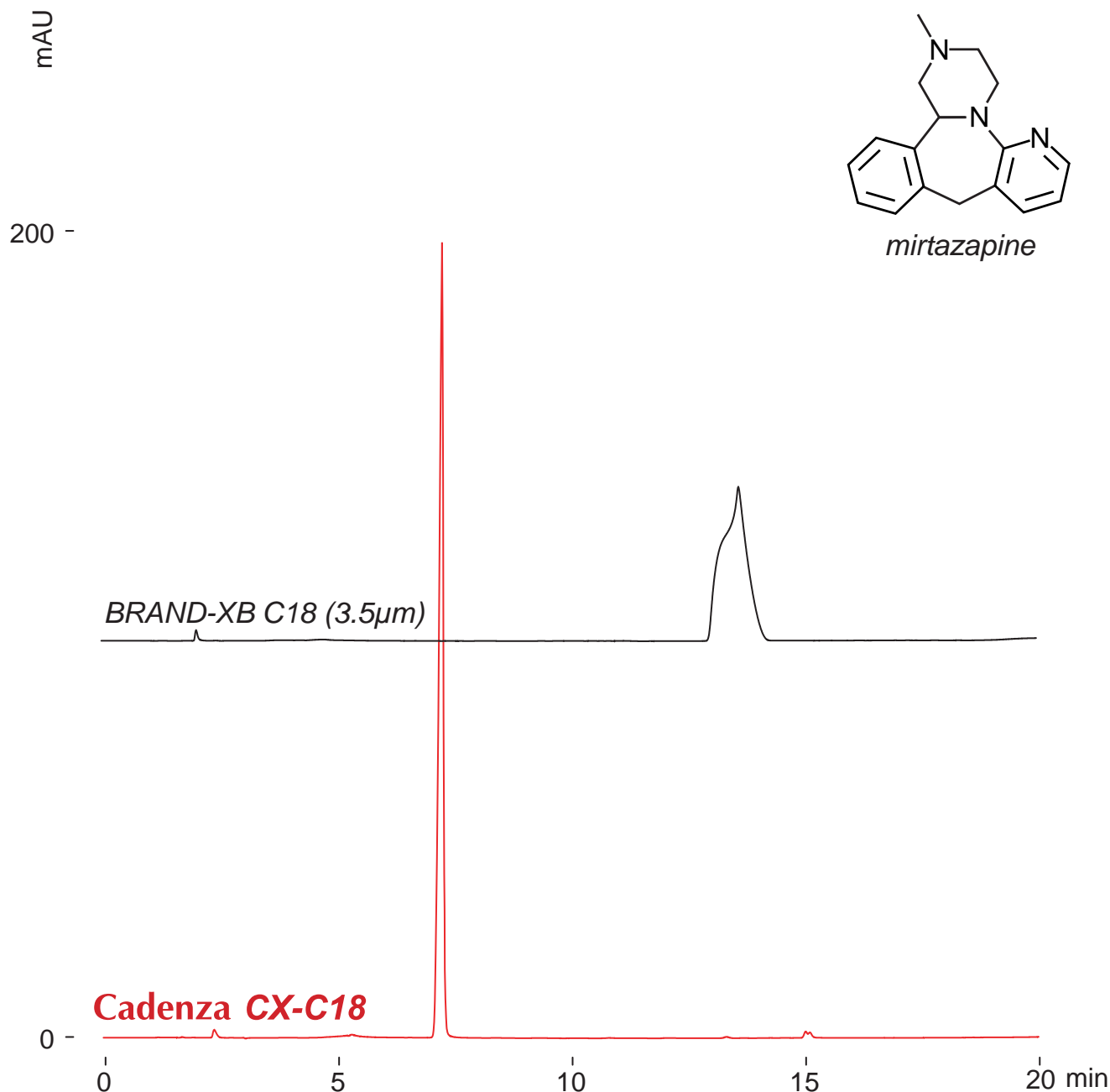

150 x 3 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 70 / 30 / 0.1

0%B (0-1min), 0-100%B (1-20min)

0.4mL/min (9MPa), 37°C, 260nm

2µL (0.7mg/mL, water / formic acid = 100 / 1)