


Nardis ND-NX MF

150 x 2 mm

HPLC COLUMN APPLICATION

Nucleoside triphosphate

ヌクレオシド三リン酸の分析

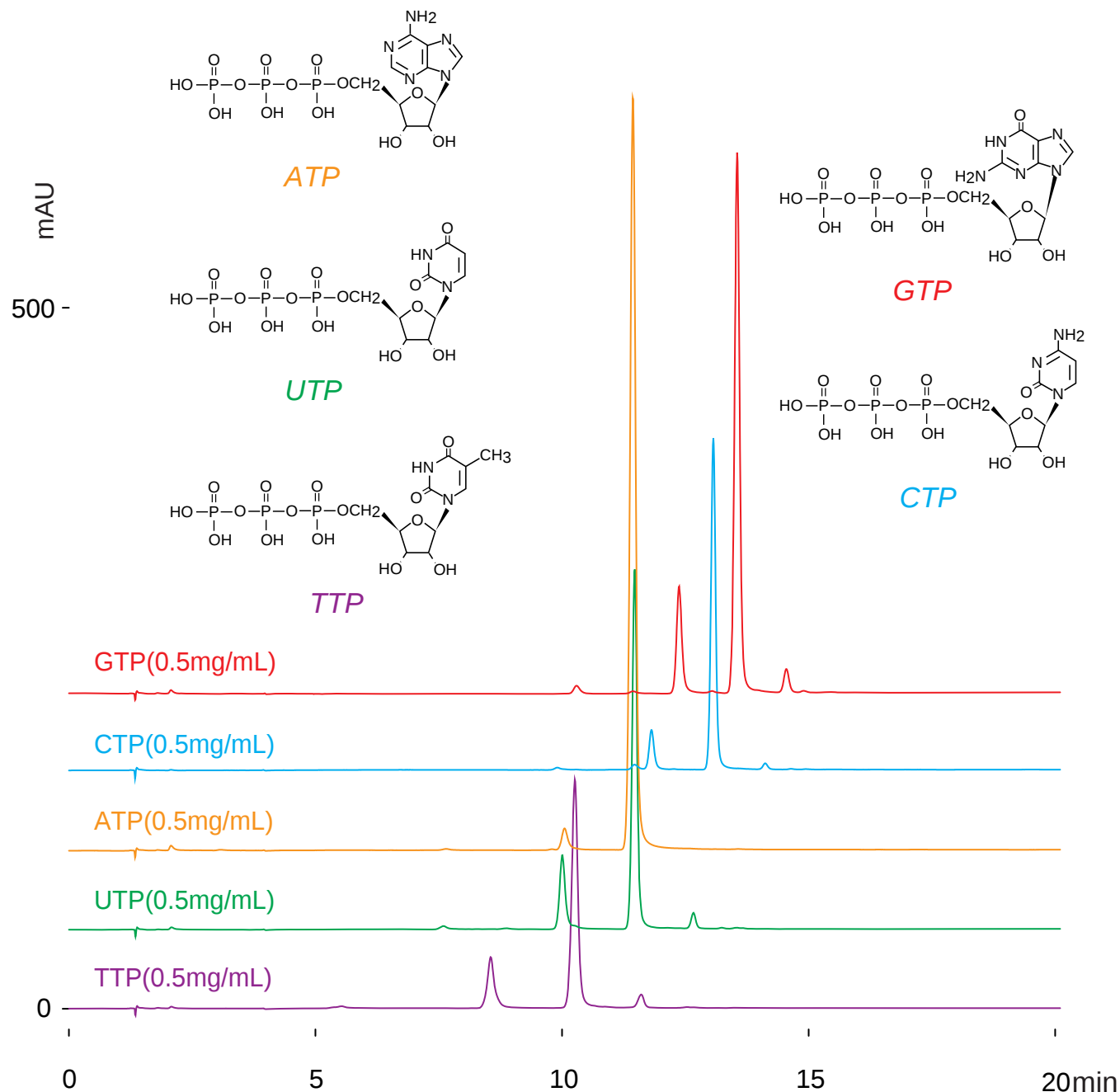

Nardis ND-NX MF, 150 x 2 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium formate = 50 / 50

50 -100%B (0 -15 min), 100%B (15 -20 min)

0.2 mL/min (6MPa), 37°C, 260nm

3 μ L (0.5mg/mL, water/HCOOH=100/1)