

Scherzo SM-C18

Scherzo SW-C18


Nardis ND-RX

100 x 4.6 mm

HPLC COLUMN APPLICATION

Dioctyl slufosuccinate (Duble-tail anionic surfactant)

ジオクチルスルホコハク酸 (二鎖型アニオン界面活性剤) の分析

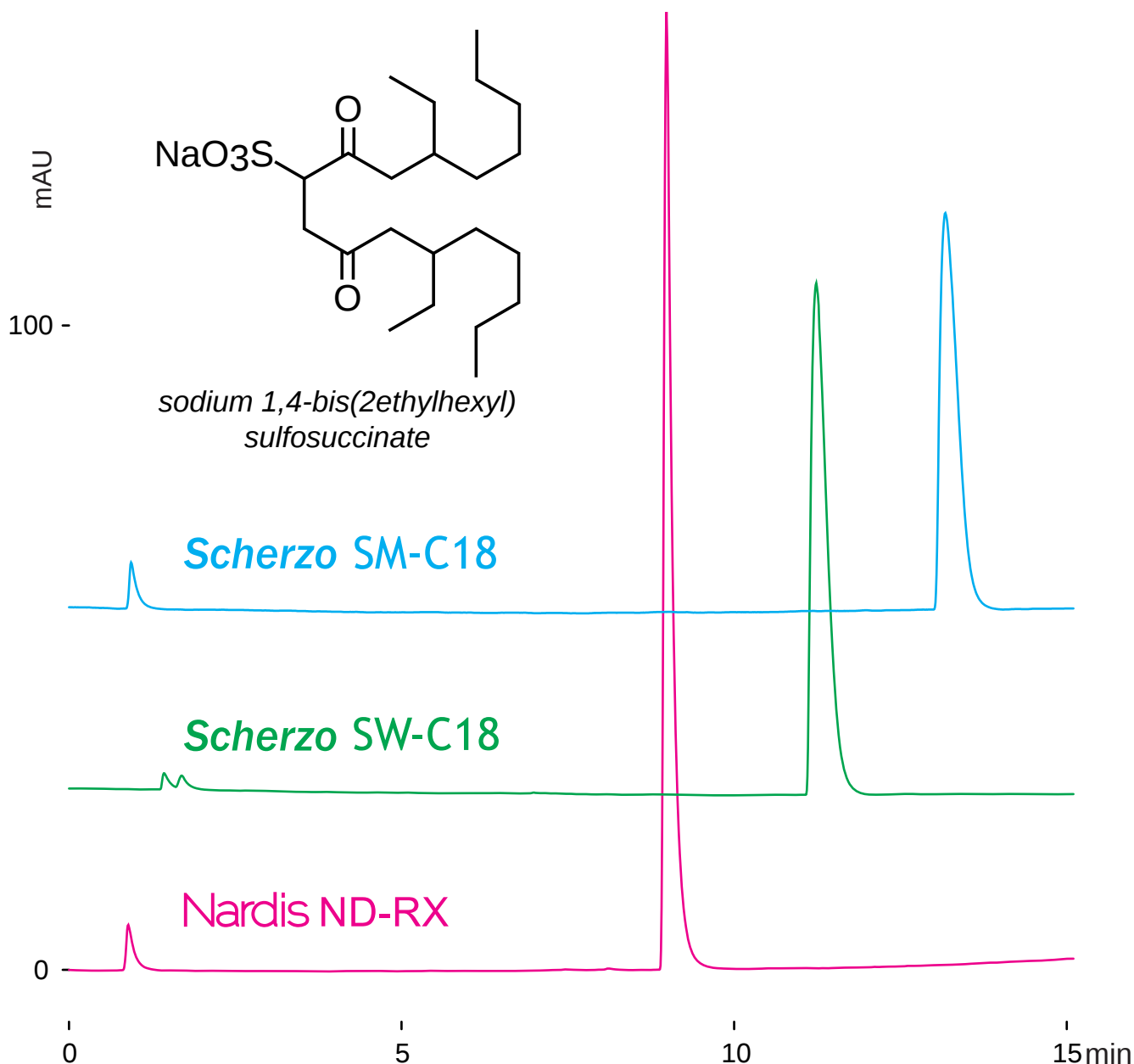

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: 100mM ammonium formate / acetonitrile = 25 / 75

0-100%B (0-10min), 100%B (10-15min)

1.0mL/min (8-10MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10μL (2mg/mL, water/acetonitrile/formic acid=50/50/1)