

Unison UK-Amino


Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

Glucuronolactone, glucuronic acid

グルクロノラクトンとグルクロン酸の分析

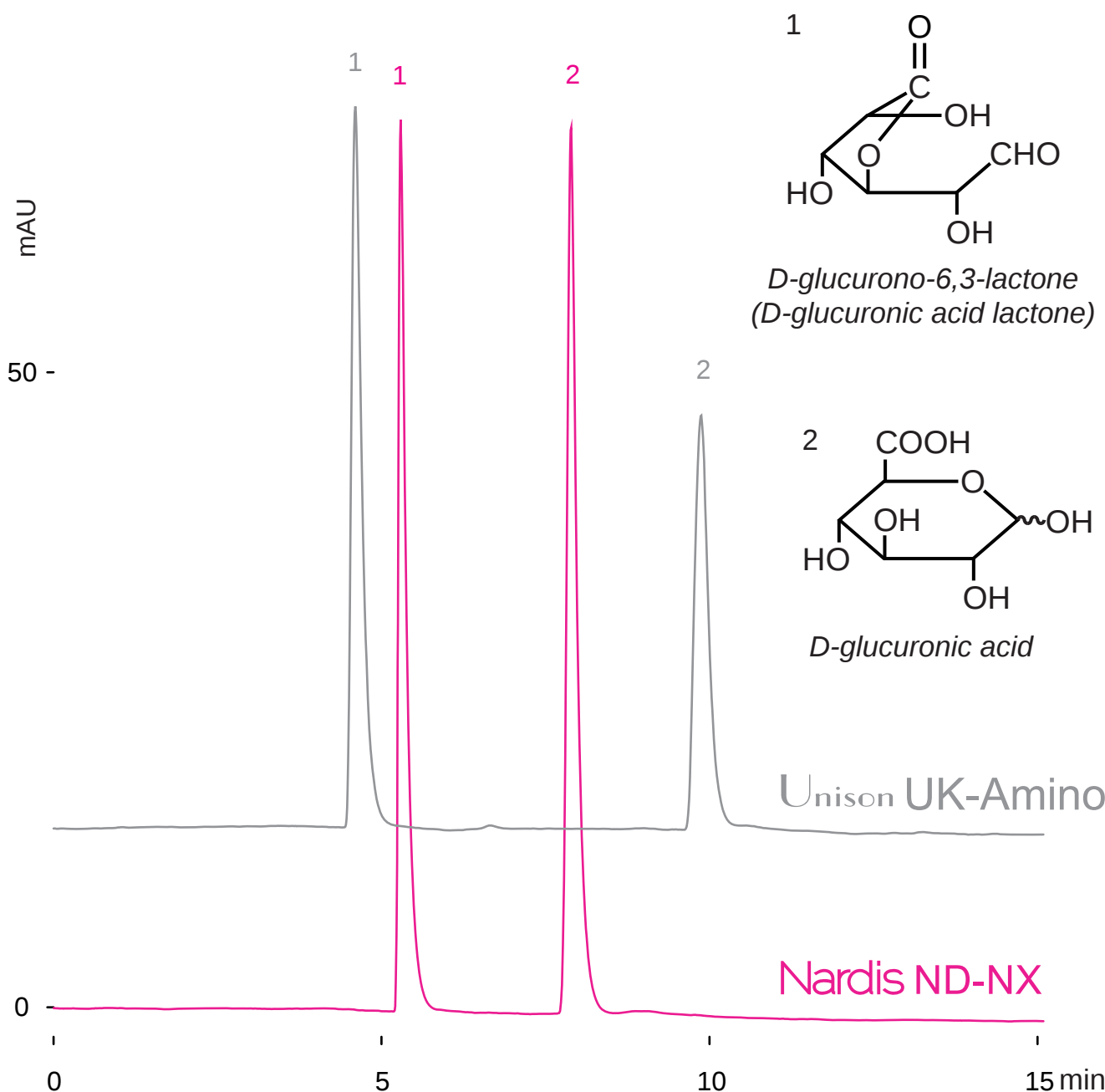

150 x 4.6 mm

A: acetonitrile

B: 100mM ammonium acetate

0 -100%B (0 -15 min), 1.0 mL/min (5 -6 MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10 µL (1-2mg/mL, water/AcOH=100/0.1)