



Unison UK-Phenyl
Unison UK-C8
Nardis ND-RX
Cadenza CD-C18
Unison UK-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Valacyclovir

バラシクロビルの分析

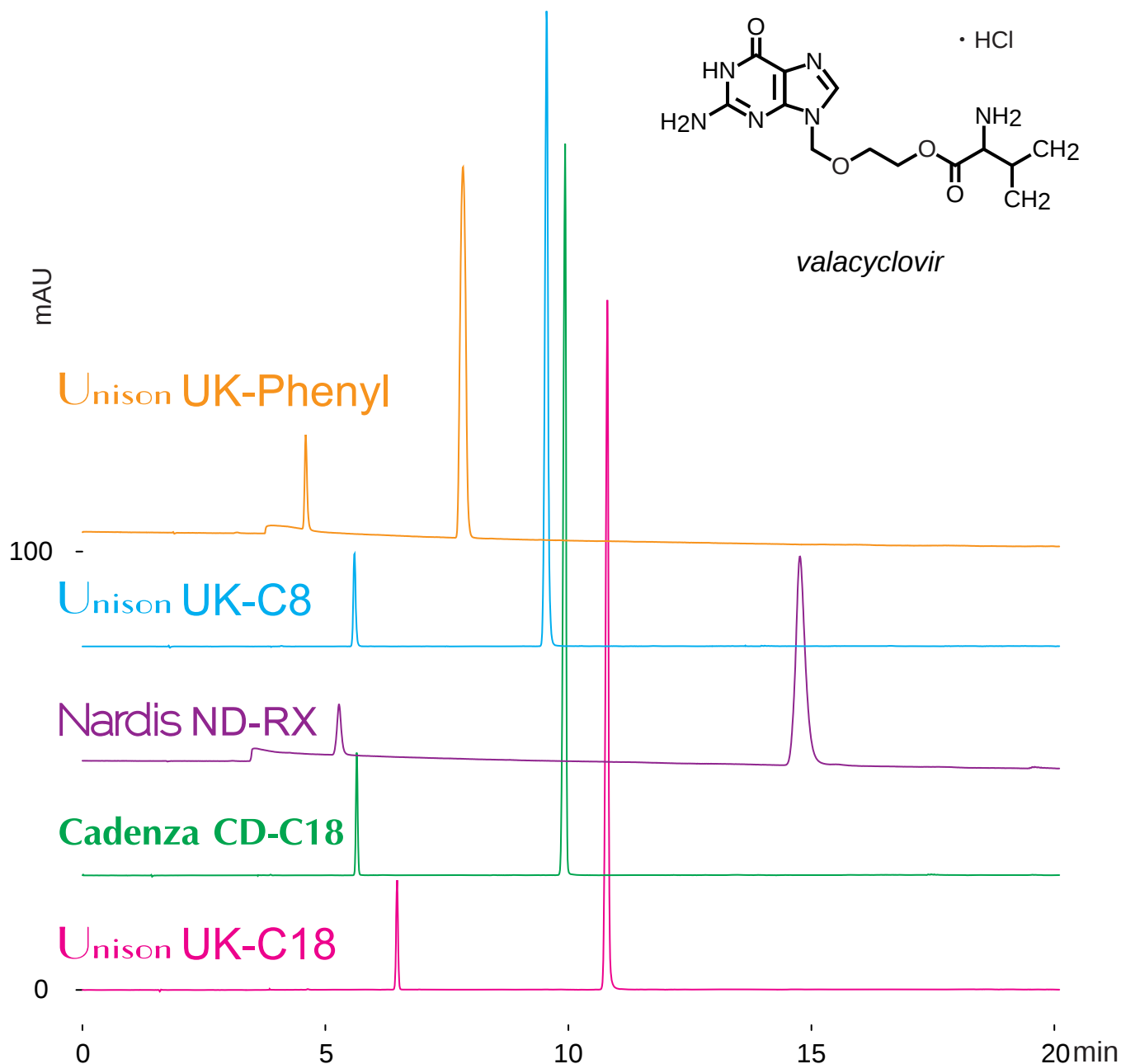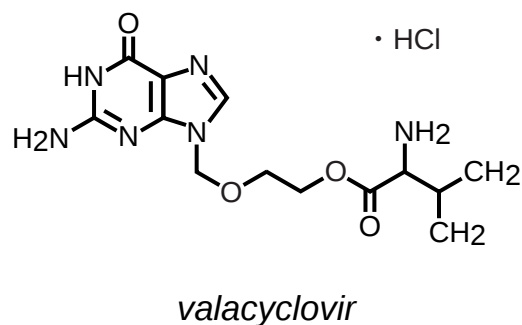

150 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-30%B (0-15min), 30%B (15-20min)

1.0mL/min (9-10MPa), 37°C, 260nm

10μL (0.1mg/mL, water)