


Unison UK-Phenyl
Unison UK-C8
Nardis ND-RX
Cadenza CD-C18
Unison UK-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Vidarabine


ビダラビンの分析

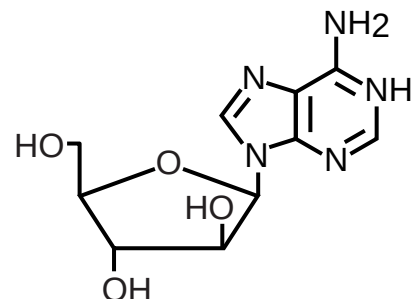

vidarabine

(9-beta-D

-arabinofuranosyladenine)

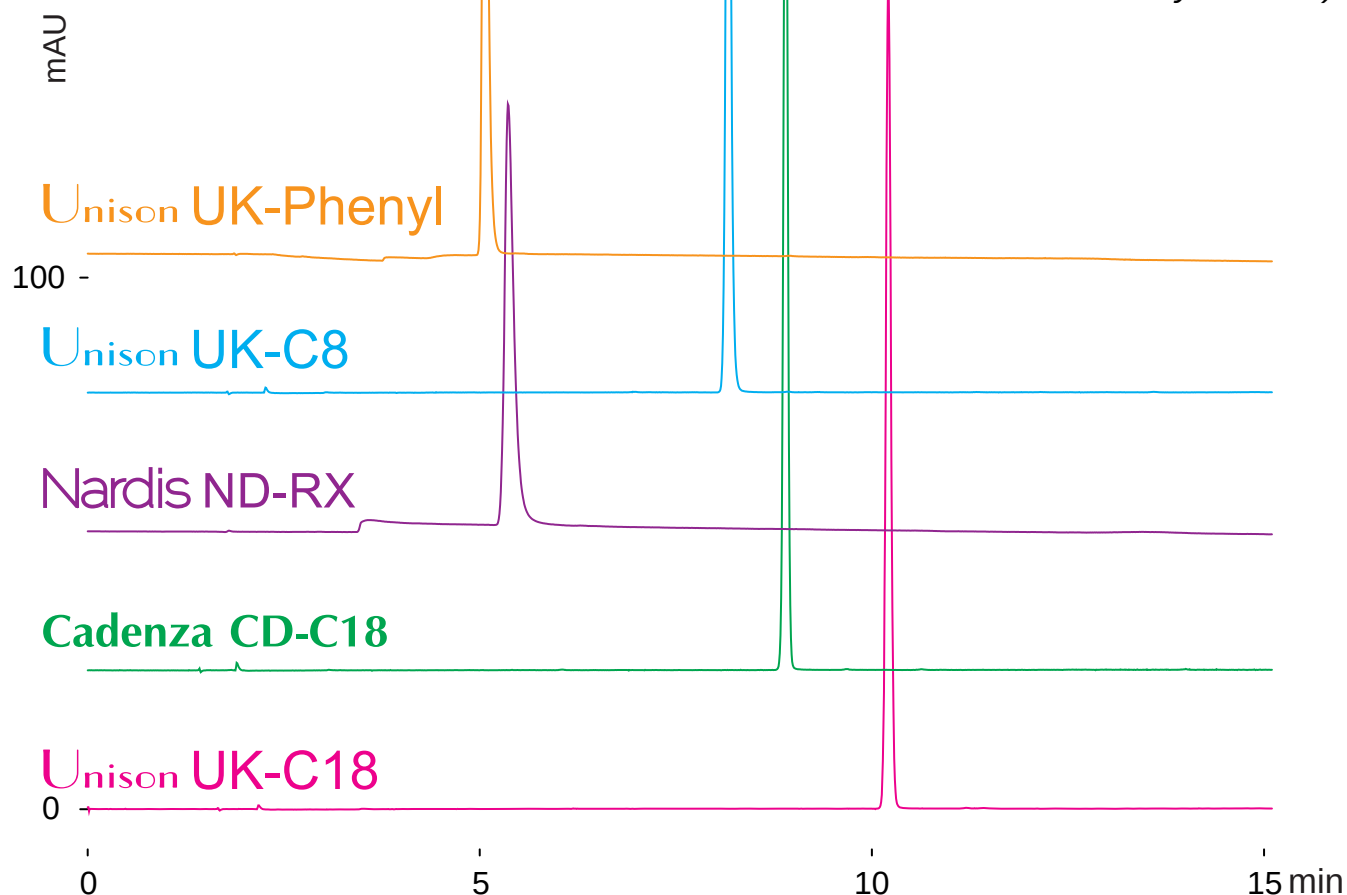

150 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-20%B (0-15min)

1.0mL/min (9-10MPa), 37°C, 260nm

10µL (50µg/mL, water/formic acid=100/1)