

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Enzalutamide (Antineoplastic) and synthetic intermediate

エンザルタミド（抗悪性腫瘍薬）と合成中間体

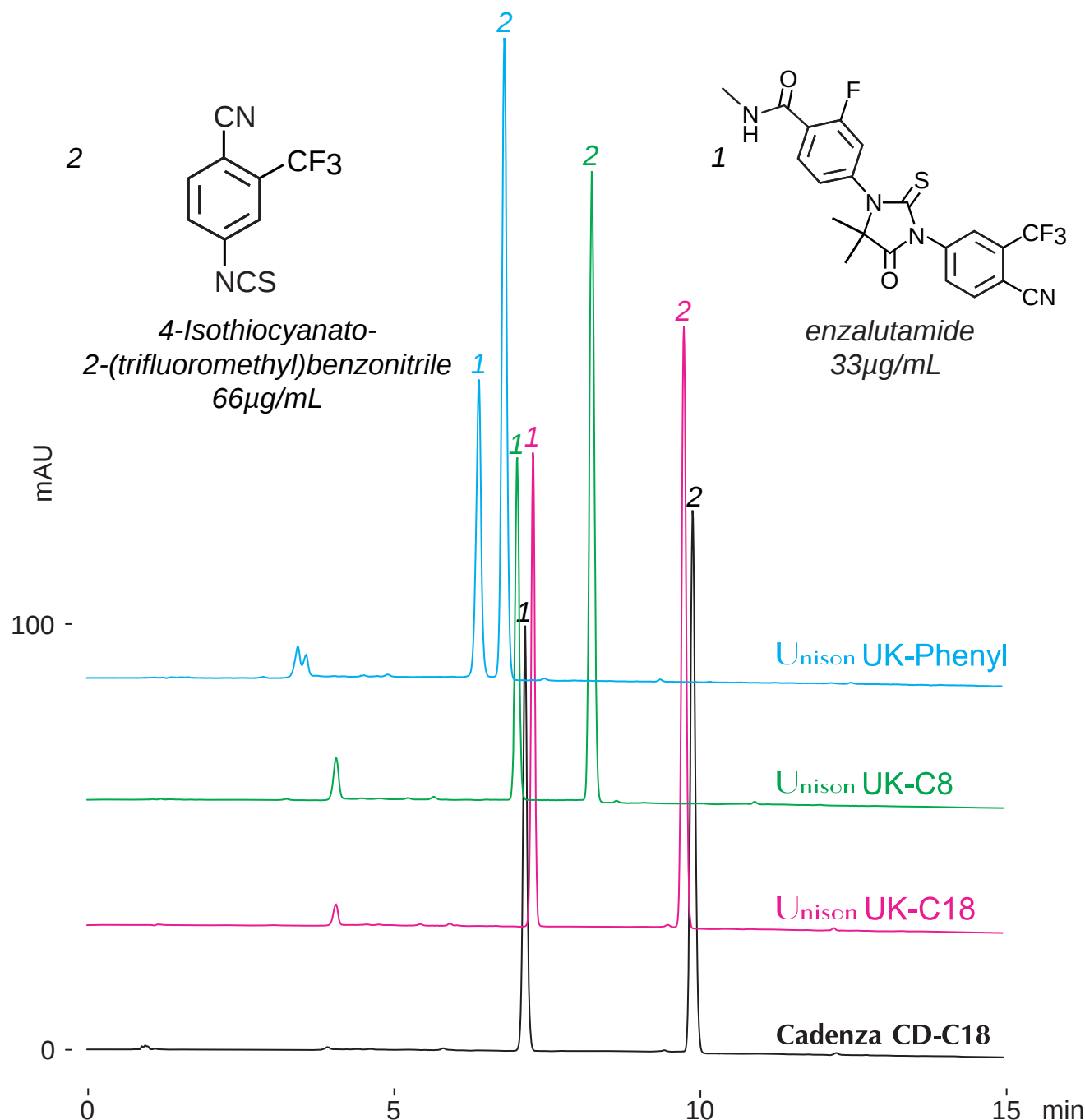

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

30-100%B (0-15min)

1.0mL/min (11MPa), 37°C, 260nm

10 µL (33-66 µg/mL, acetonitrile)