


Unison UK-Amino  
Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

## Biological buffer components

生化学用緩衝液成分の分析

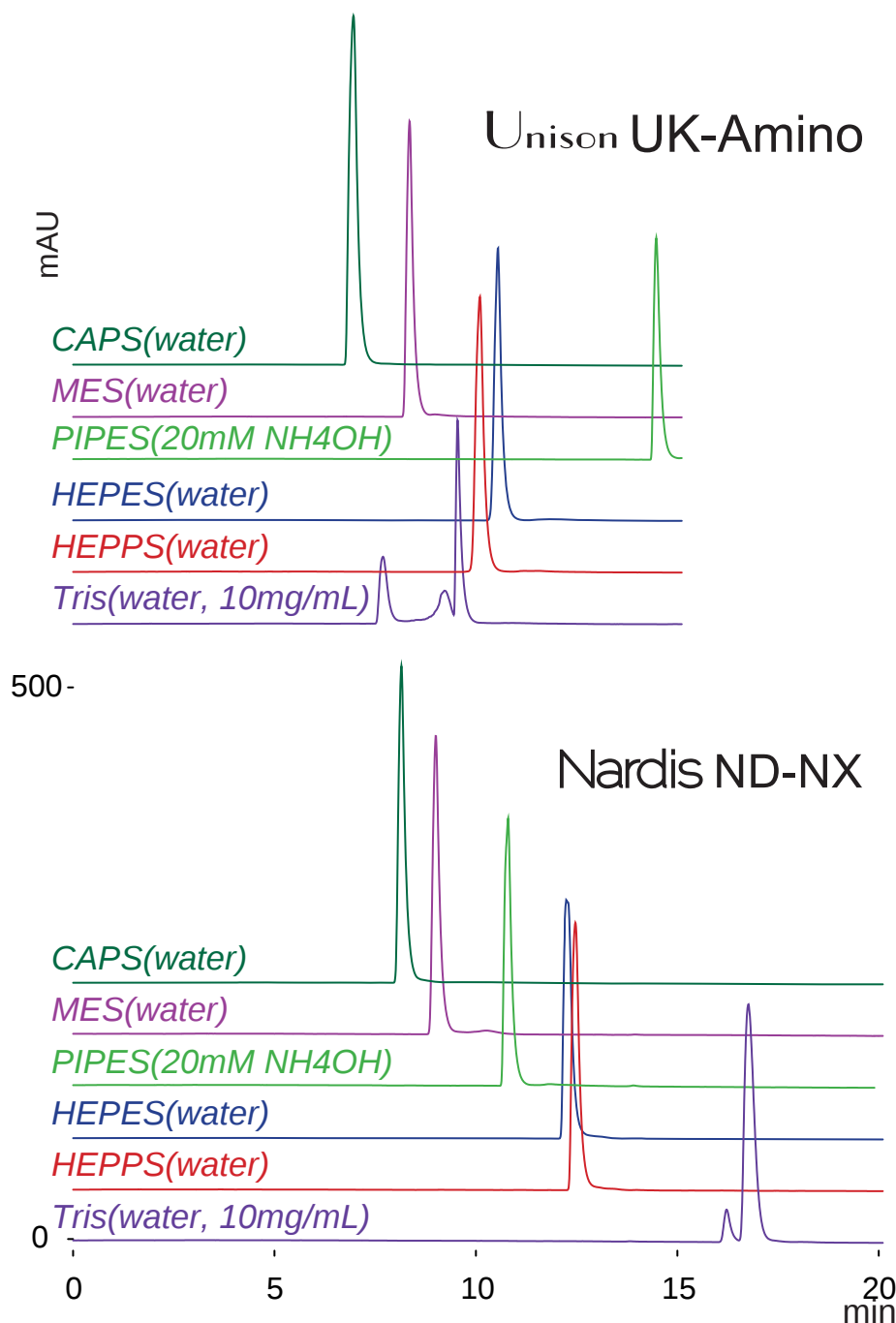

150 x 4.6 mm

A: acetonitrile / formic acid = 100 / 0.1

B: acetonitrile / 100mM ammonium formate = 50 / 50

0 -100%B (0 -15 min), 100%B (15 -20 min)

1.0 mL/min (6-8MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10 µL (5-10mg/mL)

