

Cadenza CW-C18


Presto FF-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

Span 20 (sorbitan monolaurate)

スパン 20(ソルビタンモノラウレート) の分析

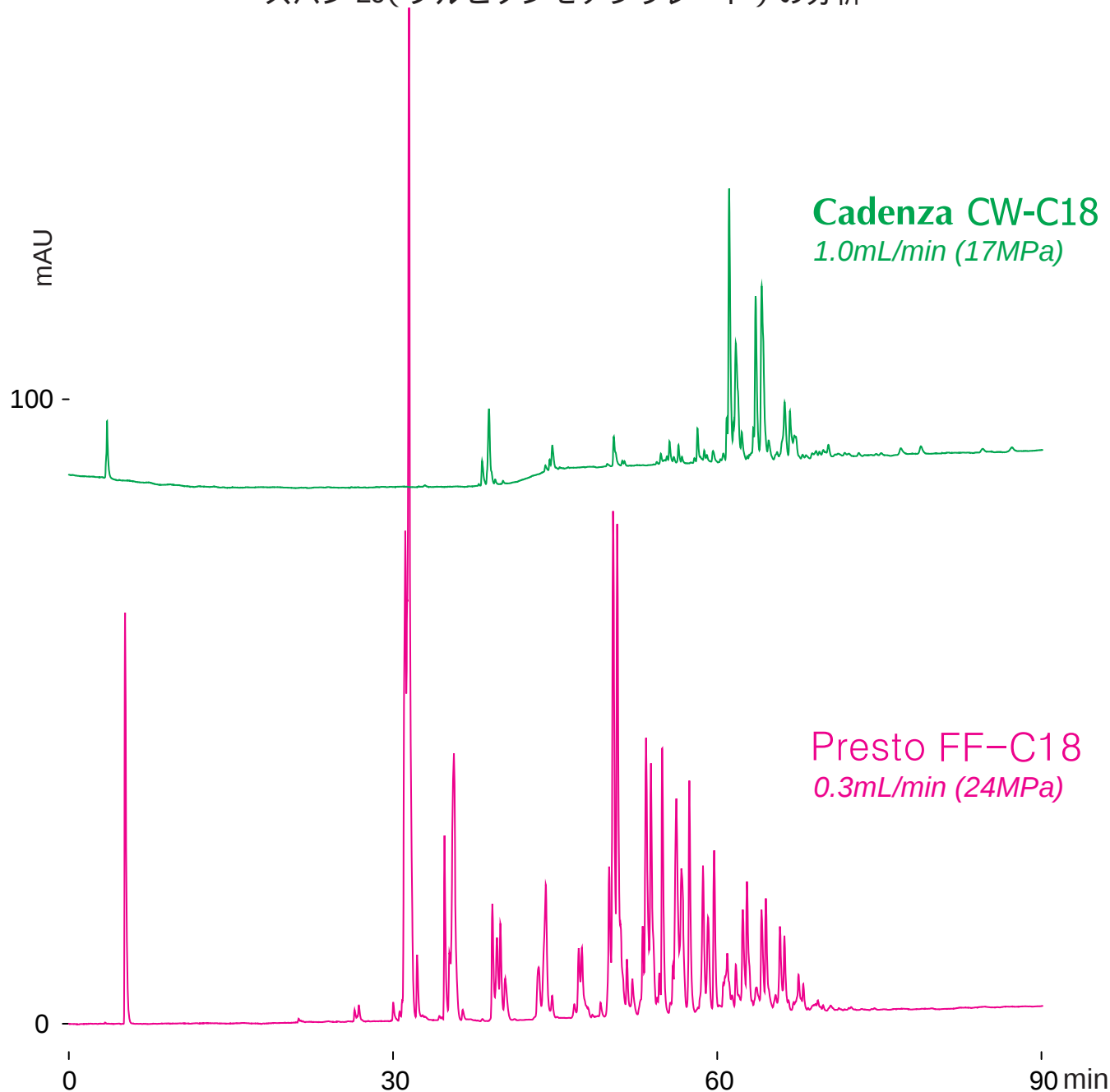

250 x 4.6 mm

A: water / acetic acid = 100 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

0-100%B (0-60min), 100%B (60-90min)

37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10μL (20mg/mL, methanol)