


Unison UK-Amino

150 x 4.6 mm

HPLC COLUMN APPLICATION

Hydroxypicolinic acid isomers

ヒドロキシピコリン酸異性体

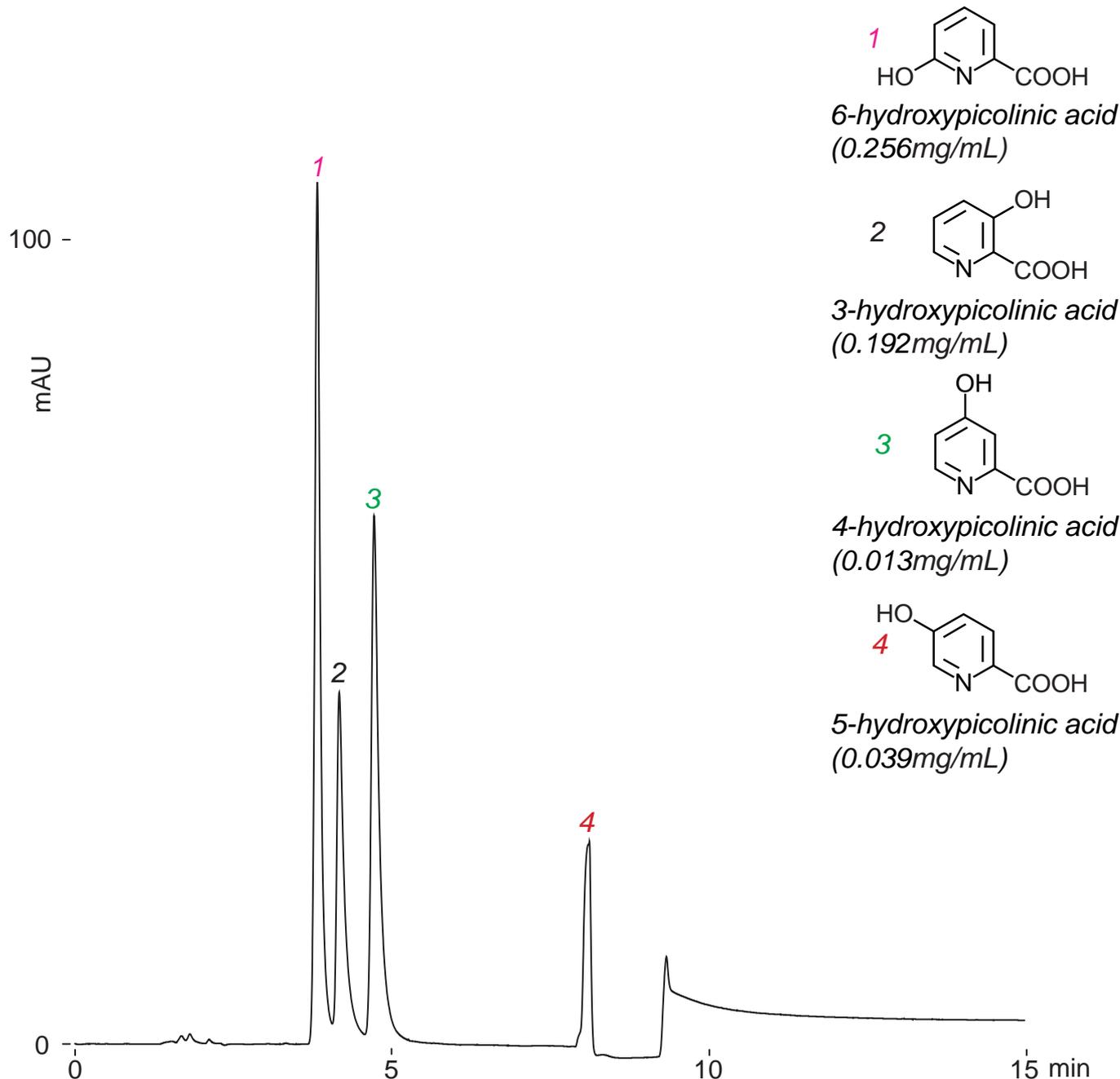

Unison UK-Amino, 150 x 4.6 mm

A: acetonitrile / 100mM ammonium acetate = 60 / 40

B: acetonitrile / water / acetic acid = 10 / 90 / 0.5

0%B (0-6min), 100%B (6-15min)

1.0mL/min (7MPa), 37°C, 260nm

10µL (0.013-0.256mg/mL, 100mM ammonium acetate / acetonitrile = 50 / 50)