

Unison UK-C18


Unison UK-C8

100 x 3 mm

HPLC COLUMN APPLICATION


LC-MS: Cedrol (Tricyclic sesquiterpene alcohol)

セドロール (三環式セスキテルペンアルコール) の LC-MS 分析

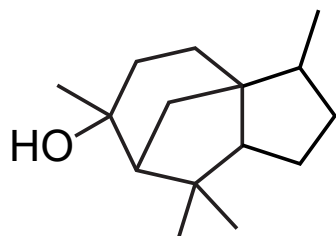

cedrol

Unison UK-C18

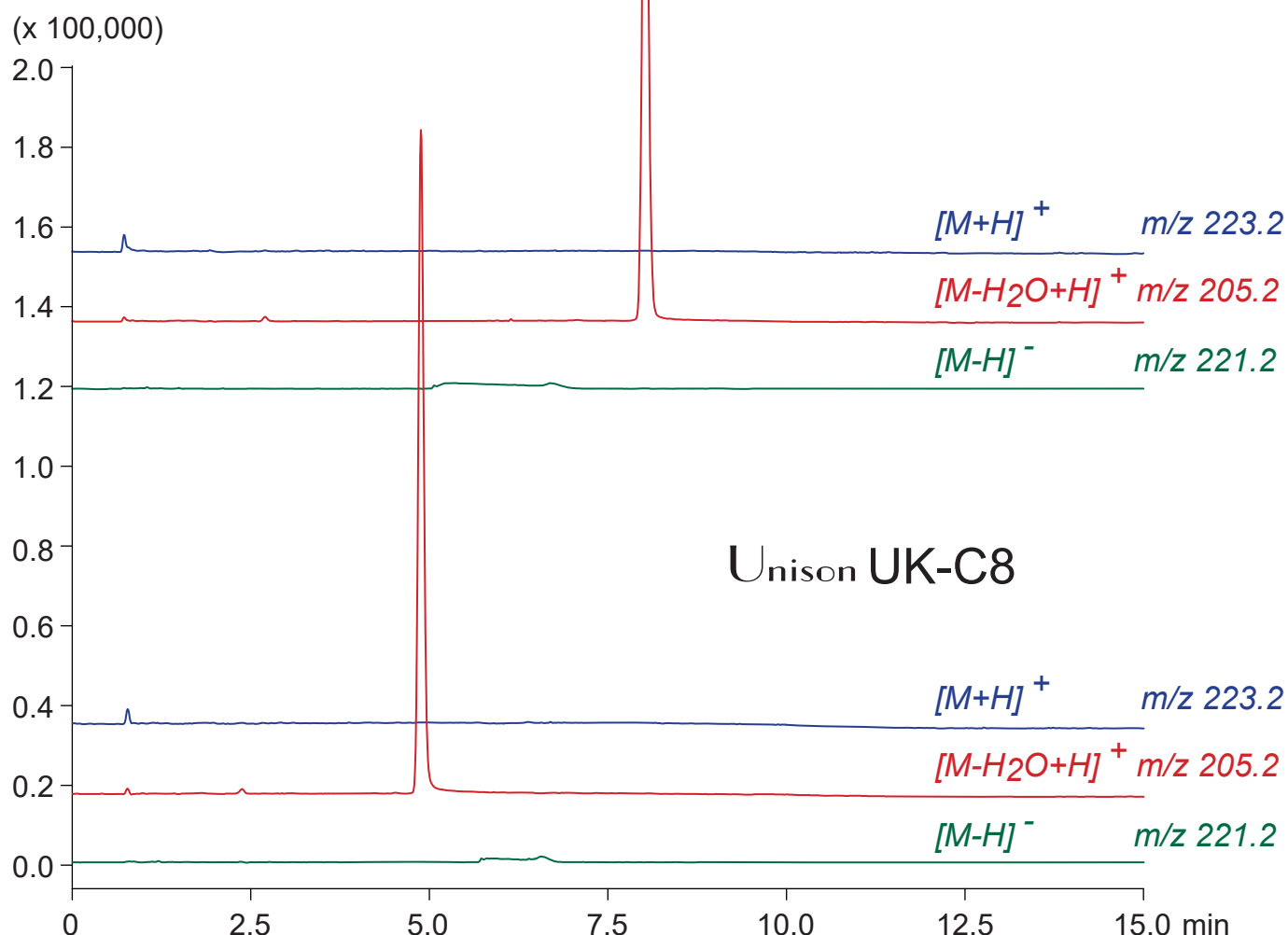

100 x 3 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

50-100%B (0-10min), 100%B (10-15min)

0.4mL/min (6MPa), 37°C, Single Quad. MS (ESI)

10μL (5μg/mL, ethanol / water / formic acid = 50 / 50 / 0.1)