


Unison UK-C18

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Cedrol (Tricyclic sesquiterpene alcohol)


セドロール（三環式セスキテルペンアルコール）のLC-MS分析

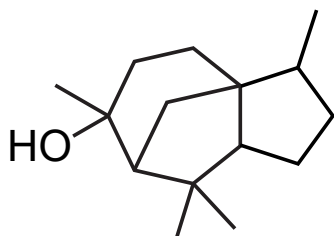

cedrol

A: ACN / water / HCOOH = 10 / 90 / 0.1
 B: ACN / HCOOH = 100 / 0.1
 10 μL (5 $\mu\text{g/mL}$, EtOH / water / HCOOH
 = 50 / 50 / 0.1)

(x 100,000)

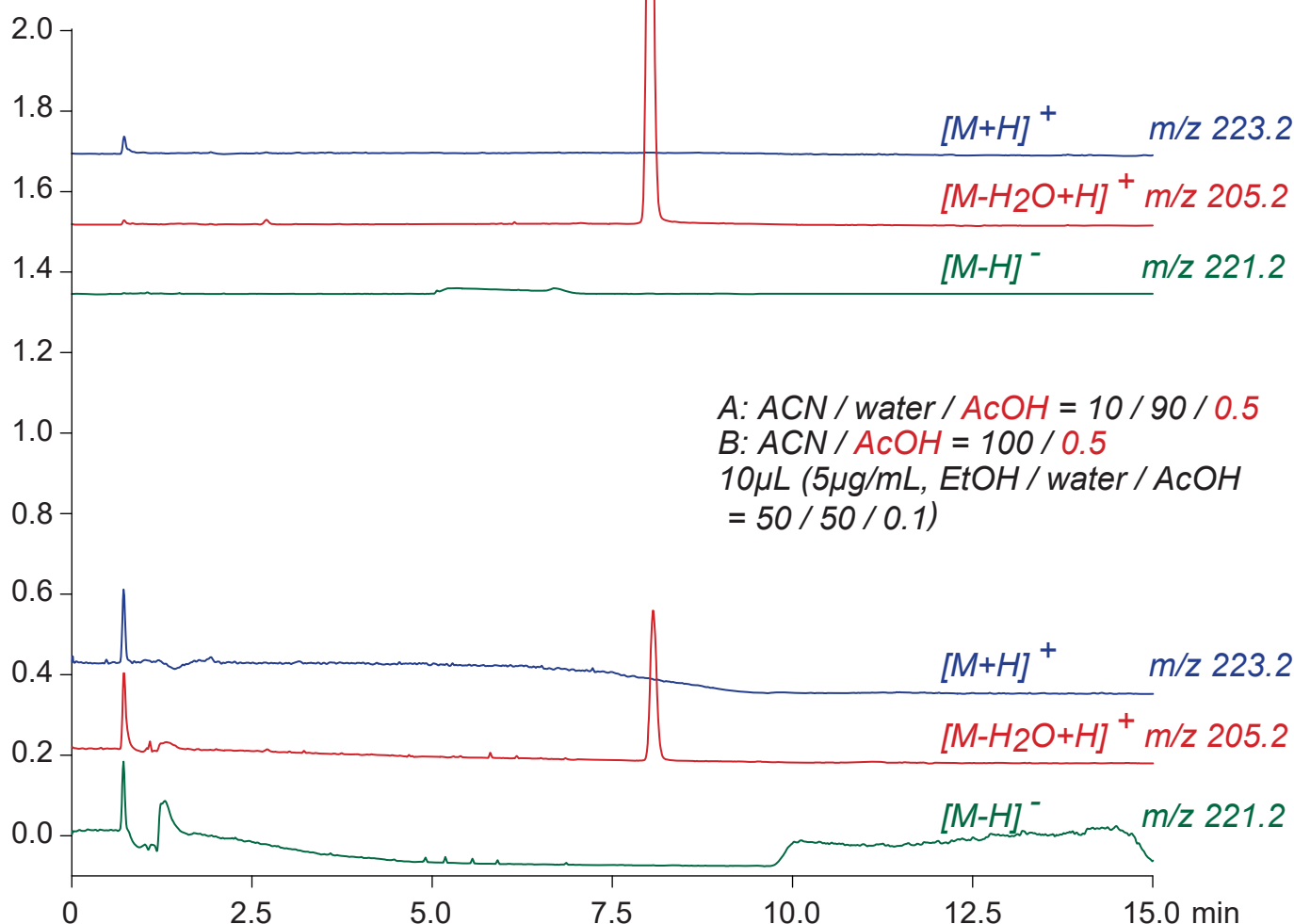

100 x 3 mm
 50-100%B (0-10min), 100%B (10-15min)
 0.4mL/min (6MPa), 37°C, Single Quad. MS (ESI)