

Unison UK-C8
Unison UK-C18


Cadenza CX-C18
Cadenza CW-C18
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Galantamine (Dementia therapeutic agent)

ガラントミン (認知症治療薬)

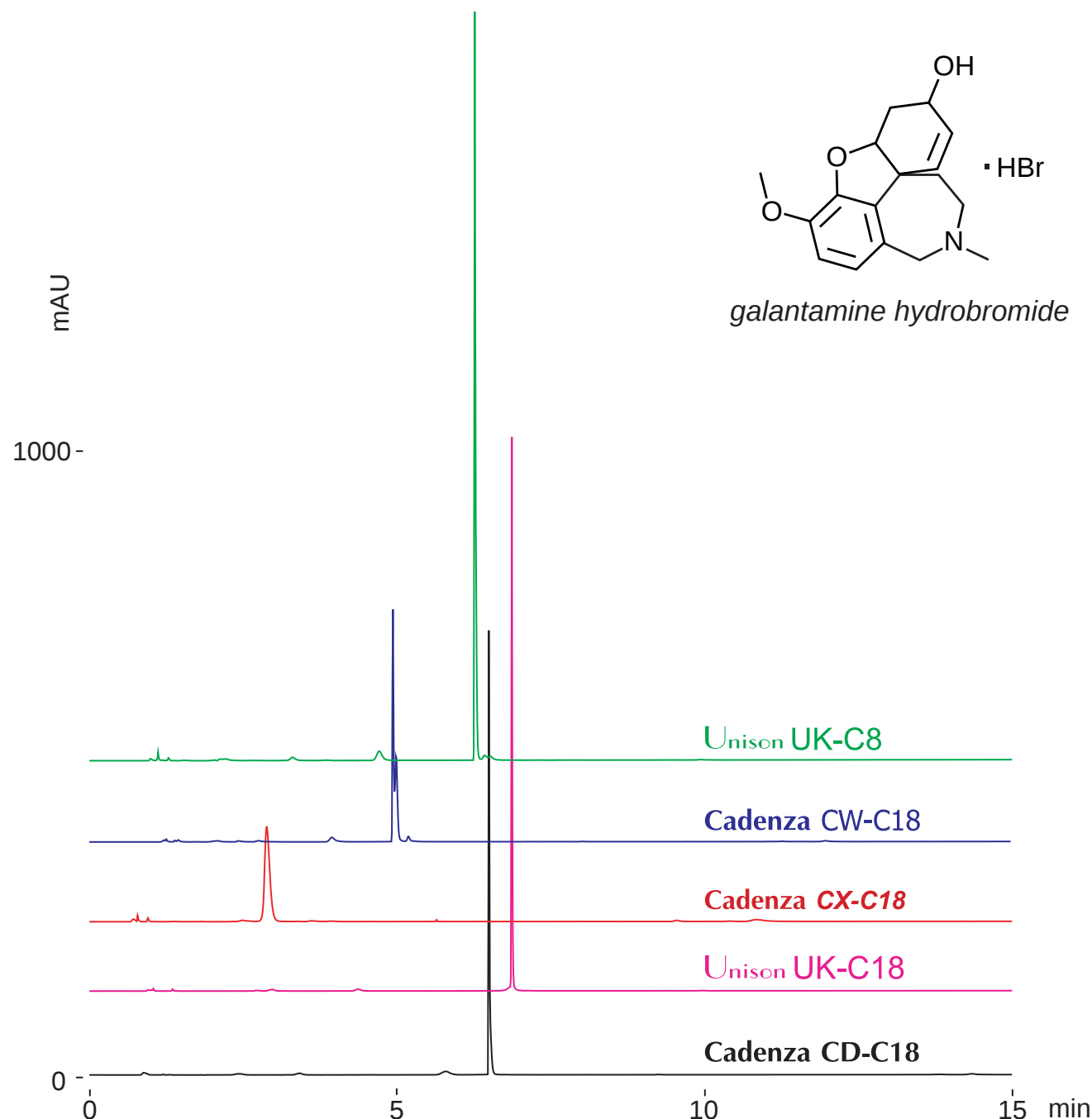

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / THF / formic acid = 90 / 10 / 0.5

0-50%B (0-15min)

1.0mL/min (10MPa), 37°C, 290nm

10μL (0.5mg/mL, water / formic acid = 100 / 1)