

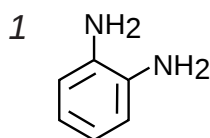
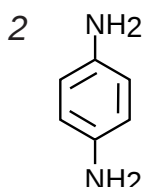
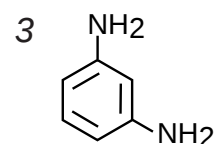
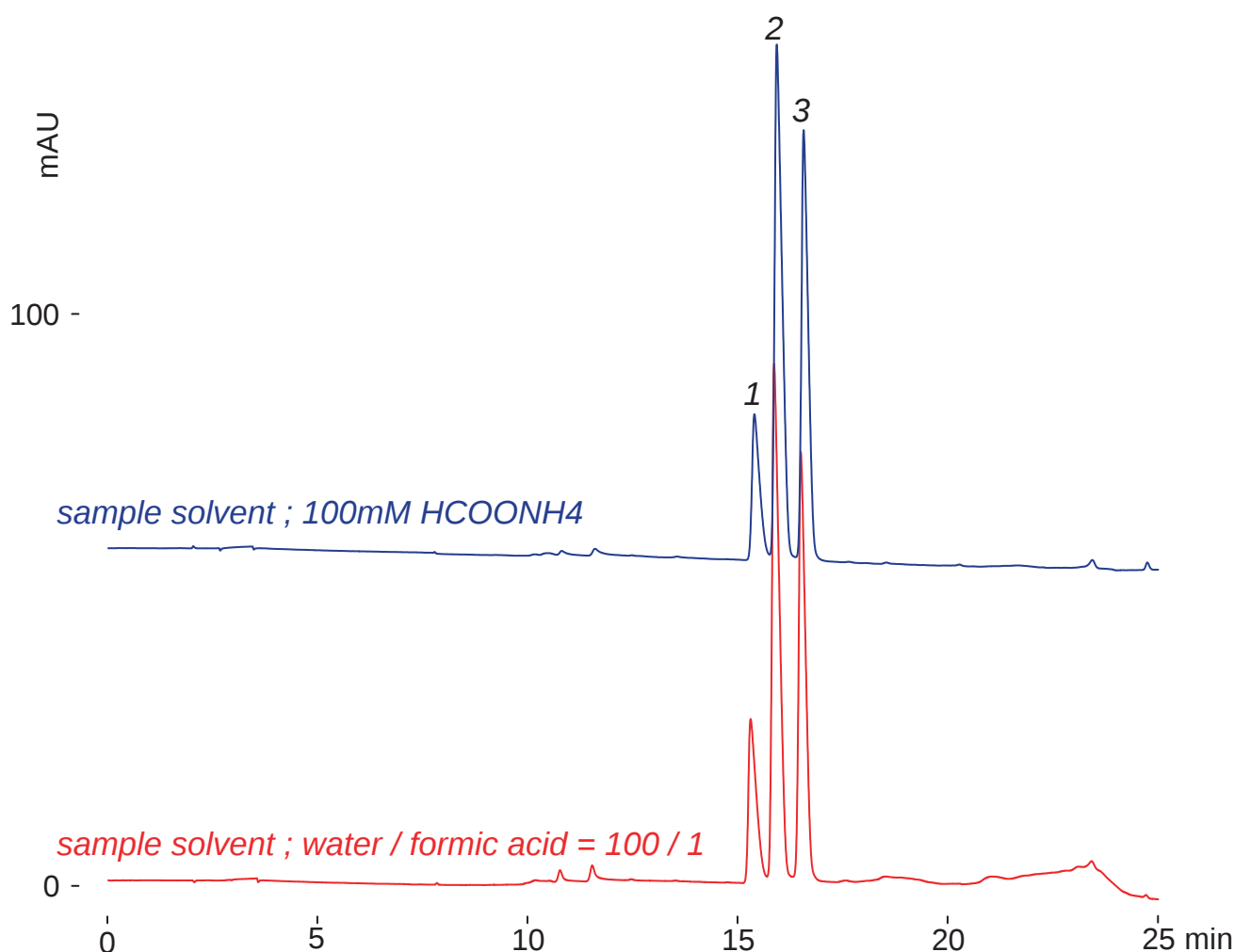
Scherzo SS-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Phenylenediamine isomers

フェニレンジアミン異性体


o-phenylenediamine

p-phenylenediamine

m-phenylenediamine


Scherzo SS-C18, 100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-100%B (0-25min)

1.0mL/min (12MPa), 37°C, 260nm

10μL (0.25mg/mL)