

Nardis ND-NX


Unison UK-Amino

250 x 4.6 mm

HPLC COLUMN APPLICATION

n-octyl-beta-D-thioglucofuranoside,
n-octyl beta-D-glucopyranoside(Membrane protein solubilizer)

n-オクチル-β-D-チオグルコピラノシド、n-オクチル-β-D-グルコピラノシド
(膜タンパク質可溶化剤)

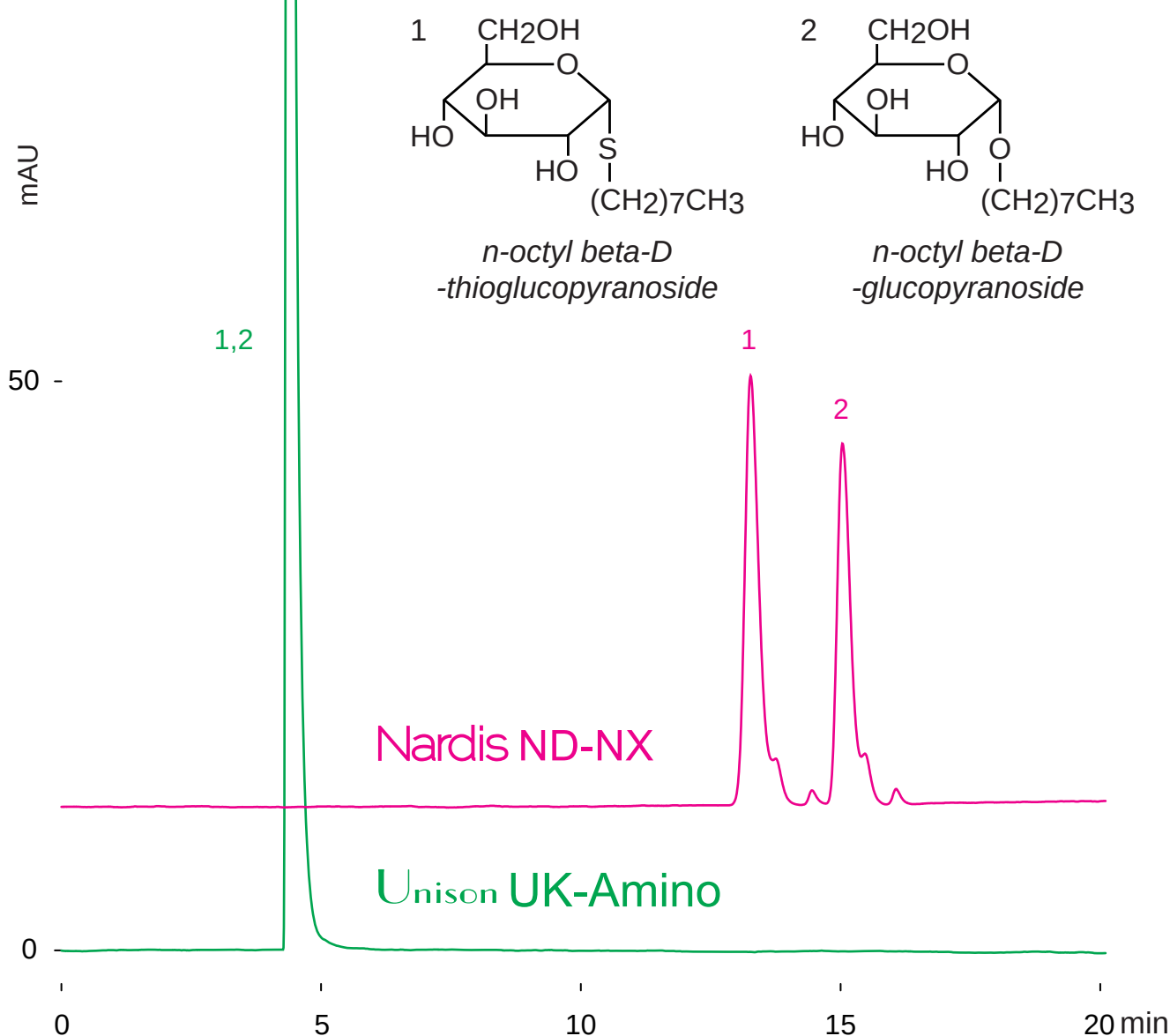

250 x 4.6 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium acetate = 50 / 50

10-20%B (0-20min)

1.0mL/min (8MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10μL (1mg/mL, water)